

(19) World Intellectual Property Organization
International Bureau



(43) International Publication Date
16 June 2011 (16.06.2011)

PCT

(10) International Publication Number
WO 2011/072168 A1

(51) International Patent Classification:
C12Q 1/68 (2006.01) *G06F 19/00* (2011.01)

(74) Agent: **KUHLMANN, Anthony C.**; Foley & Lardner LLP, 3579 Valley Centre Drive, Suite 300, San Diego, California 92130 (US).

(21) International Application Number:
PCT/US2010/059771

(81) Designated States (unless otherwise indicated, for every kind of national protection available): AE, AG, AL, AM, AO, AT, AU, AZ, BA, BB, BG, BH, BR, BW, BY, BZ, CA, CH, CL, CN, CO, CR, CU, CZ, DE, DK, DM, DO, DZ, EC, EE, EG, ES, FI, GB, GD, GE, GH, GM, GT, HN, HR, HU, ID, IL, IN, IS, JP, KE, KG, KM, KN, KP, KR, KZ, LA, LC, LK, LR, LS, LT, LU, LY, MA, MD, ME, MG, MK, MN, MW, MX, MY, MZ, NA, NG, NI, NO, NZ, OM, PE, PG, PH, PL, PT, RO, RS, RU, SC, SD, SE, SG, SK, SL, SM, ST, SV, SY, TH, TJ, TM, TN, TR, TT, TZ, UA, UG, US, UZ, VC, VN, ZA, ZM, ZW.

(22) International Filing Date:
9 December 2010 (09.12.2010)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:
61/285,938 11 December 2009 (11.12.2009) US
61/285,943 11 December 2009 (11.12.2009) US

(71) Applicant (for all designated States except US): **QUEST DIAGNOSTICS INVESTMENTS INCORPORATED** [US/US]; 300 Delaware Avenue, Wilmington, DE 19899 (US).

(84) Designated States (unless otherwise indicated, for every kind of regional protection available): ARIPO (BW, GH, GM, KE, LR, LS, MW, MZ, NA, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European (AL, AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI, FR, GB, GR, HR, HU, IE, IS, IT, LT, LU, LV, MC, MK, MT, NL, NO, PL, PT, RO, RS, SE, SI, SK, SM, TR), OAPI (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG).


(72) Inventors; and

(75) Inventors/Applicants (for US only): **HOLMQUIST, Brett** [US/US]; 21858 Bahamas, Mission Viejo, California 92692 (US). **CLARKE, Nigel J.** [GB/US]; 1689 Avenida Guillermo, Oceanside, California 92056 (US).

[Continued on next page]

(54) Title: MASS SPECTROMETRIC DETERMINATION OF COOKSON-DERIVATIZED, NON-METABOLIZED VITAMIN D

Figure 1

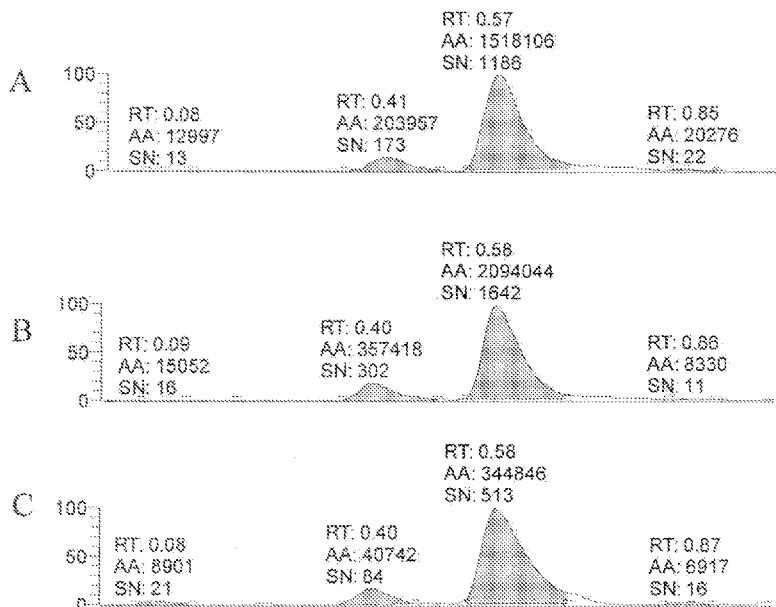

(57) Abstract: The invention relates to the detection of vitamin D. In a particular aspect, the invention relates to methods for detecting derivatized vitamin D by mass spectrometry.



WO 2011/072168 A1

Published:

— *with international search report (Art. 21(3))*

MASS SPECTROMETRIC DETERMINATION OF COOKSON-DERIVATIZED, NON-METABOLIZED VITAMIN D

FIELD OF THE INVENTION

[0001] The invention relates to the quantitative measurement of non-metabolized vitamin D. In a particular aspect, the invention relates to methods for quantitative measurement of non-metabolized vitamin D by tandem mass spectrometry.

BACKGROUND OF THE INVENTION

[0002] Vitamin D is an essential nutrient with important physiological roles in the positive regulation of calcium (Ca^{2+}) homeostasis. Vitamin D can be made *de novo* in the skin by exposure to sunlight or it can be absorbed from the diet. There are two forms of vitamin D; vitamin D₂ (ergocalciferol) and vitamin D₃ (cholecalciferol). Vitamin D₃ is the form synthesized *de novo* by animals. It is also a common supplement added to milk products and certain food products produced in the United States. Both dietary and intrinsically synthesized vitamin D₃ must undergo metabolic activation to generate the bioactive metabolites. In humans, the initial step of vitamin D₃ activation occurs primarily in the liver and involves hydroxylation to form the intermediate metabolite 25-hydroxycholecalciferol (calcifediol; 25OHD₃). Calcifediol is the major form of Vitamin D₃ in circulation. Circulating 25OHD₃ is then converted by the kidney to form 1,25-dihydroxyvitamin D₃ (calcitriol; 1,25(OH)₂D₃), which is generally believed to be the metabolite of Vitamin D₃ with the highest biological activity.

[0003] Vitamin D₂ is derived from fungal and plant sources. Many over-the-counter dietary supplements contain ergocalciferol (vitamin D₂) rather than cholecalciferol (vitamin D₃). Drisdol, the only high-potency prescription form of vitamin D available in the United States, is formulated with ergocalciferol. Vitamin D₂ undergoes a similar pathway of metabolic activation in humans as Vitamin D₃, forming the metabolites 25OHD₂ and 1,25(OH)₂D₂. Vitamin D₂ and vitamin D₃ have long been assumed to be biologically equivalent in humans, however recent reports suggest that there may be differences in the bioactivity and bioavailability of these two forms of vitamin D (Armas *et. al.*, (2004) J. Clin. Endocrinol. Metab. 89:5387-5391).

[0004] Measurement of vitamin D, the inactive vitamin D precursor, is rare in clinical settings. Rather, serum levels of 25-hydroxyvitamin D₃, 25-hydroxyvitamin D₂, and total 25-hydroxyvitamin D (“25OHD”) are useful indices of vitamin D nutritional status and the efficacy of certain vitamin D analogs. The measurement of 25OHD is commonly used in the diagnosis and management of disorders of calcium metabolism. In this respect, low levels of 25OHD are indicative of vitamin D deficiency associated with diseases such as hypocalcemia, hypophosphatemia, secondary hyperparathyroidism, elevated alkaline phosphatase, osteomalacia in adults and rickets in children. In patients suspected of vitamin D intoxication, elevated levels of 25OHD distinguishes this disorder from other disorders that cause hypercalcemia.

[0005] Measurement of 1,25(OH)₂D is also used in clinical settings. Certain disease states can be reflected by circulating levels of 1,25(OH)₂D, for example kidney disease and kidney failure often result in low levels of 1,25(OH)₂D. Elevated levels of 1,25(OH)₂D may be indicative of excess parathyroid hormone or can be indicative of certain diseases such as sarcoidosis or certain types of lymphomas.

[0006] Detection of vitamin D metabolites has been accomplished by radioimmunoassay with antibodies co-specific for 25OHD₂ and 25OHD₃. Because the current immunologically-based assays do not separately resolve 25OHD₂ and 25OHD₃, the source of any nutritional deficiency of vitamin D cannot be determined without resorting to other tests. Reports have been published that disclose methods for detecting specific vitamin D metabolites using mass spectrometry. In some of the reports, the vitamin D metabolites are derivatized prior to mass spectrometry, but in others, they are not. For example Holmquist, et al., U.S. Patent Application Serial No. 11/946765, filed December 28, 2007; Yeung B, *et al.*, J Chromatogr. 1993, 645(1):115-23; Higashi T, et al., Steroids. 2000, 65(5):281-94; Higashi T, *et al.*, Biol Pharm Bull. 2001, 24(7):738-43; Higashi T, *et al.*, J Pharm Biomed Anal. 2002, 29(5):947-55; Higashi T, et al., Anal. Biochem Chem, 2008, 391:229-38; and Aronov, et al., Anal Bioanal Chem, 2008, 391:1917-30 disclose methods for detecting various vitamin D metabolites by derivatizing the metabolites prior to mass spectrometry. Methods to detect underivatized vitamin D metabolites are reported in Clarke, *et al.*, in U.S. Patent Application Serial Nos. 11/101,166, filed April 6, 2005, and 11/386,215, filed March 21, 2006, and Singh, et al., in U.S. Patent Application Serial No. 10/977,121, filed October 24, 2004. Reports have also been published

that disclose derivatization of vitamin D₃ with Cookson-type reagents, specifically 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD) and 4-[2-(6,7-dimethoxy-4-methyl-3-oxo-3,4-dihydroquinoxalyl)ethyl]-1,2,4-triazoline-3,5-dione (DMEQ-TAD). See Aberhart, J, *et al.*, J. Org. Chem. 1976, 41(12):2098-2102, and Kamao, M, *et al.*, J Chromatogr. B 2007, 859:192-200.

SUMMARY OF THE INVENTION

[0007] The present invention provides methods for detecting the amount of one or more (non-metabolized) forms of vitamin D in a sample by mass spectrometry, including tandem mass spectrometry. In one aspect, vitamin D from the sample is derivatized prior to analysis by mass spectrometry. In a second aspect, vitamin D from the sample is not derivatized prior to analysis by mass spectrometry.

[0008] In some embodiments of the first aspect, the methods include the steps of: (i) subjecting Cookson-type-derivatized vitamin D in the sample to an ionization source under conditions suitable to generate one or more precursor ions detectable by mass spectrometry; (ii) fragmenting at least one of said precursor ions to generate one or more fragment ions detectable by mass spectrometry; (iii) determining the amount of one or more of the precursor and fragment ions by mass spectrometry; and (iv) relating the amount of ions determined in step (iii) to the amount of a vitamin D in the sample. In these methods, the sample is subjected to a Cookson-type derivatization reagent under conditions sufficient to generate Cookson-type-derivatized vitamin D prior to step (i). In some embodiments, the Cookson-type-derivatized-vitamin D is subjected to an extraction column and an analytical column prior to ionization. In related embodiments, the analytical column may be a high performance liquid chromatography (HPLC) column.

[0009] In some embodiments of the first aspect, the methods include the steps of: (i) subjecting the sample to turbulent flow liquid chromatography (TFLC); (ii) subjecting Cookson-type-derivatized vitamin D from the sample an ionization source under conditions suitable to generate one or more ions detectable by mass spectrometry; (iii) determining the amount of one or more of the Cookson-type-derivatized vitamin D ions by mass spectrometry; and (iv) relating the amount of Cookson-type-derivatized vitamin D ions determined in step (iii) to the amount of a vitamin D in the sample. In these embodiments, the sample is subjected to a Cookson-type derivatizing reagent under conditions sufficient to generate Cookson-type-derivatized vitamin D

in the sample prior to step (i). In some embodiments, the sample is subjected to high performance liquid chromatography (HPLC) after step (i) but prior to step (ii).

[0010] In some embodiments, the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD). In some of these embodiments, the one or more precursor ions comprise one or more ions selected from the group consisting of ions with a mass/charge ratio (m/z) of 572.4 ± 0.5 and 560.4 ± 0.5 . In some of these embodiments, one or more fragment ions comprise an ion with a mass/charge ratio (m/z) of 298.1 ± 0.5 .

[0011] In embodiments where vitamin D comprises vitamin D₂, one or more Cookson-type-derivatized vitamin D ions may comprise a precursor ion with a mass/charge ratio (m/z) of 572.3 ± 0.5 and a fragment ion with a mass/charge ratio (m/z) of 298.1 ± 0.5 . In embodiments where vitamin D comprises vitamin D₃, one or more Cookson-type-derivatized vitamin D ions may comprise a precursor ion with a mass/charge ratio (m/z) of 560.3 ± 0.5 and a fragment ion with a mass/charge ratio (m/z) of 298.1 ± 0.5 . In embodiments where vitamin D comprises vitamin D₂ and vitamin D₃, one or more precursor ions may comprise a vitamin D₂ precursor ion with a mass/charge ratio (m/z) of 572.4 ± 0.5 and a vitamin D₃ precursor ion with a m/z of 560.4 ± 0.5 ; and one or more fragment ions may comprise a vitamin D₂ fragment ion with a m/z of 298.1 ± 0.5 and a vitamin D₃ fragment ion with a m/z of 298.1 ± 0.5 .

[0012] In some embodiments of the second aspect, the invention provides methods for determining the amount of vitamin D in a sample by tandem mass spectrometry without ionization and detection of derivatized vitamin D. In some embodiments, the vitamin D comprises vitamin D₂; these methods include the steps of: (i) subjecting vitamin D₂ from a sample to an ionization source under conditions suitable to generate one or more precursor ions detectable by mass spectrometry selected from the group consisting of ions with a mass to charge ratio (m/z) of 397.2 ± 0.5 or 379.2 ± 0.5 ; (ii) fragmenting at least one of said precursor ions to generate one or more fragment ions detectable by mass spectrometry; (iii) determining the amount of one or more of the ions generated in steps (i) and (ii) by mass spectrometry; and (iv) relating the presence of vitamin D₂ ions determined in step (iii) to the presence of vitamin D₂ in the sample. In these methods, if the fragmented precursor ions comprise an ion with m/z of 397.2 ± 0.5 , the fragment ions comprise one or more ions selected from the group consisting of ions with m/z of 159.0 ± 0.5 , 146.9 ± 0.5 , 133.1 ± 0.5 , and 121.0 ± 0.5 . In these methods, if the

fragmented precursor ions comprise an ion with m/z of 379.2 ± 0.5 , the fragment ions comprise one or more ions selected from the group consisting of ions with m/z of 283.2 ± 0.5 , 187.3 ± 0.5 , 175.2 ± 0.5 , and 159.0 ± 0.5 . In some embodiments, vitamin D₂ from the sample may be subjected to an extraction column, such as a solid phase extraction (SPE) column or a turbulent flow liquid chromatography (THLC) column, prior to ionization. In some embodiments, vitamin D₂ from the sample is further subjected to an analytical column, such as a high performance liquid chromatography (HPLC) column, prior to ionization.

[0013] In other embodiments of the second aspect, vitamin D comprises vitamin D₃; these methods include the steps of: (i) subjecting vitamin D₃ from a sample to an ionization source under conditions suitable to generate one or more precursor ions detectable by mass spectrometry selected from the group consisting of ions with a mass to charge ratio (m/z) of 385.2 ± 0.5 or 367.2 ± 0.5 ; (ii) fragmenting at least one of the precursor ions to generate one or more fragment ions detectable by mass spectrometry; (iii) determining the amount of one or more of the ions generated in steps (i) and (ii) by mass spectrometry; and (iv) relating the presence of vitamin D₃ ions determined in step (iii) to the presence of vitamin D₃ in the sample. In these methods, if the fragmented precursor ions comprise an ion with m/z of 385.2 ± 0.5 , the fragment ions comprise one or more ions selected from the group consisting of ions with m/z of 159.0 ± 0.5 , 147.0 ± 0.5 , 133.1 ± 0.5 , and 107.1 ± 0.5 . If the fragmented precursor ions comprise an ion with m/z of 367.2 ± 0.5 , the fragment ions comprise one or more ions selected from the group consisting of ions with m/z of 172.2 ± 0.5 , 145.0 ± 0.5 , and 119.1 ± 0.5 . In some embodiments, vitamin D₃ from the sample may be subjected to an extraction column, such as a solid phase extraction (SPE) column or a turbulent flow liquid chromatography (THLC) column, prior to ionization. In some embodiments, vitamin D₃ from the sample is further subjected to an analytical column, such as a high performance liquid chromatography (HPLC) column, prior to ionization.

[0014] In the methods described herein, mass spectrometry may be tandem mass spectrometry. In embodiments utilizing tandem mass spectrometry, tandem mass spectrometry may be conducted by any method known in the art, including for example, multiple reaction monitoring, precursor ion scanning, or product ion scanning.

[0015] In embodiments which utilize an extraction column, the extraction column may be a solid phase extraction (SPE) column; such as a turbulent flow liquid chromatography (TFLC) column. In some embodiments, which utilize two or more of an extraction column, an analytical column, and an ionization source, the utilized components may be connected in an on-line fashion to allow for automated sample processing and analysis.

[0016] Cookson-type derivatizing reagents useful for certain embodiments may be selected from the group consisting of 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD), 4-methyl-1,2,4-triazoline-3,5-dione (MTAD), 4-[2-(6,7-dimethoxy-4-methyl-3-oxo-3,4-dihydroquinoxalyl)ethyl]-1,2,4-triazoline-3,5-dione (DMEQTAD), 4-(4-nitrophenyl)-1,2,4-triazoline-3,5-dione (NPTAD), and 4-ferrocenylmethyl-1,2,4-triazoline-3,5-dione (FMTAD), and isotopically labeled variants thereof. In preferred embodiments, the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD) or an isotopically labeled variant thereof. In certain preferred embodiments, the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD) or an isotopically labeled variant thereof.

[0017] As used herein, the term “vitamin D” refers to any one or more naturally occurring or synthetic analog of vitamin D which are non-metabolized. This is in contrast to vitamin D metabolites which are identified by specific chemical modification (e.g., 25-hydroxyvitamin D and 1 α ,25-dihydroxyvitamin D) which occurs during metabolism. Non-metabolized vitamin D may also be referred to as “nutritional” vitamin D to distinguish from metabolized forms. Reference to vitamin D without specifying a metabolized form is a reference to non-metabolized forms.

[0018] As used herein, “derivatizing” means reacting two molecules to form a new molecule. Thus, a derivatizing agent is an agent that may be reacted with another substance to derivatize the substance. For example, 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD) is a derivatizing reagent that may be reacted with vitamin D to form a PTAD-derivatized vitamin D.

[0019] As used here, the names of derivatized forms vitamin D include an indication as to the nature of derivatization. For example, the PTAD derivative of vitamin D₂ is indicated as PTAD-vitamin D₂ (or PTAD-derivatized vitamin D₂).

[0020] As used herein, a “Cookson-type derivatizing agent” is a 4-substituted 1,2,4-triazoline-3,5-dione compound. Exemplary Cookson-type derivatizing agents include 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD), 4-methyl-1,2,4-triazoline-3,5-dione (MTAD), 4-[2-(6,7-dimethoxy-4-methyl-3-oxo-3,4-dihydroquinoxalyl)ethyl]-1,2,4-triazoline-3,5-dione (DMEQTAD), 4-(4-nitrophenyl)-1,2,4-triazoline-3,5-dione (NPTAD), and 4-ferrocenylmethyl-1,2,4-triazoline-3,5-dione (FMTAD). Additionally, isotopically labeled variants of Cookson-type derivatizing agents may be used in some embodiments. For example, the $^{13}\text{C}_6$ -PTAD isotopic variant is 6 mass units heavier than normal PTAD and may be used in some embodiments. Derivatization of vitamin D metabolites by Cookson-type reagents can be conducted by any appropriate method. See, e.g., Holmquist, et al., U.S. Patent Application Serial No. 11/946765, filed December 28, 2007; Yeung B, *et al.*, J Chromatogr. 1993, 645(1):115-23; Higashi T, et al., Steroids. 2000, 65(5):281-94; Higashi T, *et al.*, Biol Pharm Bull. 2001, 24(7):738-43; Higashi T, *et al.*, J Pharm Biomed Anal. 2002, 29(5):947-55; Higashi T, et al., Anal. Bioanal Chem, 2008, 391:229-38; and Aronov, et al., Anal Bioanal Chem, 2008, 391:1917-30.

[0021] Vitamin D may refer to one or more forms of vitamin D, such as vitamin D₂ and/or vitamin D₃. In embodiments where a sample comprises a plurality of vitamin D forms, the plurality of vitamin D forms may be ionized simultaneously. For example, in some embodiments, the amounts of vitamin D₂ and vitamin D₃ are determined in the same sample. In these embodiments, (derivatized or underivatized) vitamin D₂ and vitamin D₃ may be ionized simultaneously.

[0022] The term “simultaneous” as applied to simultaneously detecting the amount of two or more analytes from a sample means acquiring data reflective of the amount of the two or more analytes in the sample from the same sample injection. The data for each analyte may be acquired sequentially or in parallel, depending on the instrumental techniques employed. For example, a single sample containing two analytes may be injected into a HPLC column, which may then elute each analyte one after the other, resulting in introduction of the analytes into a mass spectrometer sequentially. Determining the amount of each of these two analytes is simultaneous for the purposes herein, as both analytes result from the same sample injection into the HPLC.

[0023] Vitamin D may be found in the circulation of an animal and/or may be generated by a biological organism, such as an animal. As such, samples may be obtained, for example, from a patient; that is, a living person, male or female, presenting oneself in a clinical setting for diagnosis, prognosis, or treatment of a disease or condition. Preferred samples may be biological samples; particularly biological fluid samples such as serum or plasma. The methods presented herein may be used to determine the amount of vitamin D present in a sample when taken from the human.

[0024] In certain preferred embodiments of the methods disclosed herein, mass spectrometry is performed in positive ion mode. Alternatively, mass spectrometry is performed in negative ion mode. Various ionization sources, including for example atmospheric pressure chemical ionization (APCI), laser diode thermal desorption (LDTD), or electrospray ionization (ESI), may be used in embodiments of the present invention. In certain preferred embodiments, vitamin D metabolites are measured using APCI or LDTD in positive ion mode.

[0025] In preferred embodiments, one or more separately detectable internal standards are provided in the sample, the amount of which are also determined in the sample. In these embodiments, all or a portion of both the analyte(s) of interest and the internal standard(s) present in the sample are ionized to produce a plurality of ions detectable in a mass spectrometer, and one or more ions produced from each are detected by mass spectrometry. Preferably, the internal standard(s) are one or more of vitamin D₂-[6, 19, 19]-²H₃, vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆, vitamin D₃-[6, 19, 19]-²H₃, and vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆.

[0026] One or more separately detectable internal standards may be provided in the sample prior to treatment with a Cookson-type derivatizing reagent (if applicable), or any purification of the analyte(s) from the sample. In these embodiments, the one or more internal standards may undergo derivatization and/or purification along with the endogenous vitamin D, in which case ions of the derivatized and/or purified internal standards are detected by mass spectrometry. In these embodiments, the presence or amount of ions generated from the analyte of interest may be related to the presence or amount of analyte of interest in the sample. In some embodiments, the internal standards may be isotopically labeled versions of vitamin D, such as vitamin D₂-[6, 19, 19]-²H₃, vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆, vitamin D₃-[6, 19, 19]-²H₃, and vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆.

[0027] Ions detectable in a mass spectrometer may be generated for each of the exemplary internal standards listed above. Exemplary spectra generated for several exemplary internal standards are discussed in Examples 8 and 9, and shown in Figures 6-7, 9-10, 12-13, and 15-16.

[0028] As used herein, an “isotopic label” produces a mass shift in the labeled molecule relative to the unlabeled molecule when analyzed by mass spectrometric techniques. Examples of suitable labels include deuterium (^2H), ^{13}C , and ^{15}N . For example, vitamin D_2 -[6, 19, 19]- $^2\text{H}_3$ and vitamin D_3 -[6, 19, 19]- $^2\text{H}_3$ has masses of about 3 mass units higher than vitamin D_2 and vitamin D_3 , respectively. The isotopic label can be incorporated at one or more positions in the molecule and one or more kinds of isotopic labels can be used on the same isotopically labeled molecule.

[0029] In other embodiments, the amount of vitamin D ion or ions may be determined by comparison to one or more external reference standards. Exemplary external reference standards include blank plasma or serum spiked with one or more of vitamin D_2 -[6, 19, 19]- $^2\text{H}_3$, vitamin D_2 -[26, 26, 26, 27, 27, 27]- $^2\text{H}_6$, vitamin D_3 -[6, 19, 19]- $^2\text{H}_3$, and vitamin D_3 -[26, 26, 26, 27, 27, 27]- $^2\text{H}_6$. External standards typically will undergo the same treatment and analysis as any other sample to be analyzed, including treatment with one or more Cookson-type reagents prior to mass spectrometry in embodiments wherein derivatized vitamin D is detected.

[0030] In certain embodiments, the lower limit of quantitation (LLOQ) of vitamin D_2 and vitamin D_3 are less than 10 ng/mL; preferably less than 5 ng/mL; preferably less than 2 ng/mL.

[0031] As used herein, unless otherwise stated, the singular forms “a,” “an,” and “the” include plural reference. Thus, for example, a reference to “a protein” includes a plurality of protein molecules.

[0032] As used herein, the term “purification” or “purifying” does not refer to removing all materials from the sample other than the analyte(s) of interest. Instead, purification refers to a procedure that enriches the amount of one or more analytes of interest relative to other components in the sample that may interfere with detection of the analyte of interest. Purification of the sample by various means may allow relative reduction of one or more interfering substances, e.g., one or more substances that may or may not interfere with the detection of selected parent or daughter ions by mass spectrometry. Relative reduction as this

term is used does not require that any substance, present with the analyte of interest in the material to be purified, is entirely removed by purification.

[0033] As used herein, the term “solid phase extraction” or “SPE” refers to a process in which a chemical mixture is separated into components as a result of the affinity of components dissolved or suspended in a solution (i.e., mobile phase) for a solid through or around which the solution is passed (i.e., solid phase). In some instances, as the mobile phase passes through or around the solid phase, undesired components of the mobile phase may be retained by the solid phase resulting in a purification of the analyte in the mobile phase. In other instances, the analyte may be retained by the solid phase, allowing undesired components of the mobile phase to pass through or around the solid phase. In these instances, a second mobile phase is then used to elute the retained analyte off of the solid phase for further processing or analysis. SPE, including TFLC, may operate via a unitary or mixed mode mechanism. Mixed mode mechanisms utilize ion exchange and hydrophobic retention in the same column; for example, the solid phase of a mixed-mode SPE column may exhibit strong anion exchange and hydrophobic retention; or may exhibit column exhibit strong cation exchange and hydrophobic retention.

[0034] As used herein, the term “chromatography” refers to a process in which a chemical mixture carried by a liquid or gas is separated into components as a result of differential distribution of the chemical entities as they flow around or over a stationary liquid or solid phase.

[0035] As used herein, the term “liquid chromatography” or “LC” means a process of selective retardation of one or more components of a fluid solution as the fluid uniformly percolates through a column of a finely divided substance, or through capillary passageways. The retardation results from the distribution of the components of the mixture between one or more stationary phases and the bulk fluid, (i.e., mobile phase), as this fluid moves relative to the stationary phase(s). Examples of “liquid chromatography” include reverse phase liquid chromatography (RPLC), high performance liquid chromatography (HPLC), and turbulent flow liquid chromatography (TFLC) (sometimes known as high turbulence liquid chromatography (HTLC) or high throughput liquid chromatography).

[0036] As used herein, the term “high performance liquid chromatography” or “HPLC” (sometimes known as “high pressure liquid chromatography”) refers to liquid chromatography in

which the degree of separation is increased by forcing the mobile phase under pressure through a stationary phase, typically a densely packed column.

[0037] As used herein, the term “turbulent flow liquid chromatography” or “TFLC” (sometimes known as high turbulence liquid chromatography or high throughput liquid chromatography) refers to a form of chromatography that utilizes turbulent flow of the material being assayed through the column packing as the basis for performing the separation. TFLC has been applied in the preparation of samples containing two unnamed drugs prior to analysis by mass spectrometry. See, e.g., Zimmer *et al.*, *J Chromatogr A* 854: 23-35 (1999); see also, U.S. Patents No. 5,968,367, 5,919,368, 5,795,469, and 5,772,874, which further explain TFLC. Persons of ordinary skill in the art understand “turbulent flow”. When fluid flows slowly and smoothly, the flow is called “laminar flow”. For example, fluid moving through an HPLC column at low flow rates is laminar. In laminar flow the motion of the particles of fluid is orderly with particles moving generally in straight lines. At faster velocities, the inertia of the water overcomes fluid frictional forces and turbulent flow results. Fluid not in contact with the irregular boundary “outruns” that which is slowed by friction or deflected by an uneven surface. When a fluid is flowing turbulently, it flows in eddies and whirls (or vortices), with more “drag” than when the flow is laminar. Many references are available for assisting in determining when fluid flow is laminar or turbulent (e.g., Turbulent Flow Analysis: Measurement and Prediction, P.S. Bernard & J.M. Wallace, John Wiley & Sons, Inc., (2000); An Introduction to Turbulent Flow, Jean Mathieu & Julian Scott, Cambridge University Press (2001)).

[0038] As used herein, the term “gas chromatography” or “GC” refers to chromatography in which the sample mixture is vaporized and injected into a stream of carrier gas (as nitrogen or helium) moving through a column containing a stationary phase composed of a liquid or a particulate solid and is separated into its component compounds according to the affinity of the compounds for the stationary phase.

[0039] As used herein, the term “large particle column” or “extraction column” refers to a chromatography column containing an average particle diameter greater than about 50 μm .

[0040] As used herein, the term “analytical column” refers to a chromatography column having sufficient chromatographic plates to effect a separation of materials in a sample that elute from the column sufficient to allow a determination of the presence or amount of an analyte. In

a preferred embodiment the analytical column contains particles of about 5 μm in diameter. Such columns are often distinguished from “extraction columns”, which have the general purpose of separating or extracting retained material from non-retained materials in order to obtain a purified sample for further analysis.

[0041] As used herein, the terms “on-line” and “inline”, for example as used in “on-line automated fashion” or “on-line extraction” refers to a procedure performed without the need for operator intervention. In contrast, the term “off-line” as used herein refers to a procedure requiring manual intervention of an operator. Thus, if samples are subjected to precipitation, and the supernatants are then manually loaded into an autosampler, the precipitation and loading steps are off-line from the subsequent steps. In various embodiments of the methods, one or more steps may be performed in an on-line automated fashion.

[0042] As used herein, the term “mass spectrometry” or “MS” refers to an analytical technique to identify compounds by their mass. MS refers to methods of filtering, detecting, and measuring ions based on their mass-to-charge ratio, or “m/z”. MS technology generally includes (1) ionizing the compounds to form charged compounds; and (2) detecting the molecular weight of the charged compounds and calculating a mass-to-charge ratio. The compounds may be ionized and detected by any suitable means. A “mass spectrometer” generally includes an ionizer and an ion detector. In general, one or more molecules of interest are ionized, and the ions are subsequently introduced into a mass spectrometric instrument where, due to a combination of magnetic and electric fields, the ions follow a path in space that is dependent upon mass (“m”) and charge (“z”). See, e.g., U.S. Patent Nos. 6,204,500, entitled “Mass Spectrometry From Surfaces;” 6,107,623, entitled “Methods and Apparatus for Tandem Mass Spectrometry;” 6,268,144, entitled “DNA Diagnostics Based On Mass Spectrometry;” 6,124,137, entitled “Surface-Enhanced Photolabile Attachment And Release For Desorption And Detection Of Analytes;” Wright *et al.*, *Prostate Cancer and Prostatic Diseases* 1999, 2: 264-76; and Merchant and Weinberger, *Electrophoresis* 2000, 21: 1164-67.

[0043] As used herein, the term “operating in negative ion mode” refers to those mass spectrometry methods where negative ions are generated and detected. The term “operating in positive ion mode” as used herein, refers to those mass spectrometry methods where positive ions are generated and detected.

[0044] As used herein, the term “ionization” or “ionizing” refers to the process of generating an analyte ion having a net electrical charge equal to one or more electron units. Negative ions are those having a net negative charge of one or more electron units, while positive ions are those having a net positive charge of one or more electron units.

[0045] As used herein, the term “electron ionization” or “EI” refers to methods in which an analyte of interest in a gaseous or vapor phase interacts with a flow of electrons. Impact of the electrons with the analyte produces analyte ions, which may then be subjected to a mass spectrometry technique.

[0046] As used herein, the term “chemical ionization” or “CI” refers to methods in which a reagent gas (e.g. ammonia) is subjected to electron impact, and analyte ions are formed by the interaction of reagent gas ions and analyte molecules.

[0047] As used herein, the term “fast atom bombardment” or “FAB” refers to methods in which a beam of high energy atoms (often Xe or Ar) impacts a non-volatile sample, desorbing and ionizing molecules contained in the sample. Test samples are dissolved in a viscous liquid matrix such as glycerol, thioglycerol, m-nitrobenzyl alcohol, 18-crown-6 crown ether, 2-nitrophenyloctyl ether, sulfolane, diethanolamine, and triethanolamine. The choice of an appropriate matrix for a compound or sample is an empirical process.

[0048] As used herein, the term “matrix-assisted laser desorption ionization” or “MALDI” refers to methods in which a non-volatile sample is exposed to laser irradiation, which desorbs and ionizes analytes in the sample by various ionization pathways, including photo-ionization, protonation, deprotonation, and cluster decay. For MALDI, the sample is mixed with an energy-absorbing matrix, which facilitates desorption of analyte molecules.

[0049] As used herein, the term “surface enhanced laser desorption ionization” or “SELDI” refers to another method in which a non-volatile sample is exposed to laser irradiation, which desorbs and ionizes analytes in the sample by various ionization pathways, including photo-ionization, protonation, deprotonation, and cluster decay. For SELDI, the sample is typically bound to a surface that preferentially retains one or more analytes of interest. As in MALDI, this process may also employ an energy-absorbing material to facilitate ionization.

[0050] As used herein, the term “electrospray ionization” or “ESI,” refers to methods in which a solution is passed along a short length of capillary tube, to the end of which is applied a high positive or negative electric potential. Solution reaching the end of the tube is vaporized (nebulized) into a jet or spray of very small droplets of solution in solvent vapor. This mist of droplets flows through an evaporation chamber, which is heated slightly to prevent condensation and to evaporate solvent. As the droplets get smaller the electrical surface charge density increases until such time that the natural repulsion between like charges causes ions as well as neutral molecules to be released.

[0051] As used herein, the term “atmospheric pressure chemical ionization” or “APCI,” refers to mass spectrometry methods that are similar to ESI; however, APCI produces ions by ion-molecule reactions that occur within a plasma at atmospheric pressure. The plasma is maintained by an electric discharge between the spray capillary and a counter electrode. Then ions are typically extracted into the mass analyzer by use of a set of differentially pumped skimmer stages. A counterflow of dry and preheated N₂ gas may be used to improve removal of solvent. The gas-phase ionization in APCI can be more effective than ESI for analyzing less-polar species.

[0052] The term “atmospheric pressure photoionization” or “APPI” as used herein refers to a form of ionization where the mechanism for the ionization of molecule M is photon absorption and electron ejection to form the molecular ion M⁺. Because the energy of the absorbed photon is typically just above the ionization potential, the molecular ion is less susceptible to dissociation. In many cases it may be possible to analyze samples without the need for chromatography, thus saving significant time and expense. In the presence of water vapor or protic solvents, the molecular ion can extract H to form MH⁺. This tends to occur if M has a high proton affinity. This does not affect quantitation accuracy because the sum of M⁺ and MH⁺ is constant. Drug compounds in protic solvents are usually observed as MH⁺, whereas nonpolar compounds such as naphthalene or testosterone usually form M⁺. See, e.g., Robb *et al.*, *Anal. Chem.* 2000, 72(15): 3653-3659.

[0053] As used herein, the term “inductively coupled plasma” or “ICP” refers to methods in which a sample interacts with a partially ionized gas at a sufficiently high temperature such that most elements are atomized and ionized.

[0054] As used herein, the term “desorption” refers to the removal of an analyte from a surface and/or the entry of an analyte into a gaseous phase. Laser diode thermal desorption (LDTD) is a technique wherein a sample containing the analyte is thermally desorbed into the gas phase by a laser pulse. The laser hits the back of a specially made 96-well plate with a metal base. The laser pulse heats the base and the heat causes the sample to transfer into the gas phase. The gas phase sample may then be drawn into an ionization source, where the gas phase sample is ionized in preparation for analysis in the mass spectrometer. When using LDTD, ionization of the gas phase sample may be accomplished by any suitable technique known in the art, such as by ionization with a corona discharge (for example by APCI).

[0055] As used herein, the term “field desorption” refers to methods in which a non-volatile test sample is placed on an ionization surface, and an intense electric field is used to generate analyte ions.

[0056] As used herein, the term “selective ion monitoring” is a detection mode for a mass spectrometric instrument in which only ions within a relatively narrow mass range, typically about one mass unit, are detected.

[0057] As used herein, “multiple reaction mode,” sometimes known as “selected reaction monitoring,” is a detection mode for a mass spectrometric instrument in which a precursor ion and one or more fragment ions are selectively detected.

[0058] As used herein, the term “lower limit of quantification”, “lower limit of quantitation” or “LLOQ” refers to the point where measurements become quantitatively meaningful. The analyte response at this LOQ is identifiable, discrete and reproducible with a relative standard deviation (RSD %) of less than 20% and an accuracy of 80% to 120%.

[0059] As used herein, the term “limit of detection” or “LOD” is the point at which the measured value is larger than the uncertainty associated with it. The LOD is the point at which a value is beyond the uncertainty associated with its measurement and is defined as three times the RSD of the mean at the zero concentration.

[0060] As used herein, an “amount” of an analyte in a body fluid sample refers generally to an absolute value reflecting the mass of the analyte detectable in volume of sample. However, an amount also contemplates a relative amount in comparison to another analyte amount. For

example, an amount of an analyte in a sample can be an amount which is greater than a control or normal level of the analyte normally present in the sample.

[0061] The term “about” as used herein in reference to quantitative measurements not including the measurement of the mass of an ion, refers to the indicated value plus or minus 10%. Mass spectrometry instruments can vary slightly in determining the mass of a given analyte. The term “about” in the context of the mass of an ion or the mass/charge ratio of an ion refers to +/- 0.50 atomic mass unit.

[0062] The summary of the invention described above is non-limiting and other features and advantages of the invention will be apparent from the following detailed description of the invention, and from the claims.

BRIEF DESCRIPTION OF THE DRAWINGS

[0063] Figures 1A-C show exemplary chromatograms for PTAD-vitamin D₂, PTAD-vitamin D₃, and PTAD-vitamin D₃-[6, 19, 19]-²H₃ (internal standard), respectively. Details are discussed in Example 3.

[0064] Figures 2A and 2B show exemplary calibration curves for vitamin D₂ and vitamin D₃ in serum samples determined by methods described in Example 3.

[0065] Figure 3 shows a plots of coefficient of variation versus concentration for vitamin D₂ and vitamin D₃. Details are described in Example 5.

[0066] Figures 4 shows results of comparative studies of analysis of vitamin D₃ in different sample matrices. Details are described in Example 11.

[0067] Figure 5A shows an exemplary Q1 scan spectrum (covering the m/z range of about 300 to 450) for ionization of vitamin D₂. Figure 5B shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D₂ precursor ion with m/z of about 397.2. Figure 5C shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D₂ precursor ion with m/z of about 379.2. Details are described in Example 14.

[0068] Figure 6A shows an exemplary Q1 scan spectrum (covering the m/z range of about 300 to 450) for vitamin D₂-[6, 19, 19]-²H₃ ions. Figure 6B shows an exemplary product ion spectra

(covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_2 -[6, 19, 19]- 2H_3 precursor ion with m/z of about 400.2. Figure 6C shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_2 -[6, 19, 19]- 2H_3 precursor ion with m/z of about 382.2. Details are described in Example 14.

[0069] Figure 7A shows an exemplary Q1 scan spectrum (covering the m/z range of about 300 to 450) for vitamin D_2 -[26, 26, 26, 27, 27, 27]- 2H_6 ions. Figure 7B shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_2 -[26, 26, 26, 27, 27, 27]- 2H_6 precursor ion with m/z of about 403.2. Figure 7C shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_2 -[26, 26, 26, 27, 27, 27]- 2H_6 precursor ion with m/z of about 385.2. Details are described in Example 14.

[0070] Figure 8A shows an exemplary Q1 scan spectrum (covering the m/z range of about 300 to 450) for ionization of vitamin D_3 . Figure 8B shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_3 precursor ion with m/z of about 385.2. Figure 8C shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_3 precursor ion with m/z of about 367.2. Details are described in Example 14.

[0071] Figure 9A shows an exemplary Q1 scan spectrum (covering the m/z range of about 300 to 450) for vitamin D_3 -[6, 19, 19]- 2H_3 ions. Figure 9B shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_3 -[6, 19, 19]- 2H_3 precursor ion with m/z of about 388.2. Figure 9C shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_3 -[6, 19, 19]- 2H_3 precursor ion with m/z of about 370.2. Details are described in Example 14.

[0072] Figure 10A shows an exemplary Q1 scan spectrum (covering the m/z range of about 300 to 450) for vitamin D_3 -[26, 26, 26, 27, 27, 27]- 2H_6 ions. Figure 10B shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_3 -[26, 26, 26, 27, 27, 27]- 2H_6 precursor ion with m/z of about 391.2. Figure 10C shows an exemplary product ion spectra (covering the m/z range of about 100 to 400) for fragmentation of the vitamin D_3 -[26, 26, 26, 27, 27, 27]- 2H_6 precursor ion with m/z of about 373.2. Details are described in Example 14.

[0073] Figure 11A shows an exemplary Q1 scan spectrum (covering the m/z range of about 500 to 620) for PTAD-vitamin D₂ ions. Figure 11B shows an exemplary product ion spectra (covering the m/z range of about 250 to 350) for fragmentation of the PTAD-vitamin D₂ precursor ion with m/z of about 572.2. Details are described in Example 15.

[0074] Figure 12A shows an exemplary Q1 scan spectrum (covering the m/z range of about 500 to 620) for PTAD-vitamin D₂-[6, 19, 19]-²H₃ ions. Figure 12B shows an exemplary product ion spectra (covering the m/z range of about 250 to 350) for fragmentation of the PTAD-vitamin D₂-[6, 19, 19]-²H₃ precursor ion with m/z of about 575.2. Details are described in Example 15.

[0075] Figure 13A shows an exemplary Q1 scan spectrum (covering the m/z range of about 500 to 620) for PTAD-vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆ ions. Figure 13B shows an exemplary product ion spectra (covering the m/z range of about 250 to 350) for fragmentation of the PTAD-vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆ precursor ion with m/z of about 578.2. Details are described in Example 15.

[0076] Figure 14A shows an exemplary Q1 scan spectrum (covering the m/z range of about 500 to 620) for PTAD-vitamin D₃ ions. Figure 14B shows an exemplary product ion spectra (covering the m/z range of about 250 to 350) for fragmentation of the PTAD-vitamin D₃ precursor ion with m/z of about 560.2. Details are described in Example 15.

[0077] Figure 15A shows an exemplary Q1 scan spectrum (covering the m/z range of about 500 to 620) for PTAD-vitamin D₃-[6, 19, 19]-²H₃ ions. Figure 15B shows an exemplary product ion spectra (covering the m/z range of about 250 to 350) for fragmentation of the PTAD-vitamin D₃-[6, 19, 19]-²H₃ precursor ion with m/z of about 563.2. Details are described in Example 15.

[0078] Figure 16A shows an exemplary Q1 scan spectrum (covering the m/z range of about 500 to 620) for PTAD-vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆ ions. Figure 16B shows an exemplary product ion spectra (covering the m/z range of about 250 to 350) for fragmentation of the PTAD-vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆ precursor ion with m/z of about 566.2. Details are described in Example 15.

DETAILED DESCRIPTION OF THE INVENTION

[0079] Methods are described for measuring vitamin D in a sample. More specifically, mass spectrometric methods are described for detecting and quantifying vitamin D in a sample. The

methods may utilize Cookson-type reagents, such as PTAD, to generate derivatized vitamin D. However, in some methods, no derivatizing agent is used, and underivatized vitamin D₂ and/or vitamin D₃ are detected by mass spectrometry.

[0080] The methods may use an extraction chromatography technique, such as turbulent flow liquid chromatography (TFLC), to perform a purification of underivatized or derivatized vitamin D₂ and/or vitamin D₃, combined with methods of mass spectrometry (MS), thereby providing a high-throughput assay system for detecting and quantifying vitamin D₂ and/or vitamin D₃ in a sample. Alternatively, in some methods, no chromatography, including extraction chromatography, is necessary for sample analysis. In these methods, the underivatized or derivatized vitamin D₂ and/or vitamin D₃ is ionized with LDTD. Preferred embodiments are particularly well suited for application in large clinical laboratories for automated vitamin D quantification.

[0081] Suitable test samples for use in methods of the present invention include any test sample that may contain the analyte of interest. In some preferred embodiments, a sample is a biological sample; that is, a sample obtained from any biological source, such as an animal, a cell culture, an organ culture, etc. In certain preferred embodiments, samples are obtained from a mammalian animal, such as a dog, cat, horse, etc. Particularly preferred mammalian animals are primates, most preferably male or female humans. Preferred samples comprise bodily fluids such as blood, plasma, serum, saliva, cerebrospinal fluid, or tissue samples; preferably plasma (including EDTA and heparin plasma) and serum; most preferably serum. Such samples may be obtained, for example, from a patient; that is, a living person, male or female, presenting oneself in a clinical setting for diagnosis, prognosis, or treatment of a disease or condition.

[0082] The present invention also contemplates kits for a vitamin D quantitation assay. A kit for a vitamin D quantitation assay may include a kit comprising the compositions provided herein. For example, a kit may include packaging material and measured amounts of a Cookson-type reagent and an isotopically labeled internal standard, in amounts sufficient for at least one assay. Typically, the kits will also include instructions recorded in a tangible form (e.g., contained on paper or an electronic medium) for using the packaged reagents for use in a vitamin D quantitation assay.

[0083] Calibration and QC pools for use in embodiments of the present invention are preferably prepared using a matrix similar to the intended sample matrix.

Sample Preparation for Mass Spectrometric Analysis

[0084] In preparation for mass spectrometric analysis, vitamin D may be enriched relative to one or more other components in the sample (e.g. protein) by various methods known in the art, including for example, liquid chromatography, filtration, centrifugation, thin layer chromatography (TLC), electrophoresis including capillary electrophoresis, affinity separations including immunoaffinity separations, extraction methods including ethyl acetate or methanol extraction, and the use of chaotropic agents or any combination of the above or the like.

[0085] Protein precipitation is one method of preparing a test sample, especially a biological test sample, such as serum or plasma. Protein purification methods are well known in the art, for example, Polson *et al.*, *Journal of Chromatography B* 2003, 785:263-275, describes protein precipitation techniques suitable for use in methods of the present invention. Protein precipitation may be used to remove most of the protein from the sample leaving vitamin D in the supernatant. The samples may be centrifuged to separate the liquid supernatant from the precipitated proteins; alternatively the samples may be filtered to remove precipitated proteins. The resultant supernatant or filtrate may then be applied directly to mass spectrometry analysis; or alternatively to liquid chromatography and subsequent mass spectrometry analysis. In certain embodiments, samples, such as plasma or serum, may be purified by a hybrid protein precipitation / liquid-liquid extraction method. In these embodiments, a sample is mixed with methanol, ethyl acetate, and water, and the resulting mixture is vortexed and centrifuged. The resulting supernatant is removed, dried to completion and reconstituted in acetonitrile. The purified vitamin D may then be derivatized with any Cookson-type reagent, preferably PTAD or an isotopically labeled variant thereof.

[0086] Another method of sample purification that may be used prior to mass spectrometry is liquid chromatography (LC). Certain methods of liquid chromatography, including HPLC, rely on relatively slow, laminar flow technology. Traditional HPLC analysis relies on column packing in which laminar flow of the sample through the column is the basis for separation of the analyte of interest from the sample. The skilled artisan will understand that separation in such columns is a diffusional process and may select LC, including HPLC, instruments and columns

that are suitable for use with derivatized vitamin D. The chromatographic column typically includes a medium (i.e., a packing material) to facilitate separation of chemical moieties (i.e., fractionation). The medium may include minute particles, or may include a monolithic material with porous channels. A surface of the medium typically includes a bonded surface that interacts with the various chemical moieties to facilitate separation of the chemical moieties. One suitable bonded surface is a hydrophobic bonded surface such as an alkyl bonded, cyano bonded surface, or highly pure silica surface. Alkyl bonded surfaces may include C-4, C-8, C-12, or C-18 bonded alkyl groups. In preferred embodiments, the column is a highly pure silica column (such as a Thermo Hypersil Gold Aq column). The chromatographic column includes an inlet port for receiving a sample and an outlet port for discharging an effluent that includes the fractionated sample. The sample may be supplied to the inlet port directly, or from an extraction column, such as an on-line SPE cartridge or a TFLC extraction column.

[0087] In one embodiment, the sample may be applied to the LC column at the inlet port, eluted with a solvent or solvent mixture, and discharged at the outlet port. Different solvent modes may be selected for eluting the analyte(s) of interest. For example, liquid chromatography may be performed using a gradient mode, an isocratic mode, or a polytypic (i.e. mixed) mode. During chromatography, the separation of materials is effected by variables such as choice of eluent (also known as a “mobile phase”), elution mode, gradient conditions, temperature, etc.

[0088] In certain embodiments, an analyte may be purified by applying a sample to a column under conditions where the analyte of interest is reversibly retained by the column packing material, while one or more other materials are not retained. In these embodiments, a first mobile phase condition can be employed where the analyte of interest is retained by the column, and a second mobile phase condition can subsequently be employed to remove retained material from the column, once the non-retained materials are washed through. Alternatively, an analyte may be purified by applying a sample to a column under mobile phase conditions where the analyte of interest elutes at a differential rate in comparison to one or more other materials. Such procedures may enrich the amount of one or more analytes of interest relative to one or more other components of the sample.

[0089] In one preferred embodiment, HPLC is conducted with an alkyl bonded analytical column chromatographic system. In certain preferred embodiments, a highly pure silica column (such as a Thermo Hypersil Gold Aq column) is used. In certain preferred embodiments, HPLC and/or TFLC are performed using HPLC Grade water as mobile phase A and HPLC Grade ethanol as mobile phase B.

[0090] By careful selection of valves and connector plumbing, two or more chromatography columns may be connected as needed such that material is passed from one to the next without the need for any manual steps. In preferred embodiments, the selection of valves and plumbing is controlled by a computer pre-programmed to perform the necessary steps. Most preferably, the chromatography system is also connected in such an on-line fashion to the detector system, e.g., an MS system. Thus, an operator may place a tray of samples in an autosampler, and the remaining operations are performed under computer control, resulting in purification and analysis of all samples selected.

[0091] In some embodiments, an extraction column may be used for purification of vitamin D metabolites prior to mass spectrometry. In such embodiments, samples may be extracted using an extraction column which captures the analyte, then eluted and chromatographed on a second extraction column or on an analytical HPLC column prior to ionization. For example, sample extraction with a TFLC extraction column may be accomplished with a large particle size (50 μm) packed column. Sample eluted off of this column may then be transferred to an HPLC analytical column for further purification prior to mass spectrometry. Because the steps involved in these chromatography procedures may be linked in an automated fashion, the requirement for operator involvement during the purification of the analyte can be minimized. This feature may result in savings of time and costs, and eliminate the opportunity for operator error.

[0092] In some embodiments, protein precipitation is accomplished with a hybrid protein precipitation / liquid-liquid extraction method which includes methanol protein precipitation and ethyl acetate/water extraction from serum. The resulting vitamin D metabolites may be derivatized prior to being subjected to an extraction column. Preferably, the hybrid protein precipitation / liquid-liquid extraction method and the extraction column are connected in an on-line fashion. In preferred embodiments, the extraction column is a C-8 extraction column, such as a Cohesive Technologies C8XL online extraction column (50 μm particle size, 0.5 x 50 mm)

or equivalent. The eluent from the extraction column may then be applied to an analytical LC column, such as a HPLC column in an on-line fashion, prior to mass spectrometric analysis. Because the steps involved in these chromatography procedures may be linked in an automated fashion, the requirement for operator involvement during the purification of the analyte can be minimized. This feature may result in savings of time and costs, and eliminate the opportunity for operator error.

Detection and Quantitation by Mass Spectrometry

[0093] In various embodiments, derivatized vitamin D may be ionized by any method known to the skilled artisan. Mass spectrometry is performed using a mass spectrometer, which includes an ion source for ionizing the fractionated sample and creating charged molecules for further analysis. For example ionization of the sample may be performed by electron ionization, chemical ionization, electrospray ionization (ESI), photon ionization, atmospheric pressure chemical ionization (APCI), photoionization, atmospheric pressure photoionization (APPI), Laser diode thermal desorption (LDTD), fast atom bombardment (FAB), liquid secondary ionization (LSI), matrix assisted laser desorption ionization (MALDI), field ionization, field desorption, thermospray/plasmaspray ionization, surface enhanced laser desorption ionization (SELDI), inductively coupled plasma (ICP) and particle beam ionization. The skilled artisan will understand that the choice of ionization method may be determined based on the analyte to be measured, type of sample, the type of detector, the choice of positive versus negative mode, etc.

[0094] Derivatized vitamin D may be ionized in positive or negative mode. In preferred embodiments, derivatized vitamin D is ionized by APCI or LDTD in positive ion mode.

[0095] In mass spectrometry techniques generally, after the sample has been ionized, the positively or negatively charged ions thereby created may be analyzed to determine a mass-to-charge ratio. Suitable analyzers for determining mass-to-charge ratios include quadrupole analyzers, ion traps analyzers, and time-of-flight analyzers. Exemplary ion trap methods are described in Bartolucci, *et al.*, *Rapid Commun. Mass Spectrom.* 2000, 14:967-73.

[0096] The ions may be detected using several detection modes. For example, selected ions may be detected, i.e. using a selective ion monitoring mode (SIM), or alternatively, mass

transitions resulting from collision induced dissociation or neutral loss may be monitored, e.g., multiple reaction monitoring (MRM) or selected reaction monitoring (SRM). Preferably, the mass-to-charge ratio is determined using a quadrupole analyzer. For example, in a “quadrupole” or “quadrupole ion trap” instrument, ions in an oscillating radio frequency field experience a force proportional to the DC potential applied between electrodes, the amplitude of the RF signal, and the mass/charge ratio. The voltage and amplitude may be selected so that only ions having a particular mass/charge ratio travel the length of the quadrupole, while all other ions are deflected. Thus, quadrupole instruments may act as both a “mass filter” and as a “mass detector” for the ions injected into the instrument.

[0097] One may enhance the resolution of the MS technique by employing “tandem mass spectrometry,” or “MS/MS”. In this technique, a precursor ion (also called a parent ion) generated from a molecule of interest can be filtered in an MS instrument, and the precursor ion subsequently fragmented to yield one or more fragment ions (also called daughter ions or product ions) that are then analyzed in a second MS procedure. By careful selection of precursor ions, only ions produced by certain analytes are passed to the fragmentation chamber, where collisions with atoms of an inert gas produce the fragment ions. Because both the precursor and fragment ions are produced in a reproducible fashion under a given set of ionization/fragmentation conditions, the MS/MS technique may provide an extremely powerful analytical tool. For example, the combination of filtration/fragmentation may be used to eliminate interfering substances, and may be particularly useful in complex samples, such as biological samples.

[0098] Alternate modes of operating a tandem mass spectrometric instrument include product ion scanning and precursor ion scanning. For a description of these modes of operation, see, e.g., E. Michael Thurman, et al., *Chromatographic-Mass Spectrometric Food Analysis for Trace Determination of Pesticide Residues*, Chapter 8 (Amadeo R. Fernandez-Alba, ed., Elsevier 2005) (387).

[0099] The results of an analyte assay may be related to the amount of the analyte in the original sample by numerous methods known in the art. For example, given that sampling and analysis parameters are carefully controlled, the relative abundance of a given ion may be compared to a table that converts that relative abundance to an absolute amount of the original

molecule. Alternatively, external standards may be run with the samples, and a standard curve constructed based on ions generated from those standards. Using such a standard curve, the relative abundance of a given ion may be converted into an absolute amount of the original molecule. In certain preferred embodiments, an internal standard is used to generate a standard curve for calculating the quantity of vitamin D. Methods of generating and using such standard curves are well known in the art and one of ordinary skill is capable of selecting an appropriate internal standard. For example, in preferred embodiments one or more isotopically labeled vitamin D (e.g., vitamin D₂-[6, 19, 19]-²H₃ and vitamin D₃-[6, 19, 19]-²H₃) may be used as internal standards. Numerous other methods for relating the amount of an ion to the amount of the original molecule will be well known to those of ordinary skill in the art.

[00100] One or more steps of the methods may be performed using automated machines. In certain embodiments, one or more purification steps are performed on-line, and more preferably all of the purification and mass spectrometry steps may be performed in an on-line fashion.

[00101] In certain embodiments, such as MS/MS, where precursor ions are isolated for further fragmentation, collision activated dissociation (CAD) is often used to generate fragment ions for further detection. In CAD, precursor ions gain energy through collisions with an inert gas, and subsequently fragment by a process referred to as “unimolecular decomposition.” Sufficient energy must be deposited in the precursor ion so that certain bonds within the ion can be broken due to increased vibrational energy.

[00102] In some preferred embodiments, vitamin D in a sample is detected and/or quantified using MS/MS as follows. The samples are first purified by protein precipitation or a hybrid protein precipitation / liquid-liquid extraction. Then, vitamin D in the purified sample is optionally derivatized with a Cookson-type reagent, such as PTAD. The purified samples are then subjected to liquid chromatography, preferably on an extraction column (such as a TFLC column) followed by an analytical column (such as a HPLC column); the flow of liquid solvent from a chromatographic column enters the nebulizer interface of an MS/MS analyzer; and the solvent/analyte mixture is converted to vapor in the heated charged tubing of the interface. The analyte(s) (e.g., derivatized or underivatized vitamin D) contained in the solvent, are ionized by applying a large voltage to the solvent/analyte mixture. As the analytes exit the charged tubing of the interface, the solvent/analyte mixture nebulizes and the solvent evaporates, leaving analyte

ions. The analyte ions, e.g. precursor ions, pass through the orifice of the instrument and enter the first quadrupole. Quadrupoles 1 and 3 (Q1 and Q3) are mass filters, allowing selection of ions (i.e., selection of “precursor” and “fragment” ions in Q1 and Q3, respectively) based on their mass to charge ratio (m/z). Quadrupole 2 (Q2) is the collision cell, where ions are fragmented. The first quadrupole of the mass spectrometer (Q1) selects for ions with the mass to charge ratios of interest. Precursor ions with the correct mass/charge ratios are allowed to pass into the collision chamber (Q2), while unwanted ions with any other mass/charge ratio collide with the sides of the quadrupole and are eliminated. Precursor ions entering Q2 collide with neutral argon gas molecules and fragment. The fragment ions generated are passed into quadrupole 3 (Q3), where derivatized or underivatized vitamin D fragment ions are selected while other ions are eliminated.

[00103] The methods may involve MS/MS performed in either positive or negative ion mode; preferably positive ion mode. Using standard methods well known in the art, one of ordinary skill is capable of identifying one or more fragment ions of a particular precursor ion of a derivatized vitamin D that may be used for selection in quadrupole 3 (Q3).

[00104] As ions collide with the detector they produce a pulse of electrons that are converted to a digital signal. The acquired data is relayed to a computer, which plots counts of the ions collected versus time. The resulting mass chromatograms are similar to chromatograms generated in traditional HPLC-MS methods. The areas under the peaks corresponding to particular ions, or the amplitude of such peaks, may be measured and correlated to the amount of the analyte of interest. In certain embodiments, the area under the curves, or amplitude of the peaks, for fragment ion(s) and/or precursor ions are measured to determine the amount of vitamin D. As described above, the relative abundance of a given ion may be converted into an absolute amount of the original analyte using calibration standard curves based on peaks of one or more ions of an internal molecular standard.

EXAMPLES

Example 1: Hybrid Protein Precipitation / Liquid-Liquid Extraction and Cookson-type Derivatization

[00105] The following automated hybrid protein precipitation / liquid-liquid extraction technique was conducted on patient serum samples. Gel Barrier Serum (i.e., serum collected in Serum Separator Tubes) as well as EDTA plasma and Heparin Plasma have also been established as acceptable for this assay.

[00106] A Perkin-Elmer Janus robot and a TomTec Quadra Tower robot was used to automate the following procedure. For each sample, 50 μL of serum was added to a well of a 96 well plate. Then 25 μL of internal standard cocktail (containing isotopically labeled vitamin D_3 -[6, 19, 19]- $^2\text{H}_3$) was added to each well, and the plate vortexed. Then 75 μL of methanol was added, followed by additional vortexing. 300 μL of ethyl acetate and 75 μL of water was then added, followed by additional vortexing, centrifugation, and transfer of the resulting supernatant to a new 96-well plate.

[00107] The transferred liquid in the second 96-well plate was dried to completion under a flowing nitrogen gas manifold. Derivatization was accomplished by adding 100 μL of a 0.1 mg/mL solution of the Cookson-type derivatization agent PTAD in acetonitrile to each well. The derivatization reaction was allowed to proceed for approximately one hour, and was quenched by adding 100 μL of water to the reaction mixture.

Example 2: Extraction of Vitamin D with Liquid Chromatography

[00108] Sample injection was performed with a Cohesive Technologies Aria TX-4 TFLC system using Aria OS V 1.5.1 or newer software.

[00109] The TFLC system automatically injected an aliquot of the above prepared samples into a Cohesive Technologies C8XL online extraction column (50 μm particle size, 005 x 50 mm, from Cohesive Technologies, Inc.) packed with large particles. The samples were loaded at a high flow rate to create turbulence inside the extraction column. This turbulence ensured optimized binding of derivatized vitamin D to the large particles in the column and the passage of excess derivatizing reagent and debris to waste.

[00110] Following loading, the sample was eluted off to the analytical column, a Thermo Hypersil Gold Aq analytical column (5 μm particle size, 50 x 2.1 mm), with a water/ethanol elution gradient. The HPLC gradient was applied to the analytical column, to separate vitamin D from other analytes contained in the sample. Mobile phase A was water and mobile phase B was

ethanol. The HPLC gradient started with a 35 % organic gradient which was ramped to 99 % in approximately 65 seconds.

Example 3: Detection and Quantitation of Derivatized Vitamin D by MS/MS

[00111] MS/MS was performed on the above generated samples using a Finnigan TSQ Quantum Ultra MS/MS system (Thermo Electron Corporation). The following software programs, all from Thermo Electron, were used in the Examples described herein: Quantum Tune Master V 1.5 or newer, Xcalibur V 2.07 or newer, LCQuan V 2.56 (Thermo Finnigan) or newer, and ARIA OS v1.5.1 (Cohesive Technologies) or newer. Liquid solvent/analyte exiting the analytical column flowed to the nebulizer interface of the MS/MS analyzer. The solvent/analyte mixture was converted to vapor in the tubing of the interface. Analytes in the nebulized solvent were ionized by APCI.

[00112] Ions passed to the first quadrupole (Q1), which selected vitamin D₂ and vitamin D₃ precursor ions with a mass-to-charge ratio of 572.3 ± 0.5 m/z and 560.3 ± 0.5 m/z, respectively. Ions entering quadrupole 2 (Q2) collided with argon gas to generate ion fragments, which were passed to quadrupole 3 (Q3) for further selection. Mass spectrometer settings are shown in Table 1. Simultaneously, the same process using isotope dilution mass spectrometry was carried out with internal standard, vitamin D₃-[6, 19, 19]-²H₃. The mass transitions used for detection and quantitation during validation on positive polarity and at the indicated collision energies are shown in Table 2.

Table 1. Mass Spectrometer Settings for Detection of PTAD-vitamin D₂, PTAD-vitamin D₃, and vitamin D₃-[6, 19, 19]-²H₃ (internal standard) (Positive Polarity)

Mass Spectrometric Instrument Settings	
Discharge Current	4.0 μ A
Vaporizer Temperature	300 °C
Sheath Gas Pressure	15
Ion Sweep Gas Pressure	0.0
Aux Gas Pressure	5
Capillary Temperature	300 °C
Skimmer Offset	-10 V

Collision Pressure	1.5 mTorr
Collision Cell Energy	15 V

Table 2. Mass Transitions for PTAD derivatized vitamin D₂, vitamin D₃, and vitamin D₃-[6, 19, 19]-²H₃ (internal standard) (Positive Polarity)

<i>Analyte</i>	<i>Precursor Ion (m/z)</i>	<i>Product Ions (m/z)</i>
PTAD-vitamin D ₂	572.3 ± 0.5	298.1 ± 0.5
PTAD-vitamin D ₃	560.3 ± 0.5	298.1 ± 0.5
PTAD-vitamin D ₃ - [6, 19, 19]- ² H ₃	563.3 ± 0.5	301.1 ± 0.5

[00113] Exemplary chromatograms for PTAD-vitamin D₂, PTAD-vitamin D₃, PTAD-vitamin D₃-[6, 19, 19]-²H₃ (internal standard), are shown in Figures 1A, 1B, and 1C, respectively.

[00114] Exemplary calibration curves for the determination of vitamin D₂ and vitamin D₃ in serum specimens are shown in Figures 2A and 2B, respectively.

Example 4: Linearity of Response for Detection of Derivatized Vitamin D MS/MS

[00115] Linearity was determined by diluting four pools of serum with high endogenous concentration of either vitamin D₂ or vitamin D₃ and analyzing dilutions of 75%, 50%, and 25% in duplicate according to the methods of Examples 1-3. Specimens may be diluted 1:4 with average recovery of 102%, permitting a Clinically Reportable Range (CRR) of 2-240 ng/mL within precision limits of 85%-115%CV. Measured values and percent recoveries from these studies are shown in Table 3.

Table 3. Data Demonstrating Linearity of Response over Dilution Range

Nutritional Vitamin D ₂ (ng/mL (percent recovery %))						
Concentration	Pool 1	Pool 2	Pool 3	Pool 4	Pool 5	Pool 6
100%	25.0 (100)	24.5 (100)	28.1 (100)	24.4 (100)	28.1 (100)	26.7 (100)
75%	17.8	17.5	21.9	17.3	21.2	20.1

	(95.3)	(95.4)	(103.7)	(94.5)	(100.4)	(100.5)
50%	12.8 (102.3)	13.0 (106.2)	14.0 (99.4)	13.8 (112.8)	14.8 (105.7)	11.8 (88.9)
25%	6.7 (107.9)	6.2 (101.4)	7.7 (109.6)	6.7 (109.6)	7.3 (104.1)	7.4 (110.7)
Nutritional Vitamin D ₃ (ng/mL (percent recovery %))						
Sample concentration	Pool 1	Pool 2	Pool 3	Pool 4	Pool 5	Pool 6
100%	31.2 (100)	33.5 (100)	33.4 (100)	29.1 (100)	29.6 (100)	30.6 (100)
75%	22.1 (94.2)	23.6 (93.8)	25.8 (102.8)	22.0 (100.8)	24.8 (111.9)	24.1 (104.8)
50%	15.9 (101.9)	16.1 (96.1)	16.5 (98.5)	15.7 (107.8)	15.9 (107.8)	15.0 (98.0)
25%	7.6 (97.6)	7.9 (94.3)	8.5 (101.8)	7.3 (100.4)	8.0 (108.4)	8.0 (105.0)

Example 5: Analytical Sensitivity: Lower Limit of Quantitation (LLOQ) and Limit of Detection (LOD)

[00116] The lower limit of quantitation (LLOQ) is the point where measurements become quantitatively meaningful. The analyte response at this LLOQ is identifiable, discrete and reproducible with a precision (i.e., coefficient of variation (CV)) of greater than 20% and an accuracy of 80% to 120%. The LLOQ was determined by assaying samples with known analyte concentrations (2 ng/mL, 5 ng/mL, 10 ng/mL, 15 ng/mL, 20 ng/mL, 40 ng/mL, and 60 ng/mL) in quadruplicate five times according to the methods of Examples 1-3, then evaluating the reproducibility. Analysis of the collected data indicates that samples with concentrations of less than 2 ng/mL yielded CVs of less than 20% for both analytes. Thus, the LLOQ of each analyte was determined to be < 2 ng/mL. Data generated for the determination of LLOQ of PTAD-vitamin D₂ and PTAD-vitamin D₃ are shown in Tables 4 and 5, respectively. The graphical representations of CV versus concentration for both analytes are shown in Figure 3.

Table 4. Determination of Lower Limit of Quantitation of Vitamin D₂

Vitamin D ₂ concentration (ng/mL)	Assay 1	Assay 2	Assay 3	Assay 4	Assay 5		
2	2.1	1.9	2.1	1.9	2.0	Average	2.04
	1.7	2.0	2.1	1.9	2.0	Std Dev	0.2

	1.9	2.0	2.2	2.2	2.3	CV (%)	8.5
	1.9	1.9	2.1	2.5	2.2	Accuracy (%)	101.8
5	4.0	6.2	4.3	4.7	4.0	Average	4.7
	NC	4.6	4.6	4.6	4.6	Std Dev	0.6
	4.6	5.6	5.1	4.7	4.3	CV (%)	12.8
	5.7	4.1	4.6	4.6	5.3	Accuracy (%)	94.8
10	9.8	11.2	10.8	11.3	10.3	Average	10.3
	11.3	9.6	8.6	10.0	10.9	Std Dev	0.7
	10.9	10.6	9.9	9.6	10.1	CV (%)	6.9
	10.0	9.7	9.6	10.2	10.9	Accuracy (%)	102.6
15	13.6	14.7	15.6	15.0	14.5	Average	14.6
	13.0	16.2	15.6	12.7	17.4	Std Dev	1.2
	13.6	13.4	13.6	14.1	13.4	CV (%)	8.2
	14.8	14.6	16.0	14.8	14.6	Accuracy (%)	97.1
20	19.4	19.1	20.3	19.9	19.1	Average	20.2
	22.0	18.6	20.7	21.6	19.4	Std Dev	1.0
	19.0	21.1	21.9	20.8	20.1	CV (%)	4.9
	20.0	19.6	19.3	20.7	20.4	Accuracy (%)	100.8
40	42.2	41.5	39.1	42.5	41.8	Average	40.3
	39.2	40.6	39.8	39.0	37.8	Std Dev	1.6
	40.2	40.0	41.3	40.4	37.6	CV (%)	3.9
	41.1	39.6	38.9	39.6	43.6	Accuracy (%)	100.7
60	64.4	57.8	59.1	59.3	57.0	Average	59.8
	57.4	61.4	62.3	60.8	65.2	Std Dev	2.5
	58.9	59.2	56.8	59.0	62.0	CV (%)	4.2
	58.1	61.1	61.8	59.7	55.6	Accuracy (%)	99.7

Table 5. Determination of Lower Limit of Quantitation of Vitamin D₃

Vitamin D ₃ concentration (ng/mL)	Assay 1	Assay 2	Assay 3	Assay 4	Assay 5		
2	2.1	1.6	2.2	2.0	1.8	Average	2.03
	2.1	1.8	2.2	1.8	1.8	Std Dev	0.3
	1.6	2.3	2.0	2.3	1.7	CV (%)	14.3
	2.0	2.5	2.0	2.0	2.7	Accuracy (%)	101.4
5	4.7	5.6	4.6	5.0	4.9	Average	5.0
	7.1	4.9	4.7	5.3	4.4	Std Dev	0.6
	4.6	5.0	5.4	5.5	5.0	CV (%)	12.2
	4.6	4.7	4.3	4.7	5.1	Accuracy (%)	99.7
10	9.4	10.2	10.1	9.9	11.0	Average	10.0
	11.6	9.6	9.3	10.3	10.0	Std Dev	0.7
	10.1	10.2	9.2	9.0	11.3	CV (%)	7.5
	9.5	8.6	10.0	9.5	10.2	Accuracy (%)	99.5

15	14.5	14.9	14.4	15.7	15.8	Average	15.0
	14.2	16.6	15.7	12.8	15.5	Std Dev	0.9
	15.3	13.8	14.2	15.1	14.4	CV (%)	6.3
	15.7	15.1	16.1	16.0	14.2	Accuracy (%)	100.0
20	19.8	18.4	18.2	22.2	17.5	Average	19.5
	19.0	19.1	25.0	19.9	20.9	Std Dev	1.7
	17.8	20.1	20.0	20.1	18.3	CV (%)	8.9
	17.8	19.7	19.7	18.6	18.6	Accuracy (%)	97.7
40	44.3	41.7	40.6	39.1	42.0	Average	41.0
	39.3	40.3	38.5	40.2	38.6	Std Dev	1.8
	42.6	42.6	42.0	41.7	40.1	CV (%)	4.4
	42.2	42.1	37.8	40.3	44.0	Accuracy (%)	102.5
60	65.0	57.0	60.3	60.6	55.7	Average	59.5
	57.6	61.5	61.1	59.3	62.6	Std Dev	2.7
	55.4	55.6	56.9	60.0	61.6	CV (%)	4.5
	58.2	62.4	61.7	59.3	58.1	Accuracy (%)	99.2

[00117] The limit of detection (LOD) is the point where a measured value is larger than the uncertainty associated with it and is defined arbitrarily as four standard deviations (SD) from the zero concentration. Selectivity is the ability of an analytical method to differentiate and quantify the analyte in the presence of other components in the sample. A blank was analyzed in 20 replicates according to the methods of Examples 1-3 and the resulting area ratios were statistically analyzed to determine that the LOD for both vitamin D₂ and vitamin D₃ are 0.4 ng/mL. Data collected to determine LOD for each analyte is shown in Table 6.

Table 6. Determination of Limit of Detection of Vitamin D₂ and Vitamin D₃

Replicate #	NVD ₂ (Response Ratio)	NVD ₃ (Response Ratio)
1	0.012	0.014
2	0.004	0.005
3	0.008	0.011
4	0.005	0.004
5	0.018	0.03
6	0.004	0.001
7	0.007	0.006
8	0.006	0.002
9	0.014	0.006
10	0.015	0.006
11	0.003	0.008

12	0.005	0.004
13	0.009	0.002
14	0.017	0.016
15	0.007	0.005
16	0.006	0.012
17	0.003	0.001
18	0.002	0.001
19	0.011	0.001
20	0.023	0.009
Mean	0.009	0.007
SD	0.0058	0.0070
CV	65.1	97.0
Mean + 4SD	0.032	0.035
LOD	0.4 ng/mL	0.4 ng/mL

Example 7: Specificity of Detection

[00118] Several samples were prepared with vitamin D₂, vitamin D₃, and spiked amounts of potentially interfering species (including vitamin D metabolites and related compounds) and analyzed according to the methods of Examples 1-3. The compounds tested for potential interference are listed in Table 7. None of the tested compounds demonstrated cross-reactivity with detection of vitamin D₂ or vitamin D₃ according to the methods of Examples 1-3.

Table 7. Compounds Tested for Possible Interference with Detection of Vitamin D₂ or Vitamin D₃

Compounds Tested
1,25(OH) ₂ D ₃
1,25(OH) ₂ D ₂
1,25(OH) ₂ D ₃ -[6,19,19']- ² H
1,25(OH) ₂ D ₃ -[26,26,26,27,27,27]- ² H
1,25(OH) ₂ D ₂ -[26,26,26,27,27,27]- ² H
25OHD ₃
25OHD ₂
25OHD ₃ -IS-[6,19,19']- ² H
25OHD ₂ -IS-[6,19,19']- ² H
25OHD ₃ -IS-[26,26,26,27,27,27]- ² H
25OHD ₂ -IS-[26,26,26,27,27,27]- ² H
vitamin D ₃ -[6,19,19']- ² H
vitamin D ₂ -[6,19,19']- ² H
vitamin D ₃ -[26,26,26,27,27,27]- ² H

vitamin D ₂ -[26,26,26,27,27]- ² H
1-OH-D ₃ (Alfacalcidol)
1-OH-D ₂ (Hectoral)
24,25(OH) ₂ D ₃
25,26(OH) ₂ D ₃
3-epi-25OHD ₃
3-epi-1,25(OH) ₂ D ₃
Dihydrotachysterol
1,25(OH) ₂ D ₃ -26,23-lactone
Paracalcitol (Zemplar)
Calcipotriene (Dovonex)
7-Dehydrocholesterol

Example 8: Reproducibility of Quantitation of Vitamin D₂ and Vitamin D₃

[00119] The intra-assay variation is defined as the reproducibility of a sample within an assay and was determined by assaying 20 replicates of a sample from each of three QC pools according to the methods of Examples 1-3. Data collected from these analyses are shown in Tables 8 and 9 for vitamin D₂ and vitamin D₃, respectively. The concentrations of the analytes in the QC pools were determined to be 6.6 ng/mL, 20.6 ng/mL, and 52.6 ng/mL for vitamin D₂, and 4.9 ng/mL, 20.5 ng/mL, and 48.6 ng/mL for Vitamin D₃. Statistics performed on the results yielded reproducibility for the three QC pools at 5.1%, 4.6%, and 3.9% for vitamin D₂, and 6.4%, 4.0%, and 4.5% for vitamin D₃.

Table 8. Intra-Assay Variation Determination for Vitamin D₂

	Low QC Pool (ng/mL)	Medium QC Pool (ng/mL)	High QC Pool (ng/mL)
	Lot 120709-L	Lot 120709-M	Lot 120709-H
1	6.7	20.3	52.5
2	6.6	21.0	54.7
3	7.0	21.1	51.5
4	6.2	19.5	52.1
5	6.4	20.0	52.0
6	6.2	22.4	53.3
7	7.0	20.4	48.9
8	6.7	21.5	54.6
9	6.7	21.5	52.4
10	6.9	20.9	52.5
11	7.0	19.2	50.2

12	6.2	20.8	52.3
13	7.0	19.6	57.5
14	6.6	20.6	54.7
15	6.2	19.8	53.0
16	6.2	20.6	49.2
17	6.2	22.9	52.5
18	6.2	20.9	51.4
19	7.0	19.6	52.2
20	6.5	20.3	55.2
Mean	6.6	20.6	52.6
SD	0.3	1.0	2.0
CV %	5.1	4.6	3.9

Table 9. Intra-Assay Variation Determination for Vitamin D₃

	Low QC Pool (ng/mL)	Medium QC Pool (ng/mL)	High QC Pool (ng/mL)
	Lot 120709-L	Lot 120709-M	Lot 120709-H
1	4.3	20.0	48.3
2	5.1	21.8	47.2
3	4.8	21.7	45.4
4	4.5	21.3	49.9
5	5.0	20.5	49.8
6	4.7	21.0	48.2
7	5.8	20.4	46.3
8	4.9	21.2	50.8
9	4.7	22.0	49.1
10	4.7	20.1	48.9
11	4.8	19.1	47.3
12	4.7	19.9	49.1
13	4.9	19.5	52.9
14	5.4	20.3	52.7
15	4.7	19.6	51.3
16	4.8	19.9	46.7
17	4.8	21.0	47.5
18	4.7	21.1	46.9
19	5.0	19.7	45.3
20	5.1	20.4	48.8
Mean	4.9	20.5	48.6
SD	0.3	0.8	2.2
CV %	6.4	4.0	4.5

[00120] The inter-assay variation is defined as the reproducibility (CV) of a sample between assays. Using the three QC pools covering the reportable range of the assay, evaluated over 5 assays according to the methods of Examples 1-3, the inter-assay variation (CV) for the pools was determined for vitamin D₂ and vitamin D₃. For Vitamin D₂, the CVs were determined to be 6.7%, 5.6%, and 4.0% with mean concentrations of 6.5 ng/mL, 21.1 ng/mL, and 50.5 ng/mL, respectively. For Vitamin D₃, the CVs were determined to be 6.5%, 5.9%, and 4.2% with mean concentrations of 4.7 ng/mL, 20.8 ng/mL, and 46.8 ng/mL, respectively. Data collected from these analyses are shown in Tables 10 and 11 for vitamin D₂ and vitamin D₃, respectively. All pools met with acceptable reproducibility requirements of ≤ 15% CV.

Table 10. Inter-Assay Variation Determination for Vitamin D₂

Assay	Low QC Pool (ng/mL)	Medium QC Pool (ng/mL)	High QC Pool (ng/mL)
	Lot 120709-L	Lot 120709-M	Lot 120709-H
1	6.7	20.3	52.5
	6.6	21.0	54.7
	7.0	21.1	51.5
	6.2	19.5	52.1
	6.4	20.0	52.0
2	7.3	23.3	49.4
	6.8	22.1	51.6
	6.7	20.6	48.6
	6.7	20.6	47.2
	7.0	21.1	50.4
3	6.3	22.3	49.7
	6.1	23.0	52.2
	6.6	24.0	49.5
	5.6	22.1	49.6
	6.1	19.1	51.7
4	6.1	20.5	47.5
	7.0	21.1	50.1
	6.9	22.2	50.2
	7.3	22.1	48.3
	6.6	20.2	46.8
5	5.8	21.1	52.4
	6.9	19.6	49.6
	6.0	20.8	49.1
	6.4	20.6	56.1

	6.3	20.1	50.8
6	6.3	21.8	50.0
	7.0	22.5	49.8
	6.2	21.6	51.2
	5.9	20.1	50.3
	6.9	20.1	50.9
Mean	6.5	21.1	50.5
SD	0.4	1.2	2.0
CV %	6.7	5.6	4.0

Table 11. Inter-Assay Variation Determination for Vitamin D₃

Assay	Low QC	Medium QC	High QC
	Pool (ng/mL)	Pool (ng/mL)	Pool (ng/mL)
	Lot 120709-L	Lot 120709-M	Lot 120709-H
1	4.3	20.0	48.3
	5.1	21.8	47.2
	4.8	21.7	45.4
	4.5	21.3	49.9
	5.0	20.5	49.8
2	4.5	23.0	46.3
	4.2	22.8	47.8
	4.5	19.8	45.1
	4.4	20.1	44.4
	4.7	21.1	48.7
3	4.9	21.8	44.9
	5.1	20.9	46.8
	5.2	23.9	45.5
	4.4	22.1	45.7
	4.6	19.3	47.1
4	4.5	20.1	43.7
	4.8	21.2	44.6
	5.2	21.9	46.1
	4.7	21.3	44.1
	4.7	20.5	42.8
5	4.5	20.5	49.1
	5.1	20.8	46.9
	4.5	21.3	47.3
	4.8	20.6	45.8
	4.2	19.7	46.1

6	4.7	19.9	49.6
	4.7	18.9	47.6
	4.4	20.2	48.9
	4.4	18.6	48.5
	5.2	19.4	49.4
Mean	4.7	20.8	46.8
SD	0.3	1.2	2.0
CV %	6.5	5.9	4.2

Example 9: Method Correlation Studies for Quantitation of Vitamin D₂

[00121] A method correlation study was performed for quantitation of vitamin D₂ according to the methods of Examples 1-3 by comparing 20 split samples analyzed according to the tandem mass spectrometric methods described herein with extensive off-line extraction followed by HPLC with UV detection. Specimens were analyzed in singles for each method. Data was analyzed by Linear and Deming regressions. Correlation analyses are summarized in Table 12.

Table 12. Correlation Analyses for Method Comparison

Vitamin D ₂ (n=20)	
Linear Regression	$y = 1.26x - 0.55$ $R_2 = 0.96$
Deming Regression	$y = 1.29x - 1.31$

Example 10: Interference Studies

[00122] Hemolysis Interference: The effects of hemolysis in the assay described in Examples 1-3 were evaluated by spiking hemoglobin into serum pools containing elevated vitamin D₂ and vitamin D₃. A fresh blood sample was centrifuged to yield packed red blood cells. The cells were reconstituted in deionized water and frozen to achieve cell lysis. This crude hemoglobin solution was then spiked into the pools to generate lightly (100 mg/dL) and moderately (500 mg/dL) hemolyzed samples. Specimens were analyzed in duplicate according the method in Examples 1-3 and results were compared to the control pool result and a percent difference was calculated. The data shows that none of the hemoglobin spikes was >15% different than control, for either analyte). Therefore, light to moderate hemolyzed specimens are acceptable. For raw

data see Table 13 (%Difference = (Spiked –UnSpiked) / UnSpiked x 100%).

Table 13. Hemolysis Interference Studies

		Vitamin D ₂		Vitamin D ₃	
		ng/mL	% Diff.	ng/mL	% Diff.
Pool 1	Control	22.1		31.2	
	Light Hemolysis	24.3	9.7	33.8	8.5
	Moderate Hemolysis	23.8	7.6	34.3	10.0
Pool 2	Control	22.3		34.8	
	Light Hemolysis	21.3	-4.8	32.6	-6.4
	Moderate Hemolysis	20.7	-7.2	33.3	-4.5
Pool 3	Control	27.3		34.6	
	Light Hemolysis	26.7	-2.2	34.7	0.2
	Moderate Hemolysis	24.6	-9.8	33.4	-3.5
Pool 4	Control	23.7		30.4	
	Light Hemolysis	22.8	-3.8	31.5	3.6
	Moderate Hemolysis	21.2	-10.6	29.8	-2.0
Pool 5	Control	27.7		33.9	
	Light Hemolysis	26.7	-3.4	33.7	-0.4
	Moderate Hemolysis	25.6	-7.5	33.2	-1.9
Pool 6	Control	24.2		33.6	
	Light Hemolysis	24.3	0.3	33.5	-0.4
	Moderate Hemolysis	24.1	0.5	34.1	1.3

[00123] Icteria Interference: The effects of icteria in the assay described in Examples 1-3 were evaluated by spiking bilirubin into serum pools containing elevated vitamin D₂ and vitamin D₃. A concentrated solution of bilirubin was spiked into the pools to generate lightly (10 mg/dL) and moderately (50 mg/dL) icteric specimens. Specimens were analyzed in duplicate according to the method in Examples 1-3 and results were compared to the non-icteric pool result and the accuracy was calculated. The data shows that both analytes are unaffected by icteria (all values within acceptable accuracy range of 85-115%). Therefore, icteric specimens are acceptable. For raw data see Table 14 (%Difference = (Spiked –UnSpiked) / UnSpiked x 100%).

Table 14. Icteria Interference Studies

		Vitamin D ₂		Vitamin D ₃	
		ng/mL	% Diff.	ng/mL	% Diff.
Pool 1	Control	22.1		31.2	
	Light Icteria	22.2	0.3	32.6	4.8
	Moderate Icteria	21.5	-2.6	30.7	-1.6
Pool 2	Control	22.3		34.8	
	Light Icteria	22.0	-1.3	33.4	-4.0
	Moderate Icteria	20.5	-8.2	32.4	-6.9
Pool 3	Control	27.3		34.6	
	Light Icteria	28.3	3.8	35.2	1.8
	Moderate Icteria	27.1	-0.7	35.0	1.0
Pool 4	Control	23.7		30.4	
	Light Icteria	22.2	-6.3	31.7	4.0
	Moderate Icteria	21.7	-8.7	30.7	1.0
Pool 5	Control	27.7		33.9	
	Light Icteria	26.7	-3.7	32.9	-2.8
	Moderate Icteria	27.1	-2.0	32.4	-4.4
Pool 6	Control	24.2		33.6	
	Light Icteria	22.8	-5.8	31.8	-5.5
	Moderate Icteria	22.6	-6.4	30.5	-9.4

[00124] Lipemia Interference: The effect of lipemia in the assay described in Examples 1-3 was evaluated by spiking porcine brain extract into serum pools containing elevated vitamin D₂ and vitamin D₃. Powdered lipid sample (Avanti Polar Lipids) was dissolved into each pool to generate lightly (400 mg/dL) and moderately (2000 mg/dL) lipemic specimens. Specimens were analyzed in duplicate according to the method in Examples 1-3 and results were compared to the control pool result and the accuracy was calculated. The data shows that both analytes are unaffected by lipemia (all values within acceptable accuracy range of 85-115%). For raw data see Table 15 ((%Difference = (Spiked – UnSpiked) / UnSpiked x 100%).

Table 15. Lipemia Interference (Porcine Brain Extract) Studies

		Vitamin D ₂		Vitamin D ₃	
		ng/mL	% Diff.	ng/mL	% Diff.
Pool 1	Control	22.1		31.2	
	Light Lipemia	22.0	-0.6	32.7	4.8
	Moderate Lipemia	23.0	3.8	32.7	5.0
Pool 2	Control	22.3		34.8	
	Light Lipemia	23.0	3.2	35.4	1.6
	Moderate Lipemia	24.2	8.3	34.7	-0.2
Pool 3	Control	27.3		34.6	
	Light Lipemia	28.7	5.4	35.7	2.9
	Moderate Lipemia	28.1	3.3	36.5	5.3
Pool 4	Control	23.7		30.4	
	Light Lipemia	25.5	7.4	33.2	9.0
	Moderate Lipemia	23.2	-2.5	32.7	7.3
Pool 5	Control	27.7		33.9	
	Light Lipemia	27.1	-2.2	34.4	1.5
	Moderate Lipemia	26.7	-3.4	33.9	0.1
Pool 6	Control	24.2		33.6	
	Light Lipemia	26.1	8.0	34.5	2.5
	Moderate Lipemia	24.0	-0.7	33.6	-0.2

[00125] The effects of lipemia in the assay described in Examples 1-3 were also evaluated by spiking Intralipid emulsion into serum pools containing elevated vitamin D₂ and vitamin D₃. To the serum pools, Intralipid (20% emulsion) was added to generate lightly (400 mg/dL) and moderately (2000 mg/dL) lipemic specimens. Specimens were analyzed in duplicate according to the method in Examples 1-3 and results were compared to the control pool result and the accuracy was calculated. The data shows that both analytes are unaffected by lipemia (all values within acceptable accuracy range of 85-115%). For raw data see Table 16.

Table 16. Lipemia Interference (Intralipid) Studies

		Vitamin D ₂		Vitamin D ₃	
		ng/mL	% Diff.	ng/mL	% Diff.
Pool 1	Control	22.1		31.2	
	Light Lipemia	22.5	1.6	33.3	6.8
	Moderate Lipemia	19.8	-10.4	28.8	-7.6
Pool 2	Control	22.3		34.8	
	Light Lipemia	22.4	0.4	33.5	-3.6
	Moderate Lipemia	19.2	-13.8	30.5	-12.3
Pool 3	Control	27.3		34.6	
	Light Lipemia	27.4	0.6	34.5	-0.3
	Moderate Lipemia	24.9	-8.6	32.5	-6.1
Pool 4	Control	23.7		30.4	
	Light Lipemia	23.6	-0.6	31.3	3.0
	Moderate Lipemia	19.9	-16.0	29.8	-2.0
Pool 5	Control	27.7		33.9	
	Light Lipemia	27.4	-1.1	32.0	-5.6
	Moderate Lipemia	25.8	-6.9	31.9	-5.9
Pool 6	Control	24.2		33.6	
	Light Lipemia	24.5	1.2	33.8	0.6
	Moderate Lipemia	22.8	-5.9	30.0	-10.7

[00126] Based upon the two lipemia experiments using porcine brain extract and Intralipid, lipemic specimens are acceptable.

Example 11: Specimen Type Studies

[00127] Specimens were collected from 10 sources into four different Vacutainer® containers. The Vacutainers used were Red-Top (silicon-coated serum tubes), SST (Serum Separator Tubes, which result in gel-barrier serum), EDTA tubes, and Sodium Heparin tubes.

[00128] These 40 samples were analyzed for nutritional vitamin D₃ according to the method in Examples 1-3. Comparative results are presented in Figure 4. The data demonstrates that all four sample types are suitable for analysis.

Example 12: Demonstration of Routine Range for Vitamin D₃

[00129] Serum specimens from 140 patients were analyzed according to the method described in Examples 1-3 to quantitate vitamin D₃. The results ranged from < 2 ng/mL to about 63 ng/mL vitamin D₃, with 95% of the results falling within the range of < 2 ng/mL to about 20 ng/mL. The results of these analyses are presented in Table 17.

Table 17. Routine Range Studies for Vitamin D₃ (in descending order)

Patient	Vit D ₃ (ng/mL)	Patient	Vit D ₃ (ng/mL)	Patient	Vit D ₃ (ng/mL)	Patient	Vit D ₃ (ng/mL)
1	62.8	36	4.3	71	<2	106	<2
2	35.8	37	4.2	72	<2	107	<2
3	24.3	38	4.2	73	<2	108	<2
4	20.2	39	4.0	74	<2	109	<2
5	15.6	40	3.9	75	<2	110	<2
6	15.6	41	3.9	76	<2	111	<2
7	15.5	42	3.9	77	<2	112	<2
8	14.6	43	3.7	78	<2	113	<2
9	14.3	44	3.6	79	<2	114	<2
10	12.7	45	3.4	80	<2	115	<2
11	12.7	46	3.3	81	<2	116	<2
12	12.5	47	3.2	82	<2	117	<2
13	10.5	48	3.2	83	<2	118	<2
14	10.3	49	3.2	84	<2	119	<2
15	9.8	50	3.2	85	<2	120	<2
16	9.6	51	3.2	86	<2	121	<2
17	9.6	52	3.1	87	<2	122	<2
18	7.9	53	2.8	88	<2	123	<2
19	7.4	54	2.8	89	<2	124	<2
20	6.6	55	2.7	90	<2	125	<2
21	6.5	56	2.6	91	<2	126	<2
22	6.4	57	2.6	92	<2	127	<2
23	5.9	58	2.5	93	<2	128	<2
24	5.8	59	2.3	94	<2	129	<2
25	5.5	60	2.2	95	<2	130	<2
26	5.4	61	2.2	96	<2	131	<2
27	5.3	62	2.2	97	<2	132	<2
28	5.1	63	2.2	98	<2	133	<2
29	5.1	64	2.1	99	<2	134	<2
30	5.0	65	2.0	100	<2	135	<2
31	4.9	66	2.0	101	<2	136	<2
32	4.8	67	2.0	102	<2	137	<2
33	4.6	68	2.0	103	<2	138	<2
34	4.6	69	<2	104	<2	139	<2
35	4.4	70	<2	105	<2	140	<2

Example 13: Recovery Studies

[00130] Mix recovery studies were performed by analysis of specimens with naturally elevated levels of 25-hydroxyvitamin D₂ or 25-hydroxyvitamin D₃, and therefore also had some endogenous circulating vitamin D₂ or vitamin D₃. Six pairs of specimens were selected for the studies. From each pair of specimens (generically referred to as specimen A and specimen B), five samples were prepared and analyzed in quadruplicate according to the method in Examples 1-3. The samples corresponded to 100% A, 80% A – 20% B, 50% A – 50% B, 20% A – 80% B, and 100% B. The results of the recovery studies are presented in Tables 18 and 19.

Table 18. Mixed Specimen Recovery Studies

<i>Pool</i>	<i>Spike Level</i>	Vitamin D₂			Vitamin D₃		
		Measured (ng/mL)	Expected (ng/mL)	Recovery (%)	Measured (ng/mL)	Expected (ng/mL)	Recovery (%)
1	100% A	49.0			2.1		
	80/20 A/B	37.8	39.3	103.8	3.3	3.5	106.1
	50/50 A/B	26.0	24.7	94.8	5.6	5.5	96.9
	20/80 A/B	10.5	10.1	96.5	7.6	7.5	98.5
	100% B	0.4			8.8		
2	100% A	5.0			2.6		
	80/20 A/B	4.9	4.2	85.5	6.6	6.6	99.4
	50/50 A/B	2.9	2.9	101.7	11.8	12.5	106.2
	20/80 A/B	1.6	1.7	106.2	17.7	18.5	104.3
	100% B	0.9			22.4		
3	100% A	8.8			3.0		
	80/20 A/B	7.6	7.4	96.7	6.6	6.6	101.2
	50/50 A/B	5.4	5.2	95.8	11.3	12.1	106.7
	20/80 A/B	3.1	3.0	99.1	16.3	17.5	107.6
	100% B	1.6			21.2		

*Measured values are averages of analysis of four aliquots.

Table 19. Summary of Results of Mixed Specimen Recovery Studies

Pool	Vitamin D₂	Vitamin D₃
1	103.8%	106.1%
	94.8%	96.9%
	96.5%	98.5%
2	85.5%	99.4%
	101.7%	106.2%
	106.2%	104.3%
3	96.7%	101.2%
	95.8%	106.7%
	99.1%	107.6%
Avg	97.8%	103.0%
Avg	100.4%	

Example 14: Exemplary spectra from MS/MS analysis of Vitamin D₂ and Vitamin D₃

[00131] Exemplary Q1 scan spectra from the tandem mass spectrometric analysis of vitamin D₂, vitamin D₂-[6, 19, 19]-²H₃, and vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆ are shown in Figures 5A, 6A, and 7A, respectively. These analyses were conducted by directly injecting standard solutions containing the analyte of interest into a Finnigan TSQ Quantum Ultra MS/MS system (Thermo Electron Corporation). A liquid chromatography mobile phase was simulated by passing 800 μL/min of 80% acetonitrile, 20% water with 0.1% formic acid through an HPLC column, upstream of introduction of the analyte. The analytes were ionized by APCI as described above. The spectra were collected by scanning Q1 across a m/z range of about 300 to 450.

[00132] Exemplary product ion scans generated from two different precursor ions for each of vitamin D₂, vitamin D₂-[6, 19, 19]-²H₃, and vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆ are presented in Figures 5B-C, 6B-C, and 7B-C, respectively. The precursor ions selected in Q1 and the collision energies used to generate these product ion spectra are indicated in Table 20.

[00133] Exemplary MRM transitions for the quantitation of vitamin D₂ include fragmenting a precursor ion with a m/z of about 397.2 to a product ion with a m/z of about 159.0; and fragmenting a precursor ion with a m/z of about 379.2 to a product ion with a m/z of about 158.9. Exemplary MRM transitions for the quantitation of vitamin D₂-[6, 19, 19]-²H₃ include fragmenting a precursor ion with a m/z of about 400.2 to a product ion with a m/z of about

147.0; and fragmenting a precursor ion with a m/z of about 382.2 to a product ion with a m/z of about 312.2. Exemplary MRM transitions for the quantitation of vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆ include fragmenting a precursor ion with a m/z of about 403.2 to a product ion with a m/z of about 159.1; and fragmenting a precursor ion with a m/z of about 385.2 to a product ion with a m/z of about 159.0. However, as can be seen in the product ion scans in Figures 5B-C, 6B-C, and 7B-C, several other product ions are generated upon fragmentation of the precursor ions. Additional product ions may be selected from those indicated in Figures 5B-C, 6B-C, and 7B-C to replace or augment the exemplary fragment ions. For example, additional product ions generated by fragmentation of the vitamin D₂ precursor ion with m/z of about 397.2 include ions with m/z of about 146.9, 133.1, and 121.0. Exemplary additional product ions generated by fragmentation of the vitamin D₂ precursor ion with m/z of about 379.2 include ions with m/z of about 283.2, 187.3, and 175.2.

Table 20. Precursor Ions and Collision Cell Energies for Fragmentation of vitamin D₂, vitamin D₂-[6, 19, 19]-²H₃, and vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆

<i>Analyte</i>	<i>Precursor Ion (m/z)</i>	<i>Collision Cell Energy (V)</i>
vitamin D ₂	397.2, 379.2	25
vitamin D ₂ -[6, 19, 19]- ² H ₃	400.2, 382.2	25
vitamin D ₂ -[26, 26, 26, 27, 27, 27]- ² H ₆	403.2, 385.2	25

[00134] Exemplary Q1 scan spectra from the tandem mass spectrometric analysis of vitamin D₃, vitamin D₃-[6, 19, 19]-²H₃, and vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆ are shown in Figures 8A, 9A, and 10A, respectively. These analyses were conducted by directly injecting standard solutions containing the analyte of interest into a Finnigan TSQ Quantum Ultra MS/MS system (Thermo Electron Corporation). A liquid chromatography mobile phase was simulated by passing 800 μL/min of 80% acetonitrile, 20% water with 0.1% formic acid through an HPLC column, upstream of introduction of the analyte. The spectra were collected by scanning Q1 across a m/z range of about 300 to 450.

[00135] Exemplary product ion scans generated from two different precursor ions for each of vitamin D₃, vitamin D₃-[6, 19, 19]-²H₃, and vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆ are

presented in Figures 8B-C, 9B-C, and 10B-C, respectively. The precursor ions selected in Q1 and the collision energies used to generate these product ion spectra are indicated in Table 21.

[00136] Exemplary MRM transitions for the quantitation of vitamin D₃ include fragmenting a precursor ion with a m/z of about 385.2 to a product ion with a m/z of about 147.0; and fragmenting a precursor ion with a m/z of about 367.2 to a product ion with a m/z of about 159.0. Exemplary MRM transitions for the quantitation of vitamin D₃-[6, 19, 19]-²H₃ include fragmenting a precursor ion with a m/z of about 388.2 to a product ion with a m/z of about 147.0; and fragmenting a precursor ion with a m/z of about 370.2 to a product ion with a m/z of about 162.0. Exemplary MRM transitions for the quantitation of vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆ include fragmenting a precursor ion with a m/z of about 391.2 to a product ion with a m/z of about 159.1; and fragmenting a precursor ion with a m/z of about 373.2 to a product ion with a m/z of about 159.0. However, as can be seen in the product ion scans in Figures 8B-C, 9B-C, and 10B-C, several other product ions are generated upon fragmentation of the precursor ions. Additional product ions may be selected from those indicated in Figures 8B-C, 9B-C, and 10B-C to replace or augment the exemplary fragment ions. For example, additional product ions generated by fragmentation of the vitamin D₃ precursor ion with m/z of about 385.2 include ions with m/z of about 159.0, 133.1, and 107.1. Exemplary additional product ions generated by fragmentation of the vitamin D₃ precursor ion with m/z of about 367.2 include ions with m/z of about 172.9, 145.0, and 119.1.

Table 21. Precursor Ions and Collision Cell Energies for Fragmentation of vitamin D₃, vitamin D₃-[6, 19, 19]-²H₃, and vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆

<i>Analyte</i>	<i>Precursor Ion (m/z)</i>	<i>Collision Cell Energy (V)</i>
vitamin D ₃	385.2, 367.2	25
vitamin D ₃ -[6, 19, 19]- ² H ₃	388.2, 370.2	25
vitamin D ₃ -[26, 26, 26, 27, 27, 27]- ² H ₆	391.2, 373.2	25

Example 15: Exemplary spectra from MS/MS analysis of PTAD derivatized vitamin D₂ and vitamin D₃

[00137] PTAD derivatives of vitamin D₂, vitamin D₂-[6, 19, 19]-²H₃, vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆, vitamin D₃, vitamin D₃-[6, 19, 19]-²H₃, and vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆ were prepared by treating aliquots of stock solutions of each analyte with PTAD in acetonitrile. The derivatization reactions were allowed to proceed for approximately one hour, and were quenched by adding water to the reaction mixture. The derivatized analytes were then analyzed according to the procedure outlined above in Examples 2-3.

[00138] Exemplary Q1 scan spectra from the analysis of samples containing PTAD-vitamin D₂, PTAD-vitamin D₂-[6, 19, 19]-²H₃, and PTAD-vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆ are shown in Figures 11A, 12A, and 13A, respectively. These analyses were conducted by directly injecting standard solutions containing the analyte of interest into a Finnigan TSQ Quantum Ultra MS/MS system (Thermo Electron Corporation). A liquid chromatography mobile phase was simulated by passing 800 μL/min of 80% acetonitrile, 20% water with 0.1% formic acid through an HPLC column, upstream of introduction of the analyte. The analytes were ionized by APCI as described above. The spectra were collected by scanning Q1 across a m/z range of about 500 to 620.

[00139] Exemplary product ion scans generated from precursor ions for each of PTAD-vitamin D₂, PTAD-vitamin D₂-[6, 19, 19]-²H₃, and PTAD-vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆ are presented in Figures 11B, 12B, and 13B, respectively. The precursor ions selected in Q1 and the collision energies used to generate these product ion spectra are indicated in Table 22.

[00140] An exemplary MRM transition for the quantitation of PTAD-vitamin D₂ includes fragmenting a precursor ion with a m/z of about 572.2 to a product ion with a m/z of about 297.9. An exemplary MRM transition for the quantitation of PTAD-vitamin D₂-[6, 19, 19]-²H₃ includes fragmenting a precursor ion with a m/z of about 575.2 to a product ion with a m/z of about 301.0. An exemplary MRM transition for the quantitation of PTAD-vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆ includes fragmenting a precursor ion with a m/z of about 578.2 to a product ion with a m/z of about 297.9. However, as can be seen in the product ion scans in Figures 11B, 12B, and 13B, several other product ions are generated upon fragmentation of the precursor ions. Additional product ions may be selected from those indicated in Figures 11B, 12B, and 13B to replace or

augment the exemplary fragment ions. For example, additional product ions generated by fragmentation of the PTAD-vitamin D₂ precursor ion with m/z of about 572.2 include ions with m/z of about 280.1.

Table 22. Precursor Ions and Collision Cell Energies for Fragmentation of PTAD-vitamin D₂, PTAD-vitamin D₂-[6, 19, 19]-²H₃, and PTAD-vitamin D₂-[26, 26, 26, 27, 27, 27]-²H₆

<i>Analyte</i>	<i>Precursor Ion (m/z)</i>	<i>Collision Cell Energy (V)</i>
PTAD-vitamin D ₂	572.2	15
PTAD-vitamin D ₂ -[6, 19, 19]- ² H ₃	575.2	15
PTAD-vitamin D ₂ -[26, 26, 26, 27, 27, 27]- ² H ₆	578.2	15

[00141] Exemplary Q1 scan spectra from the analysis of samples containing PTAD-vitamin D₃, PTAD-vitamin D₃-[6, 19, 19]-²H₃, and PTAD-vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆ are shown in Figures 14A, 15A, and 16A, respectively. These analyses were conducted by directly injecting standard solutions containing the analyte of interest into a Finnigan TSQ Quantum Ultra MS/MS system (Thermo Electron Corporation). A liquid chromatography mobile phase was simulated by passing 800 μL/min of 80% acetonitrile, 20% water with 0.1% formic acid through an HPLC column, upstream of introduction of the analyte. The spectra were collected by scanning Q1 across a m/z range of about 500 to 620.

[00142] Exemplary product ion scans generated from precursor ions for each of PTAD-vitamin D₃, PTAD-vitamin D₃-[6, 19, 19]-²H₃, and PTAD-vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆ are presented in Figures 14B, 15B, and 16B, respectively. The precursor ions selected in Q1 and the collision energies used to generate these product ion spectra are indicated in Table 23.

[00143] An exemplary MRM transition for the quantitation of PTAD-vitamin D₃ includes fragmenting a precursor ion with a m/z of about 560.2 to a product ion with a m/z of about 298.0. An exemplary MRM transition for the quantitation of PTAD-vitamin D₃-[6, 19, 19]-²H₃ includes fragmenting a precursor ion with a m/z of about 563.2 to a product ion with a m/z of about 301.0. An exemplary MRM transition for the quantitation of PTAD-vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆ includes fragmenting a precursor ion with a m/z of about 566.2 to a product ion with a m/z of about 298.0. However, as can be seen in the product ion scans in Figures 14B, 15B, and 16B,

several other product ions are generated upon fragmentation of the precursor ions. Additional product ions may be selected from those indicated in Figures 14B, 15B, and 16B to replace or augment the exemplary fragment ions. For example, additional product ions generated by fragmentation of the PTAD-vitamin D₃ precursor ion with m/z of about 560.2 include ions with m/z of about 280.0.

Table 23. Precursor Ions and Collision Cell Energies for Fragmentation of PTAD-vitamin D₃, PTAD-vitamin D₃-[6, 19, 19]-²H₃, and PTAD-vitamin D₃-[26, 26, 26, 27, 27, 27]-²H₆

<i>Analyte</i>	<i>Precursor Ion (m/z)</i>	<i>Collision Cell Energy (V)</i>
PTAD-vitamin D ₃	560.2	15
PTAD-vitamin D ₃ -[6, 19, 19]- ² H ₃	563.2	15
PTAD-vitamin D ₃ -[26, 26, 26, 27, 27, 27]- ² H ₆	566.2	15

[00144] The contents of the articles, patents, and patent applications, and all other documents and electronically available information mentioned or cited herein, are hereby incorporated by reference in their entirety to the same extent as if each individual publication was specifically and individually indicated to be incorporated by reference. Applicants reserve the right to physically incorporate into this application any and all materials and information from any such articles, patents, patent applications, or other physical and electronic documents.

[00145] The methods illustratively described herein may suitably be practiced in the absence of any element or elements, limitation or limitations, not specifically disclosed herein. Thus, for example, the terms “comprising”, “including,” containing”, etc. shall be read expansively and without limitation. Additionally, the terms and expressions employed herein have been used as terms of description and not of limitation, and there is no intention in the use of such terms and expressions of excluding any equivalents of the features shown and described or portions thereof. It is recognized that various modifications are possible within the scope of the invention claimed. Thus, it should be understood that although the present invention has been specifically disclosed by preferred embodiments and optional features, modification and variation of the invention embodied therein herein disclosed may be resorted to by those skilled in the art, and that such modifications and variations are considered to be within the scope of this invention.

[00146] The invention has been described broadly and generically herein. Each of the narrower species and subgeneric groupings falling within the generic disclosure also form part of the methods. This includes the generic description of the methods with a proviso or negative limitation removing any subject matter from the genus, regardless of whether or not the excised material is specifically recited herein.

[00147] Other embodiments are within the following claims. In addition, where features or aspects of the methods are described in terms of Markush groups, those skilled in the art will recognize that the invention is also thereby described in terms of any individual member or subgroup of members of the Markush group.

THAT WHICH IS CLAIMED IS:

1. A method for determining the amount of vitamin D in a sample by tandem mass spectrometry, the method comprising the steps of:
 - (i) subjecting Cookson-type-derivatized vitamin D from the sample to an ionization source to generate one or more precursor ions detectable by mass spectrometry;
 - (ii) fragmenting at least one of said precursor ions to generate one or more fragment ions detectable by mass spectrometry;
 - (iii) determining the amount of one or more precursor and fragment ions by mass spectrometry; and
 - (iv) relating the amount of ions determined in step (iii) to the amount of a vitamin D in the sample,wherein the sample is subjected to a Cookson-type derivatization reagent under conditions sufficient to generate Cookson-type-derivatized vitamin D prior to step (i).
2. The method of claim 1, wherein the Cookson-type-derivatized-vitamin D is subjected to an extraction column and an analytical column prior to ionization.
3. The method of claim 2, wherein the extraction column is a solid-phase extraction (SPE) column.
4. The method of claim 2, wherein the extraction column is a turbulent flow liquid chromatography (TFLC) column.
5. The method of any one of claims 2-4, wherein the analytical column is a high performance liquid chromatography (HPLC) column.
6. The method of any one of claims 2-5, wherein the extraction and analytical columns and the ionization source of step (i) are connected in an on-line fashion.

7. The method of any one of claims 1-6, wherein said ionization source is an atmospheric pressure chemical ionization (APCI) source.
8. The method of any one of claims 1-7, wherein said tandem mass spectrometry is conducted as multiple reaction monitoring, precursor ion scanning, or product ion scanning.
9. The method of any one of claims 1-8, wherein the Cookson-type derivatizing reagent is selected from the group consisting of 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD), 4-methyl-1,2,4-triazoline-3,5-dione (MTAD), 4-[2-(6,7-dimethoxy-4-methyl-3-oxo-3,4-dihydroquinoxalyl)ethyl]-1,2,4-triazoline-3,5-dione (DMEQTAD), 4-(4-nitrophenyl)-1,2,4-triazoline-3,5-dione (NPTAD), and 4-ferrocenylmethyl-1,2,4-triazoline-3,5-dione (FMTAD), and isotopically labeled variants thereof.
10. The method of claim 9, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD) or an isotopically labeled variant thereof.
11. The method of any one of claims 1-10, wherein said vitamin D comprises one or more of vitamin D₂ and vitamin D₃.
12. The method of any one of claims 1-8, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); and said one or more precursor ions comprise one or more ions selected from the group consisting of ions with a mass/charge ratio (m/z) of 572.4 ± 0.5 and 560.4 ± 0.5 .
13. The method of any one of claims 1-8 and 12, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); and said one or more fragment ions comprise an ion with a mass/charge ratio (m/z) of 298.1 ± 0.5 .
14. The method of any one of claims 1-13, wherein said vitamin D comprises vitamin D₂.
15. The method of claim 14, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); said one or more precursor ions comprise a vitamin D₂ precursor ion with a mass/charge ratio (m/z) of 572.4 ± 0.5 ; and said one or more fragment ions comprise a vitamin D₂ fragment ion with a m/z of 298.1 ± 0.5 .

16. The method of any one of claims 1-15, wherein said vitamin D comprises vitamin D₃.
17. The method of claim 16, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); said one or more precursor ions comprise a vitamin D₃ precursor ion with a mass/charge ratio (m/z) of 560.4 ± 0.5 ; and said one or more fragment ions comprise a vitamin D₃ fragment ion with a m/z of 298.1 ± 0.5 .
18. The method of any one of claims 1-17, wherein two or more vitamin D derivatives are ionized simultaneously.
19. The method of claim 18, wherein said vitamin D comprises vitamin D₂ and vitamin D₃.
20. The method of claim 19, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); said one or more precursor ions comprise a vitamin D₂ precursor ion with a mass/charge ratio (m/z) of 572.4 ± 0.5 and a vitamin D₃ precursor ion with a m/z of 560.4 ± 0.5 ; and said one or more fragment ions comprise a vitamin D₂ fragment ion with a m/z of 298.1 ± 0.5 and a vitamin D₃ fragment ion with a m/z of 298.1 ± 0.5 .
21. The method of any one of claims 1-20, wherein said sample comprises a biological sample.
22. The method of claim 21, wherein said biological sample is from a human, and the amount of vitamin D determined in the sample is the amount present in the sample when taken from the human.
23. The method of any one of claims 21-22, wherein said sample comprises plasma or serum.
24. A method for determining the amount of vitamin D in a sample by mass spectrometry, the method comprising the steps of:
 - (i) subjecting the sample to turbulent flow liquid chromatography (TFLC);
 - (ii) subjecting Cookson-type-derivatized vitamin D from the sample to an ionization source under conditions suitable to generate one or more ions detectable by mass spectrometry;

(iii) determining the amount of one or more of the Cookson-type-derivatized vitamin D ions by mass spectrometry; and

(iv) relating the amount of Cookson-type-derivatized vitamin D ions determined in step (iii) to the amount of a vitamin D in the sample,

wherein the sample is subjected to a Cookson-type derivatizing reagent under conditions sufficient to generate Cookson-type-derivatized vitamin D in the sample prior to ionization.

25. The method of claim 24, wherein the sample is subjected to high performance liquid chromatography (HPLC) after step (i) but prior to step (ii).

26. The method of claim 25, wherein the TFLC of step (i), the HPLC, and the ionization of step (ii) are conducted in an on-line fashion.

27. The method of any one of claims 24-26, wherein the ionization source is an atmospheric pressure chemical ionization (APCI) source.

28. The method of any one of claims 24-27, wherein said mass spectrometry is tandem mass spectrometry.

29. The method of claim 28, wherein said tandem mass spectrometry is conducted as multiple reaction monitoring, precursor ion scanning, or product ion scanning.

30. The method of any one of claims 24-29, wherein the Cookson-type derivatizing reagent is selected from the group consisting of 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD), 4-methyl-1,2,4-triazoline-3,5-dione (MTAD), 4-[2-(6,7-dimethoxy-4-methyl-3-oxo-3,4-dihydroquinoxalyl)ethyl]-1,2,4-triazoline-3,5-dione (DMEQTAD), 4-(4-nitrophenyl)-1,2,4-triazoline-3,5-dione (NPTAD), and 4-ferrocenylmethyl-1,2,4-triazoline-3,5-dione (FMTAD), and isotopically labeled variants thereof.

31. The method of claim 30, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD) or an isotopically labeled variant thereof.

32. The method of any one of claims 24-31, wherein said vitamin D comprises one or more of vitamin D₂ and vitamin D₃.
33. The method of any one of claims 24-32, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); and said one or more precursor ions comprise one or more ions selected from the group consisting of ions with a mass/charge ratio (m/z) of 572.4 ± 0.5 and 560.4 ± 0.5 .
34. The method of any one of claims 24-32, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); and said one or more fragment ions comprise an ion with a mass/charge ratio (m/z) of 298.1 ± 0.5 .
35. The method of any one of claims 24-34, wherein said vitamin D comprises vitamin D₂.
36. The method of claim 35, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); said one or more precursor ions comprise a vitamin D₂ precursor ion with a mass/charge ratio (m/z) of 572.4 ± 0.5 ; and said one or more fragment ions comprise a vitamin D₂ fragment ion with a m/z of 298.1 ± 0.5 .
37. The method of any one of claims 24-36, wherein said vitamin D comprises vitamin D₃.
38. The method of claim 37, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); said one or more precursor ions comprise a vitamin D₃ precursor ion with a mass/charge ratio (m/z) of 560.4 ± 0.5 ; and said one or more fragment ions comprise a vitamin D₃ fragment ion with a m/z of 298.1 ± 0.5 .
39. The method of any one of claims 24-38, wherein two or more vitamin D derivatives are ionized simultaneously.
40. The method of claim 39, wherein said vitamin D comprises vitamin D₂ and vitamin D₃.
41. The method of any one of claims 24-40, wherein the Cookson-type derivatizing reagent is 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD); said one or more precursor ions comprise a vitamin D₂ precursor ion with a mass/charge ratio (m/z) of 572.4 ± 0.5 and a vitamin D₃ precursor ion

with a m/z of 560.4 ± 0.5 ; and said one or more fragment ions comprise a vitamin D₂ fragment ion with a m/z of 298.1 ± 0.5 and a vitamin D₃ fragment ion with a m/z of 298.1 ± 0.5 .


42. The method of any one of claims 24-41, wherein said sample comprises a biological sample.

43. The method of claim 42, wherein said biological sample is from a human, and the amount of vitamin D determined in the sample is the amount present in the sample when taken from the human.

44. The method of any one of claims 42-43, wherein said sample comprises plasma or serum.


1/35

Figure 1

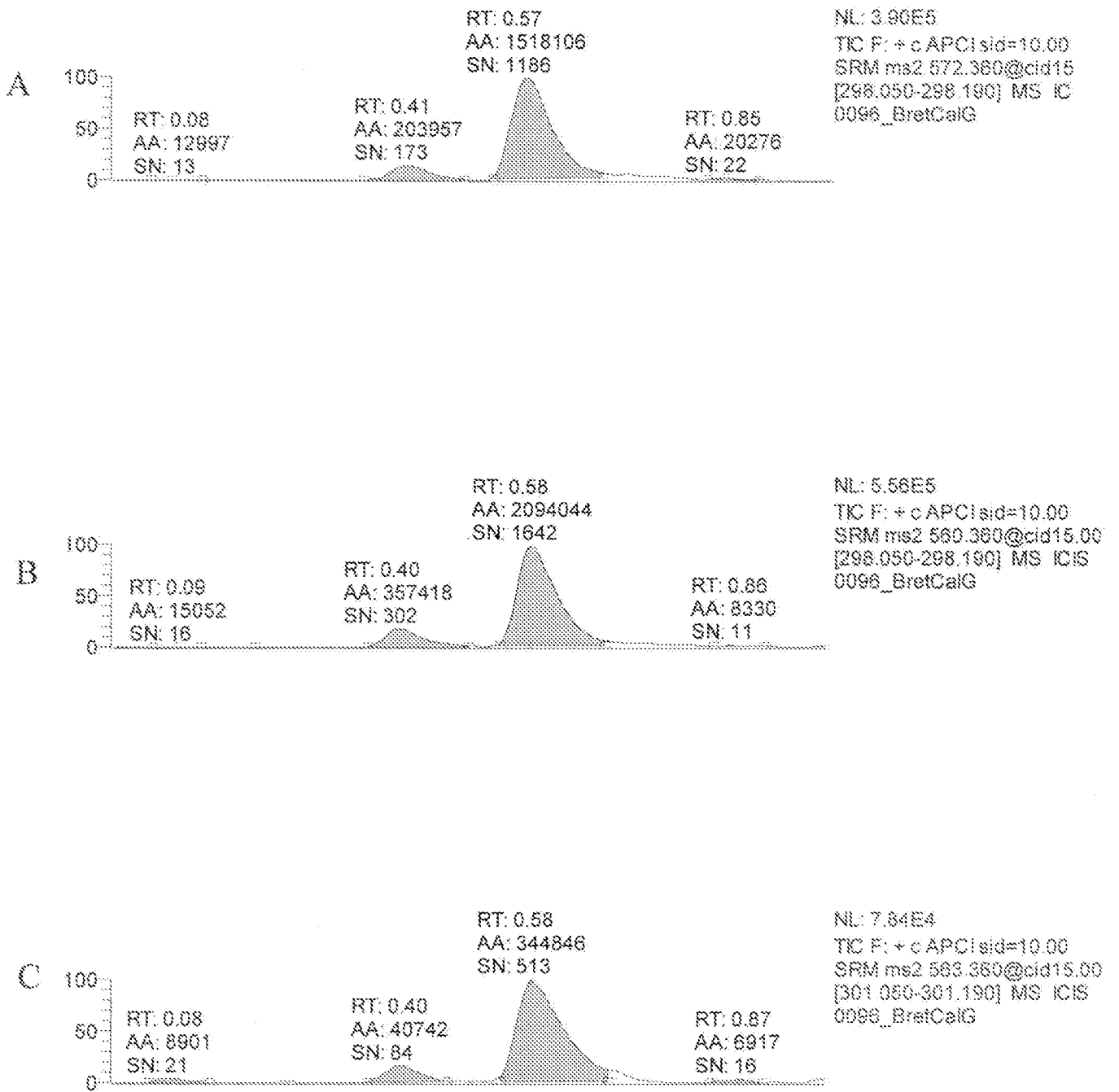


2/35

Figure 2A

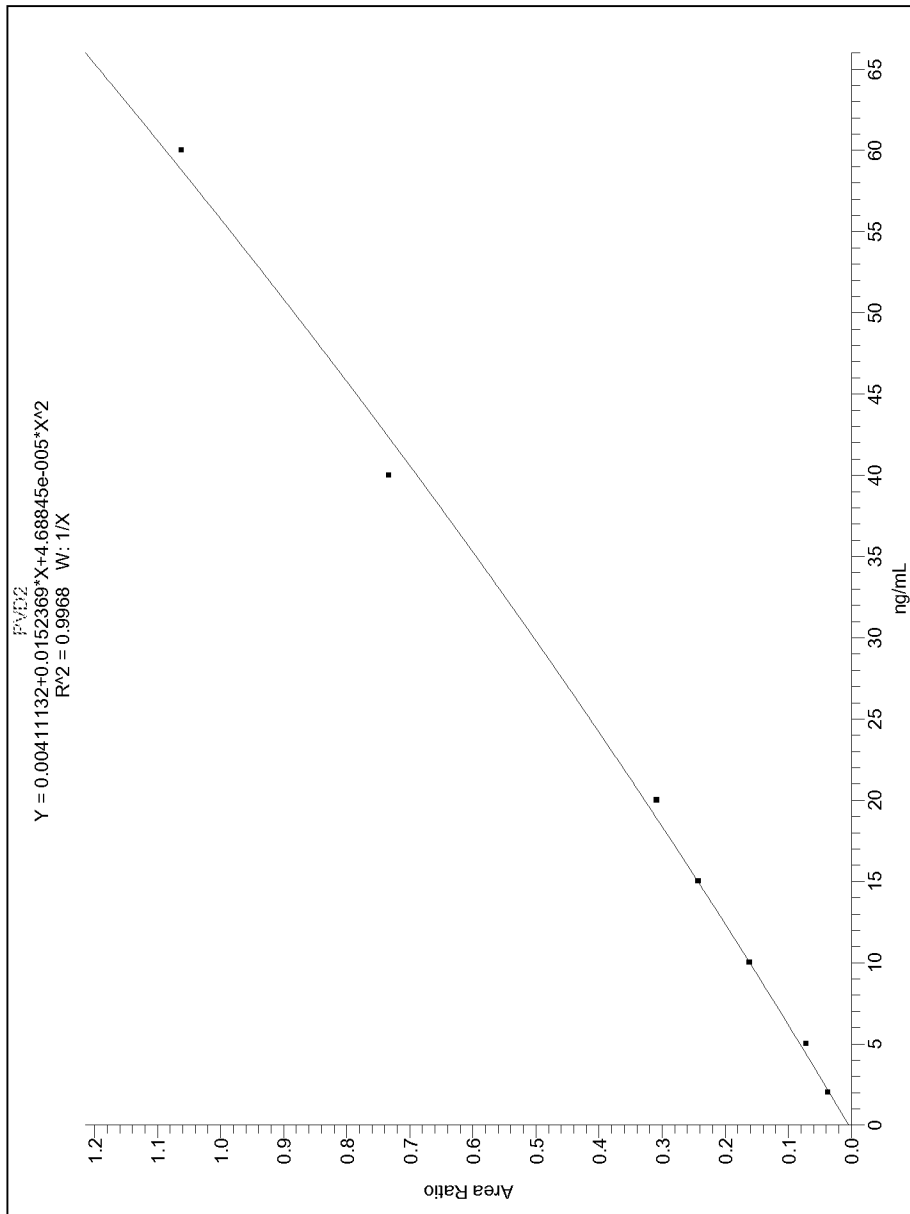

3/35

Figure 2B

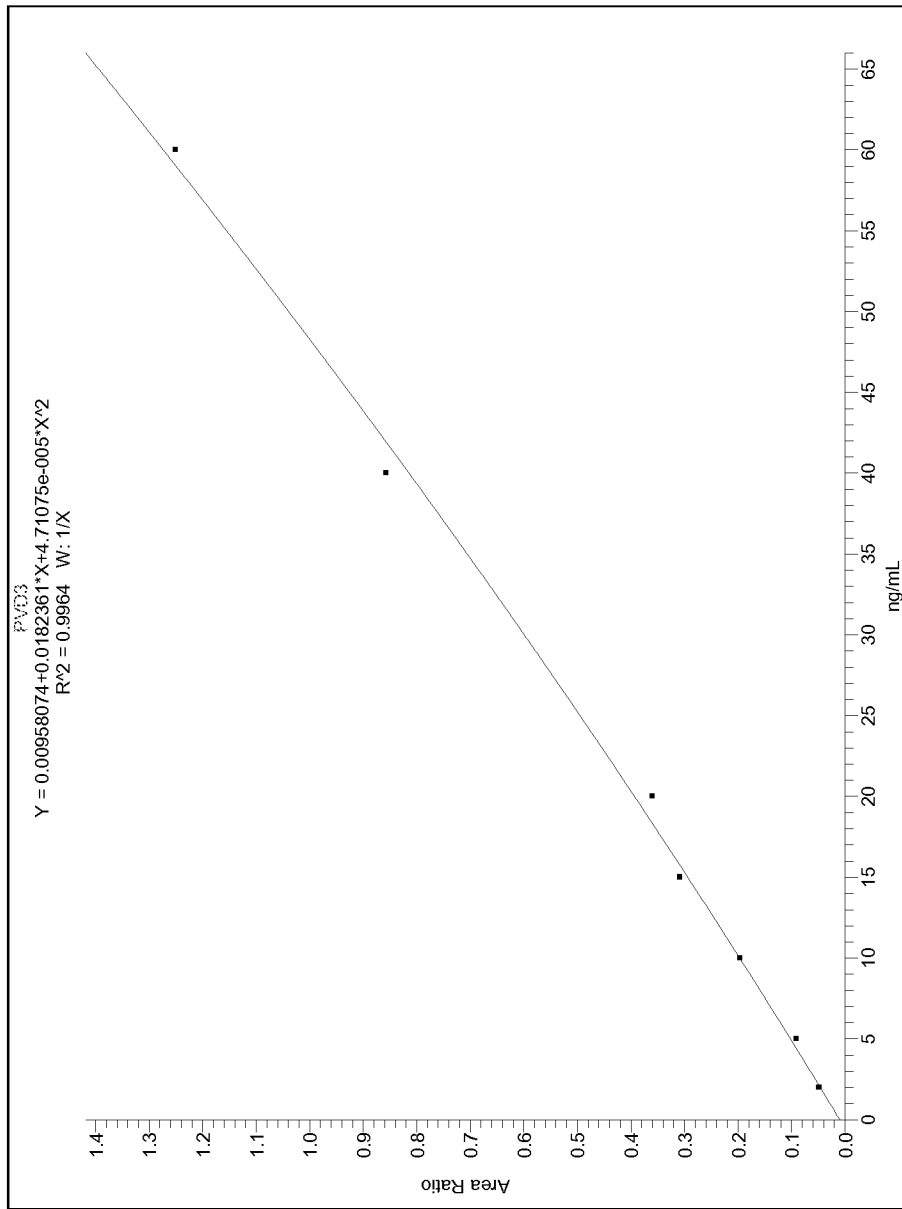

4/35

Figure 3

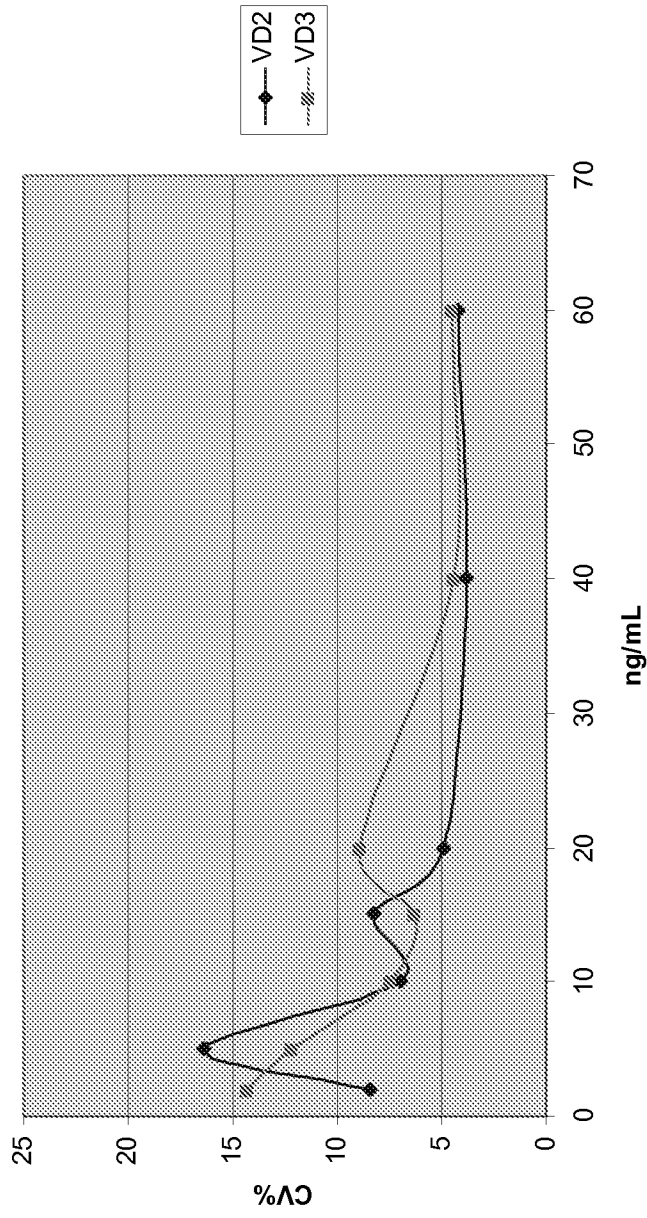


Figure 4



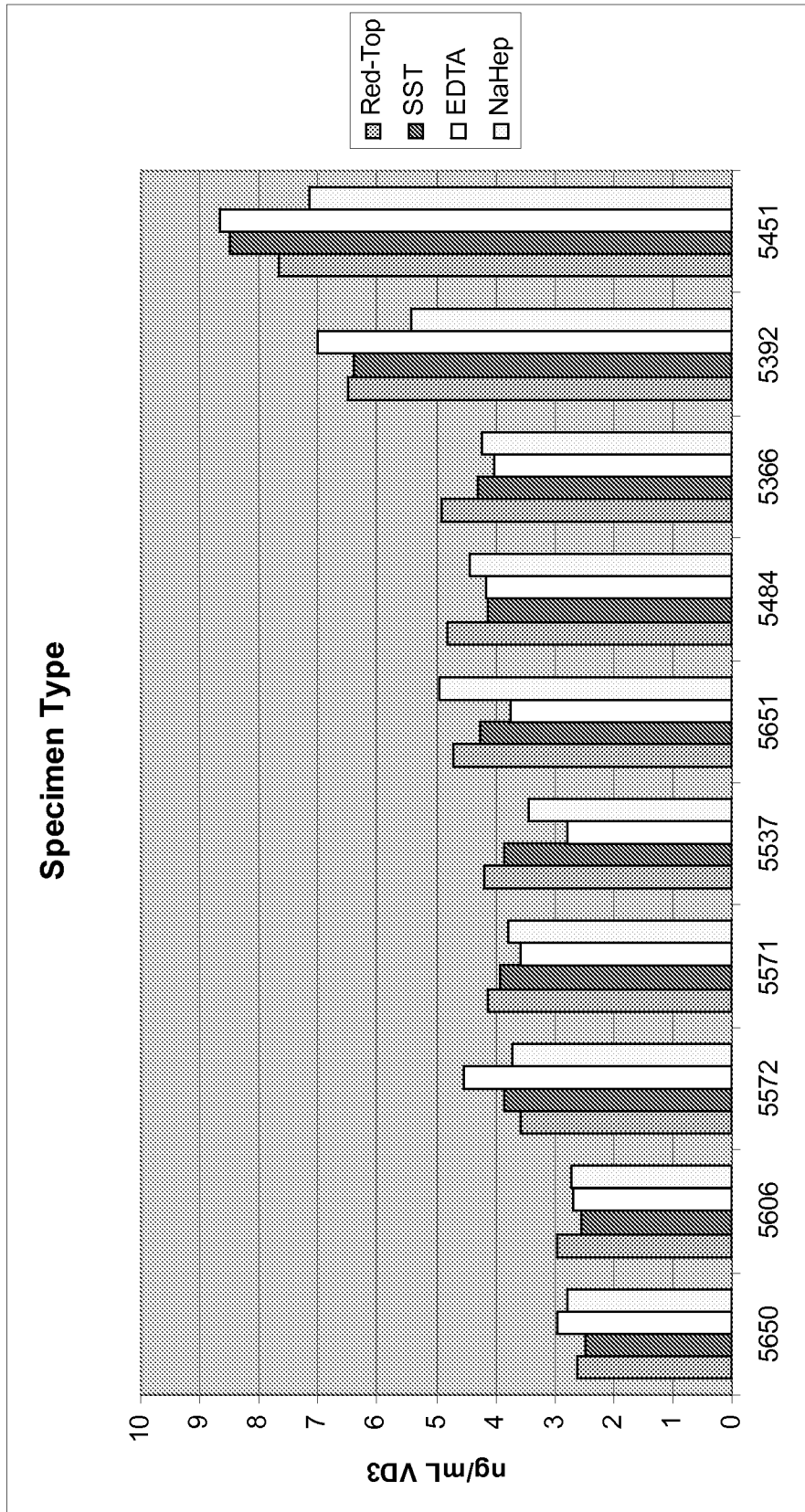


Figure 5A

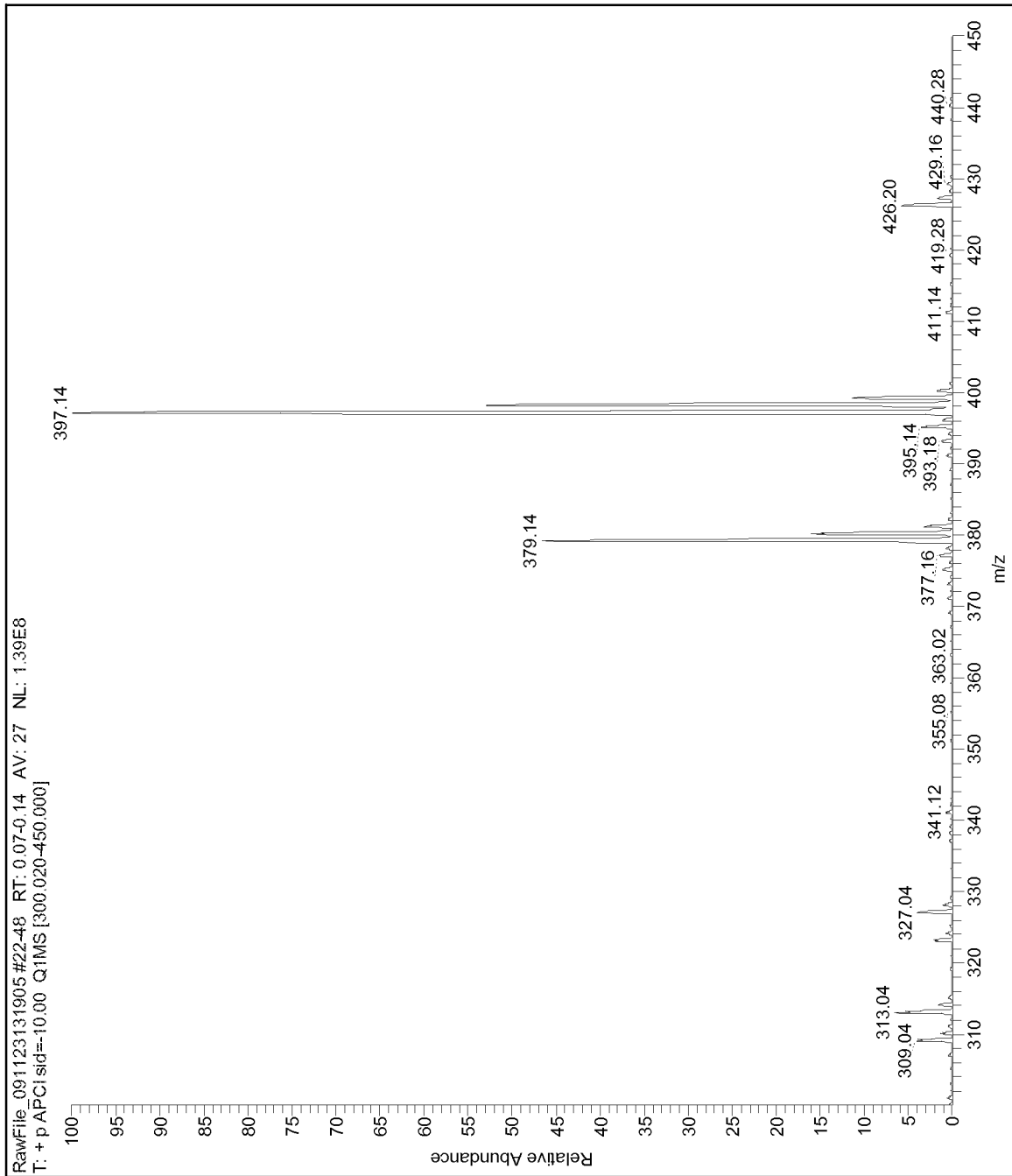


7/35

Figure 5B

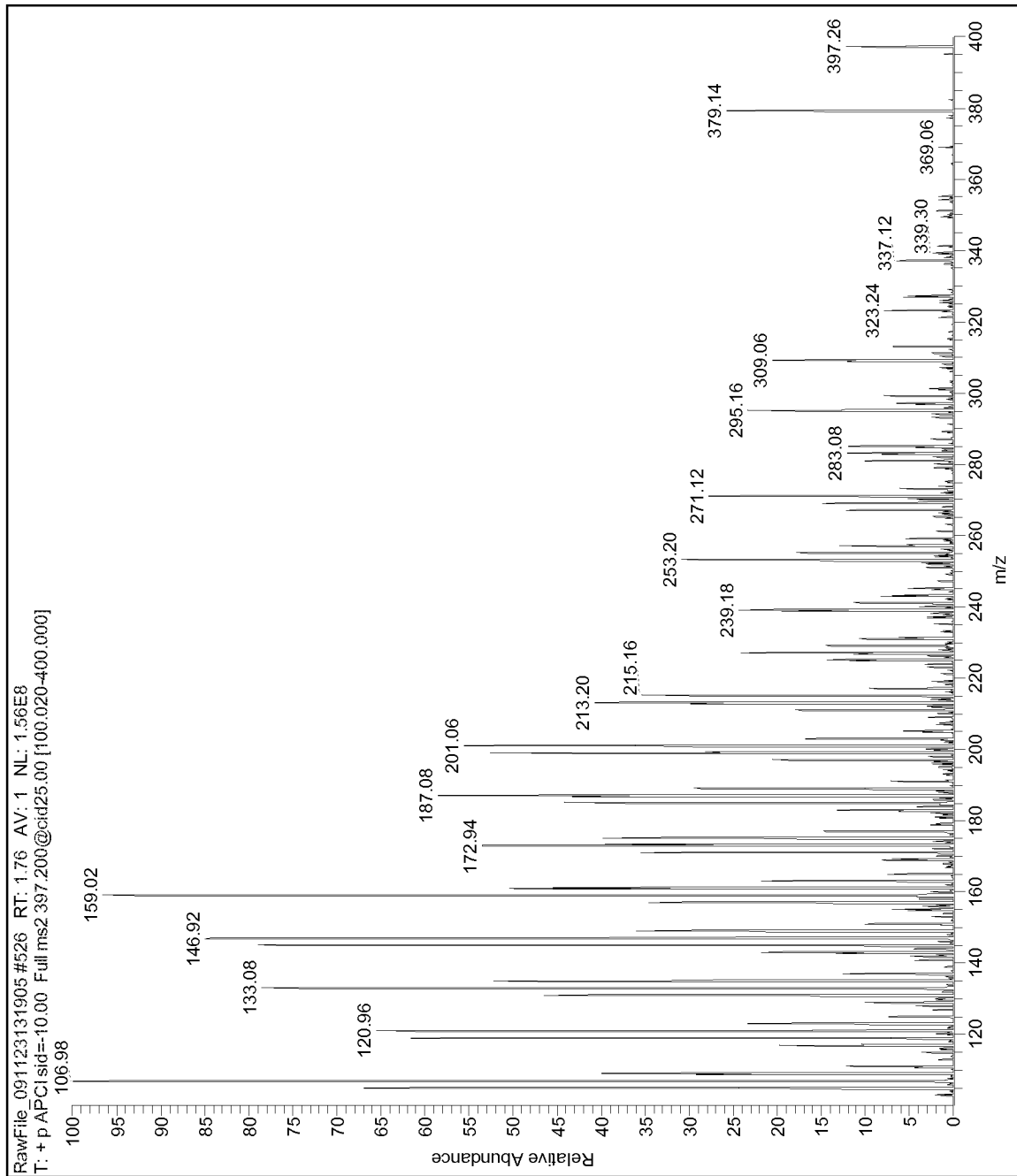


8/35

Figure 5C

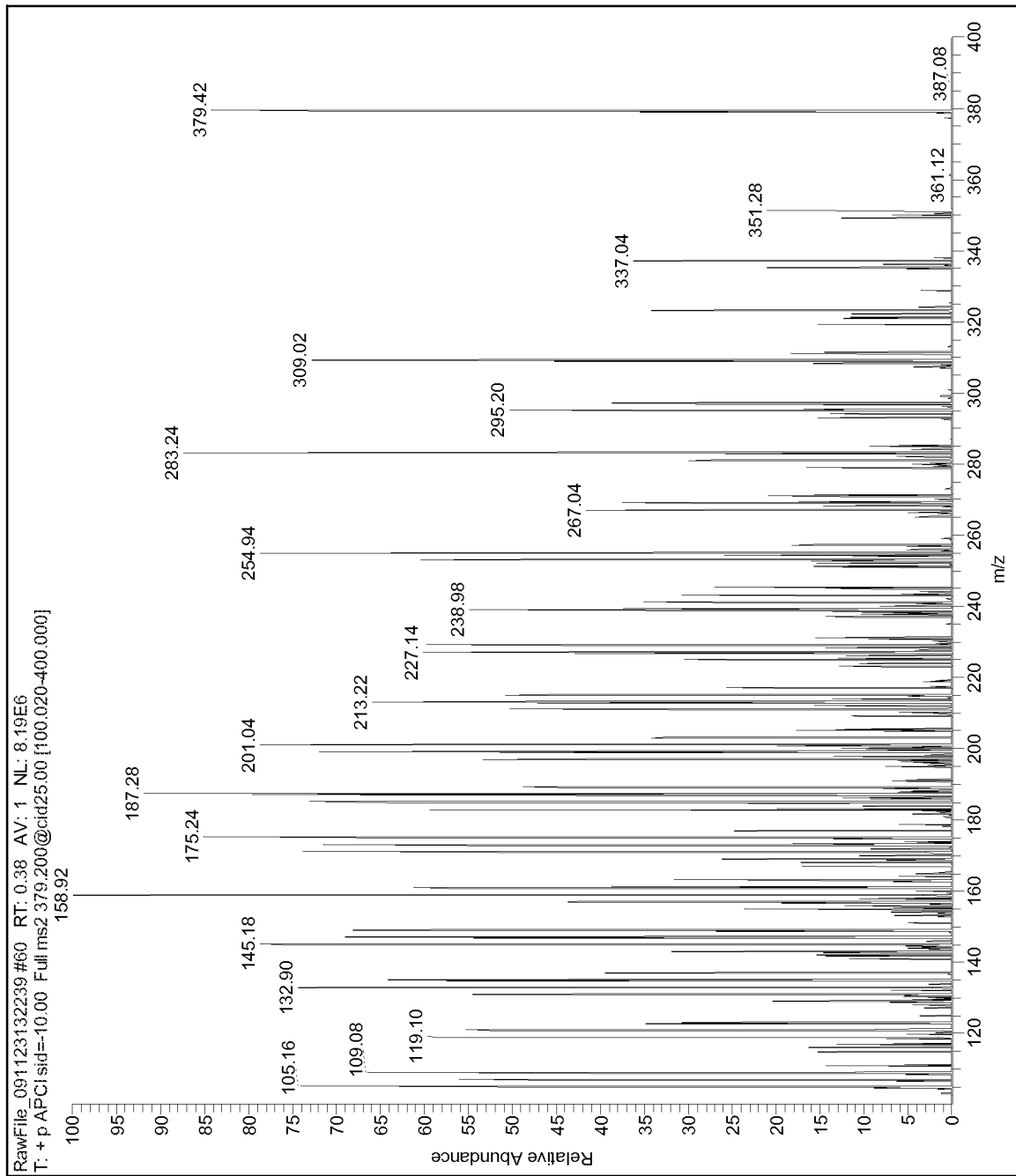


9/35

Figure 6A

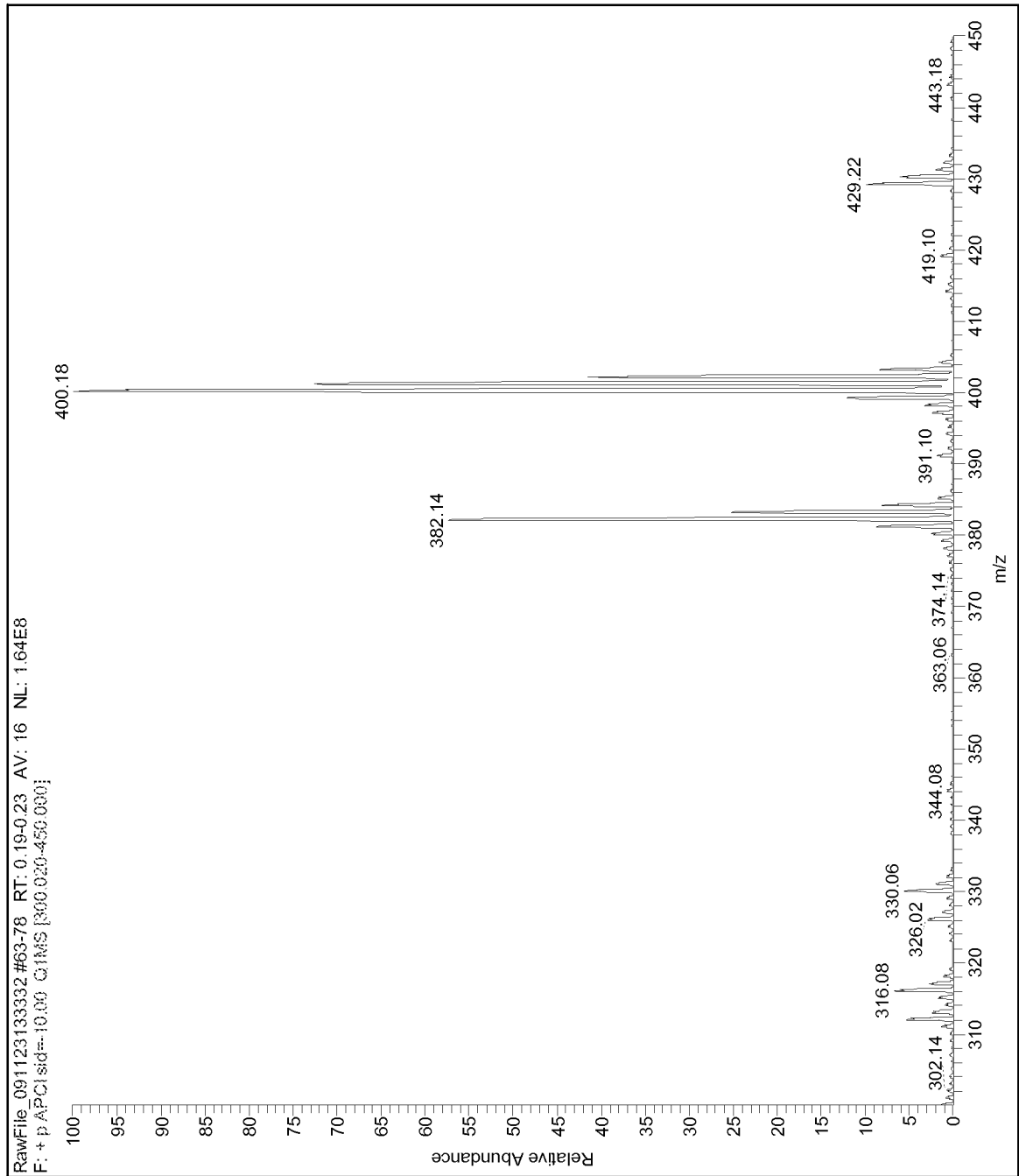


10/35

Figure 6B

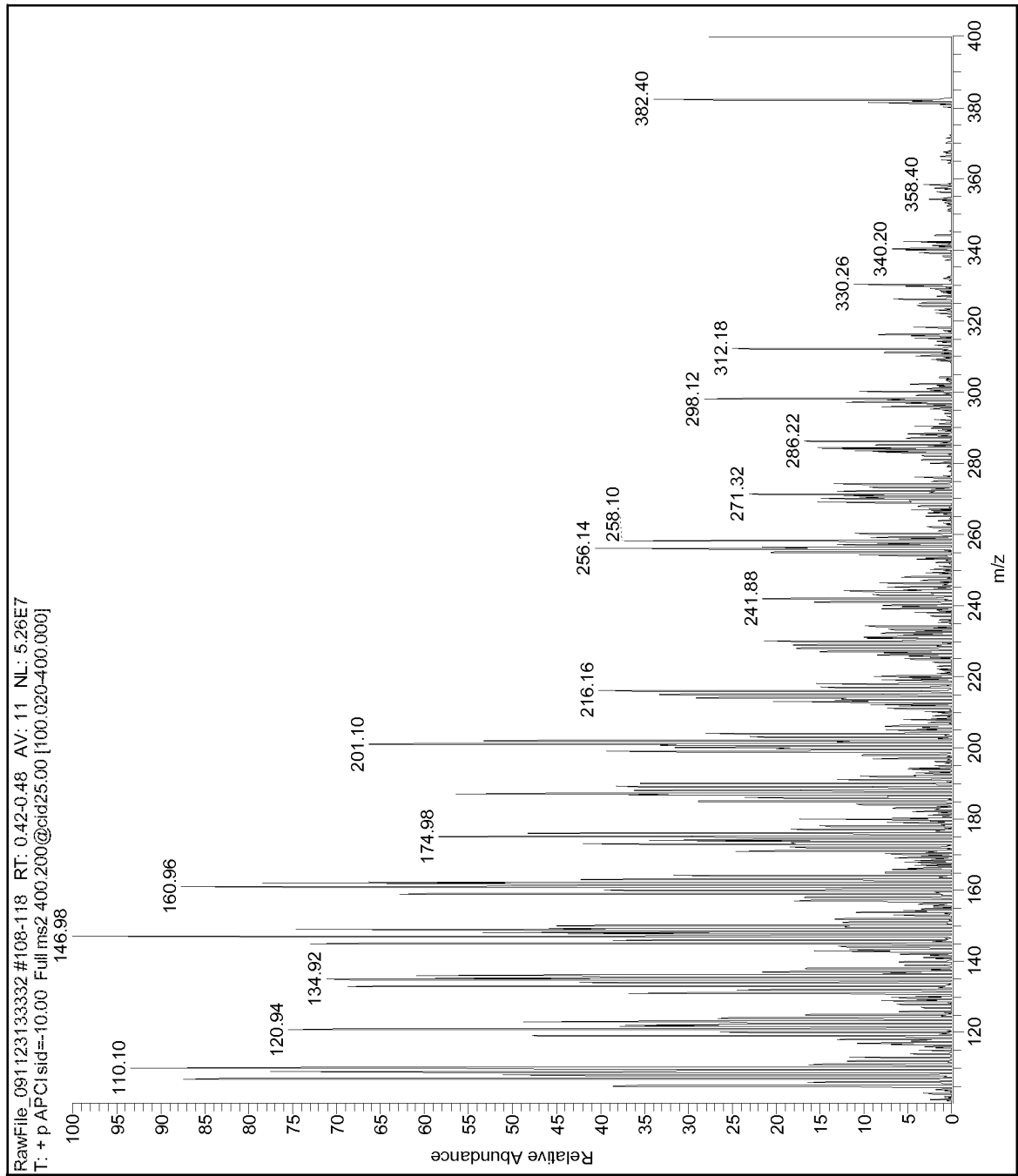


11/35

Figure 6C

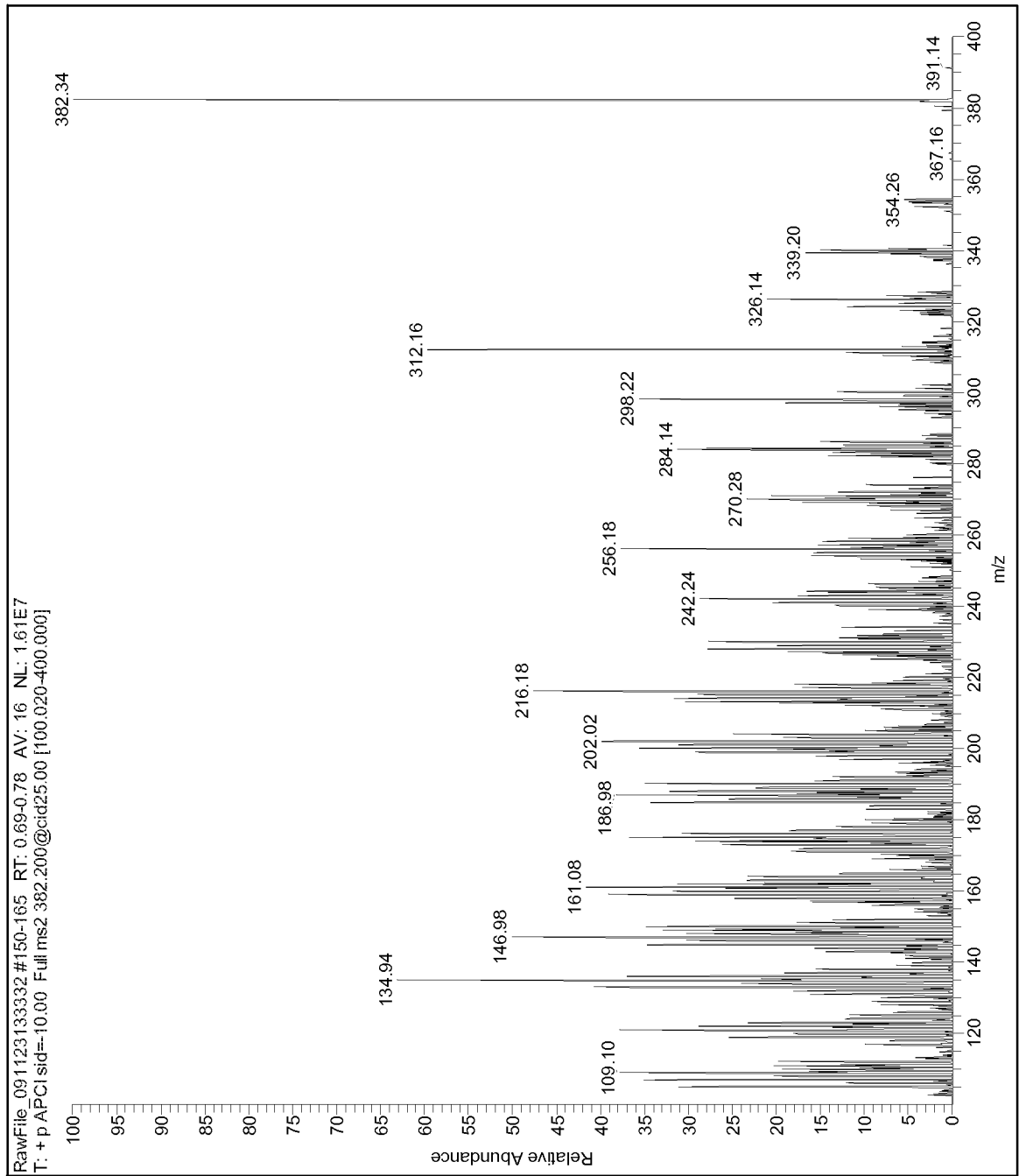


12/35

Figure 7A

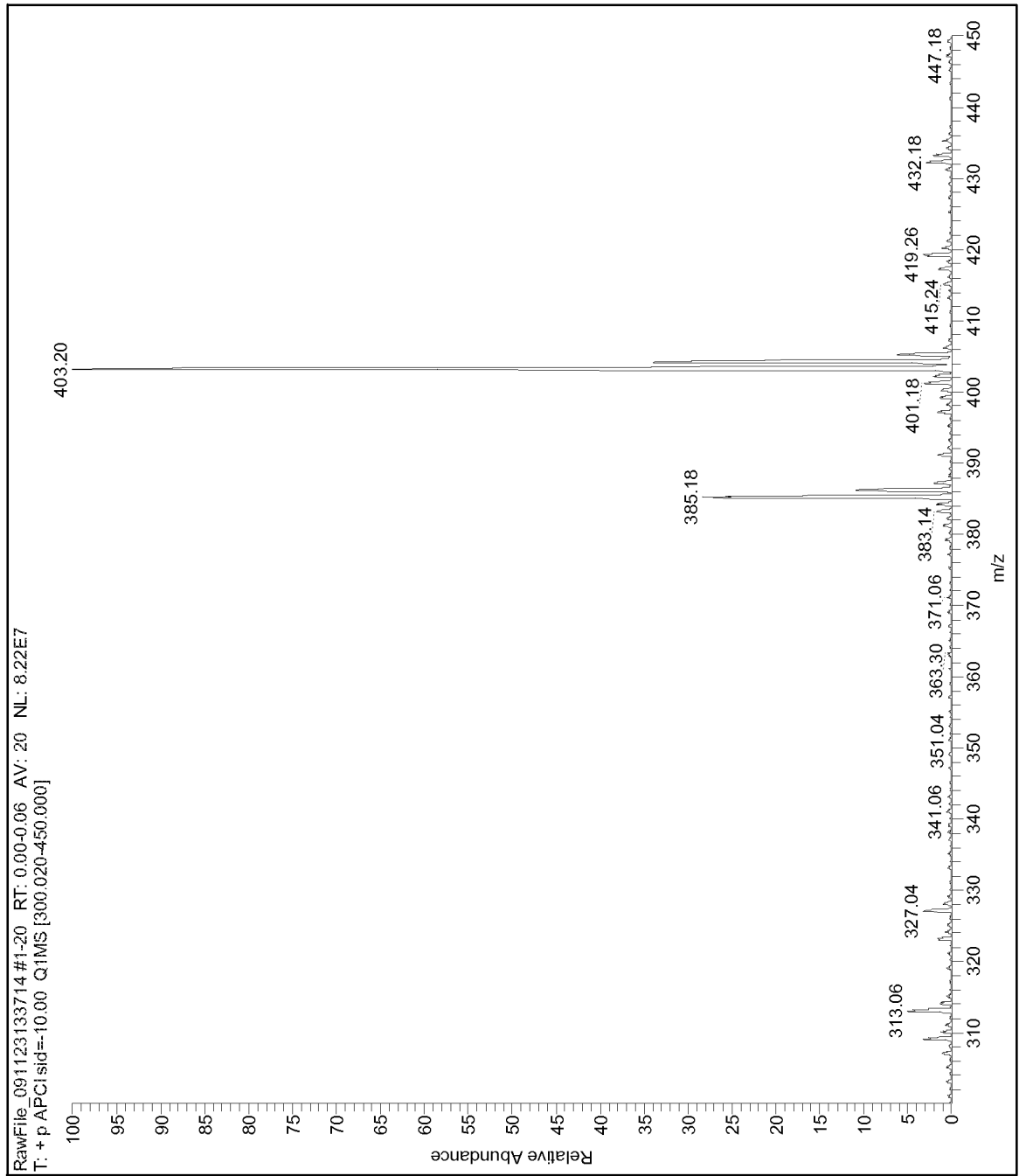


13/35

Figure 7B

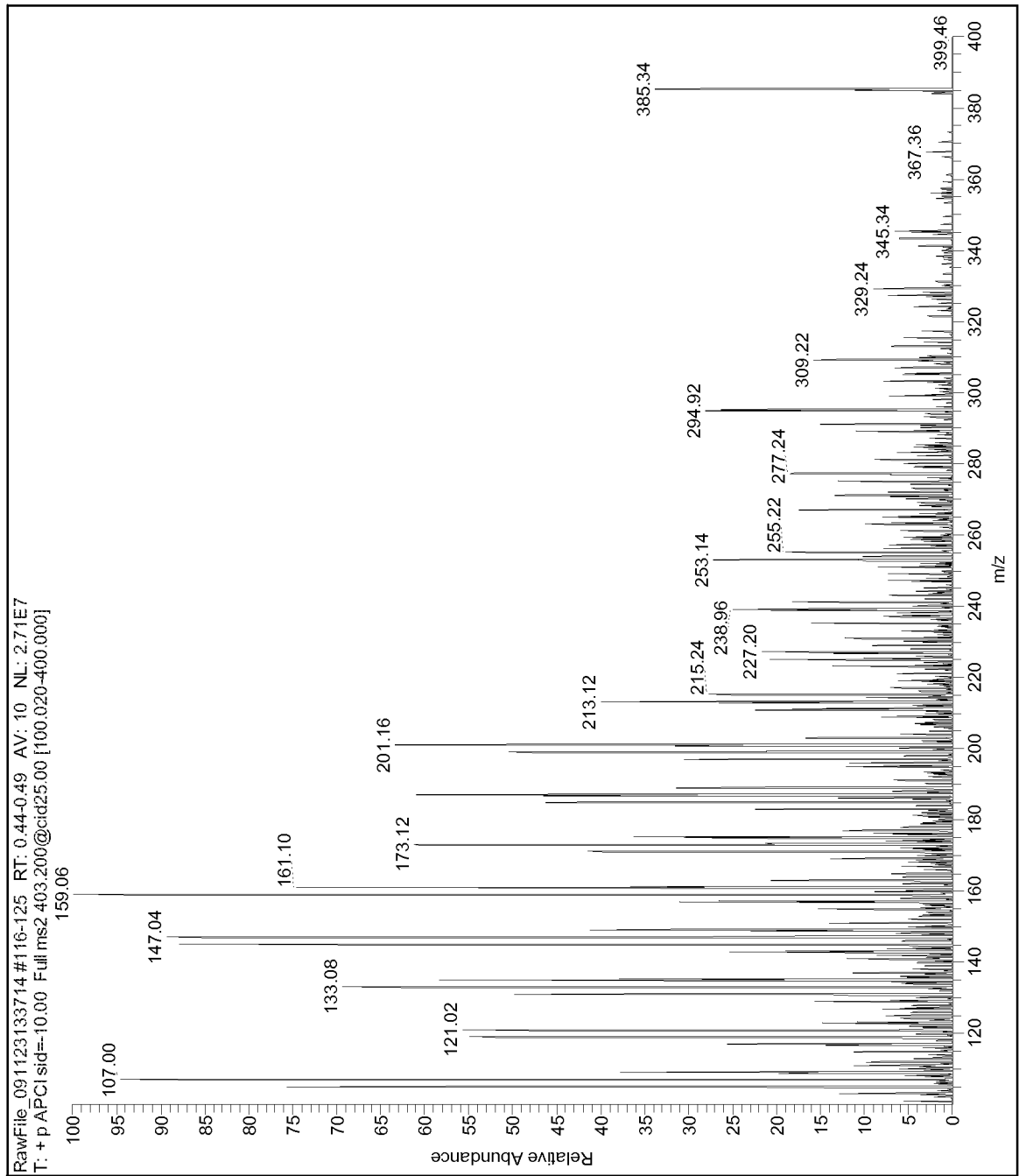


14/35

Figure 7C

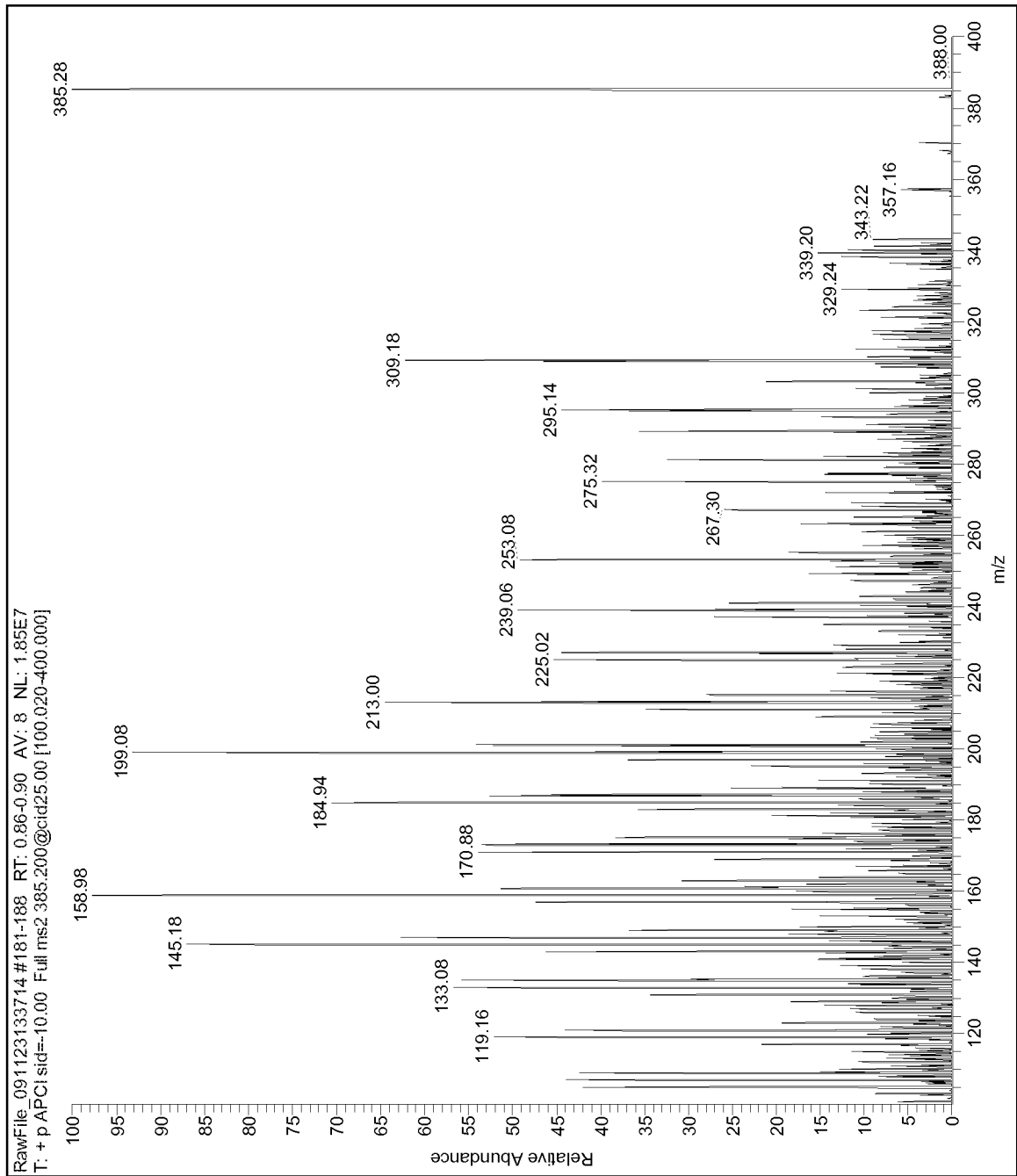


15/35

Figure 8A

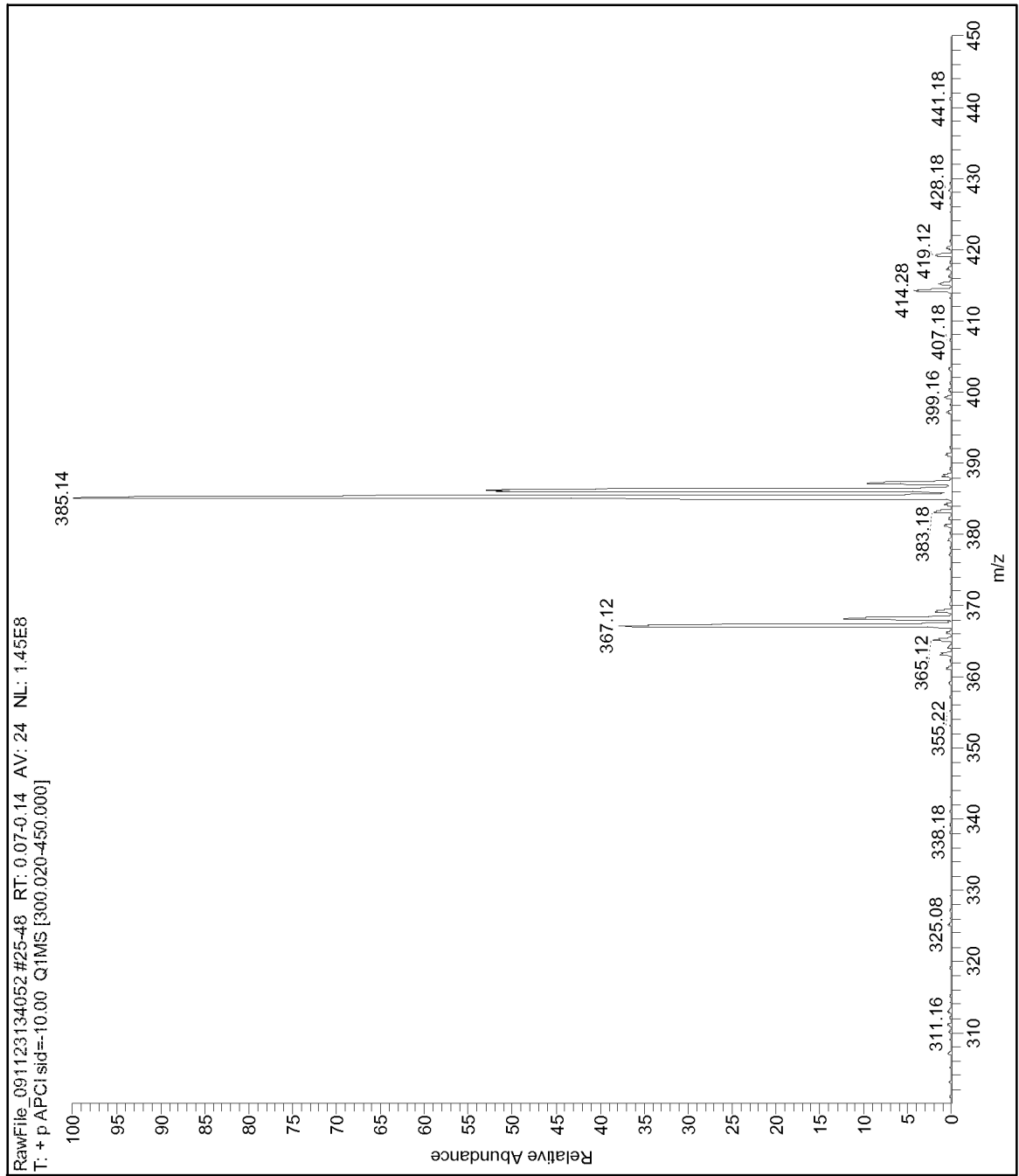


16/35

Figure 8B

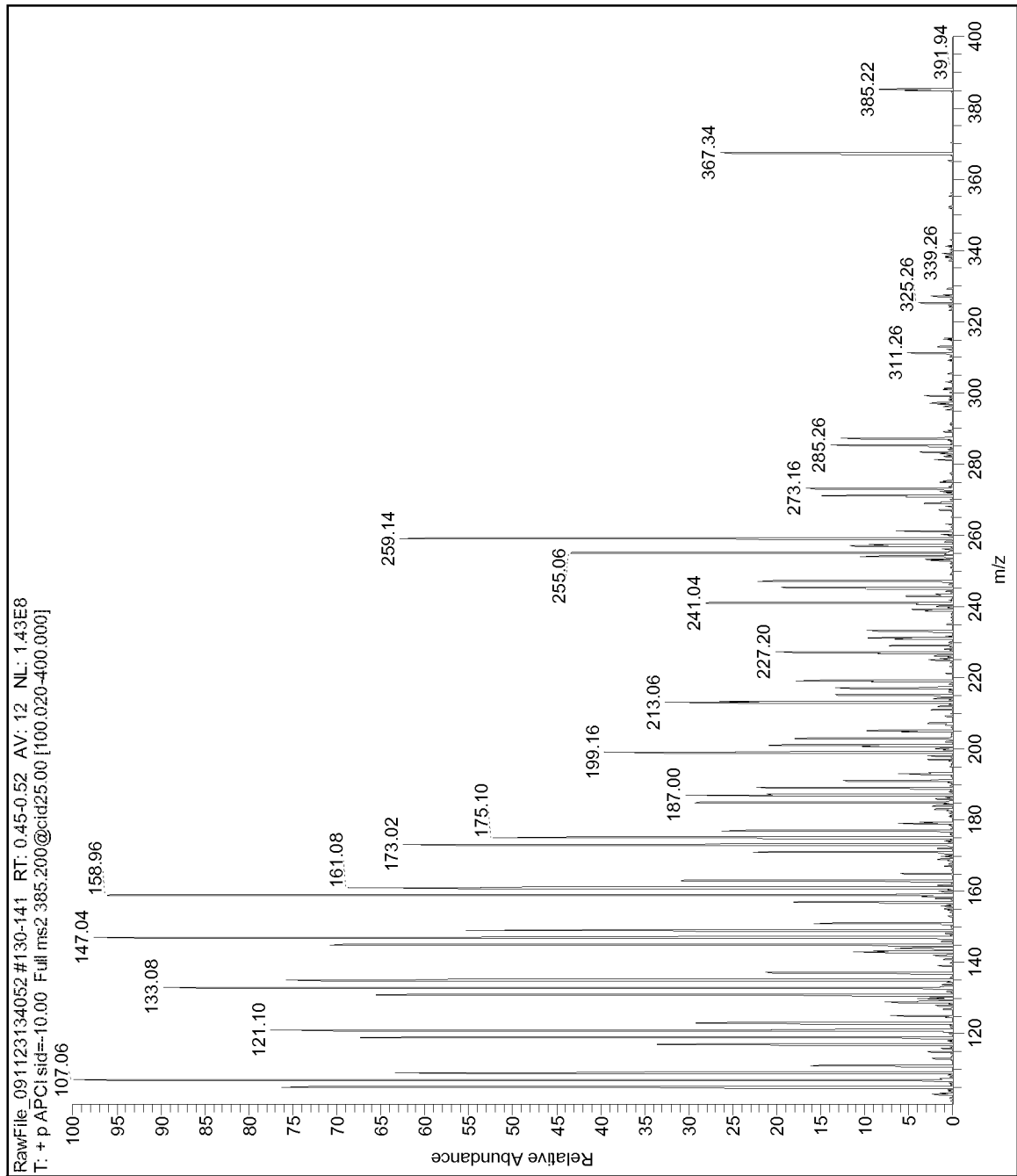


17/35

Figure 8C

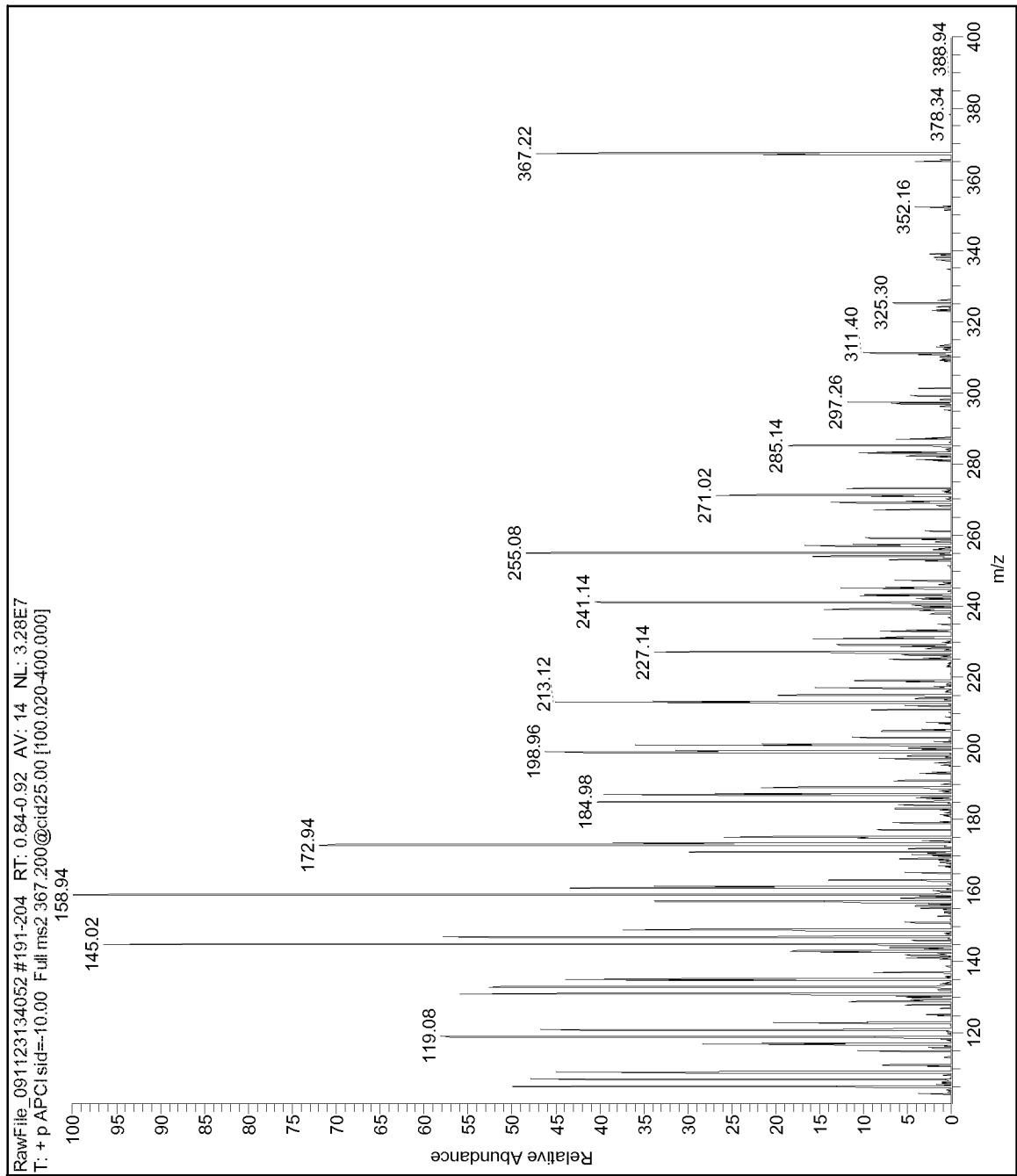


18/35

Figure 9A

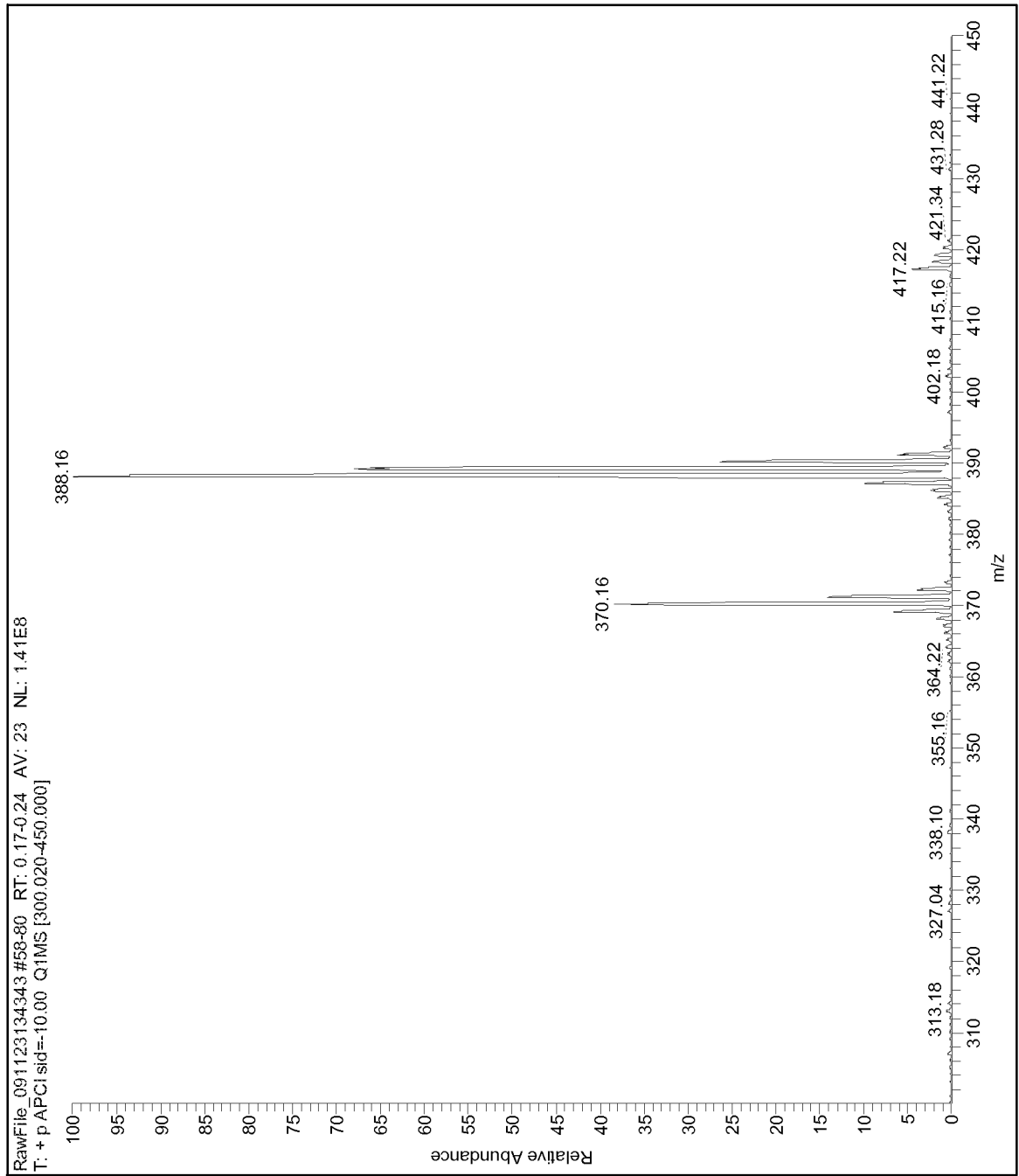


19/35

Figure 9B

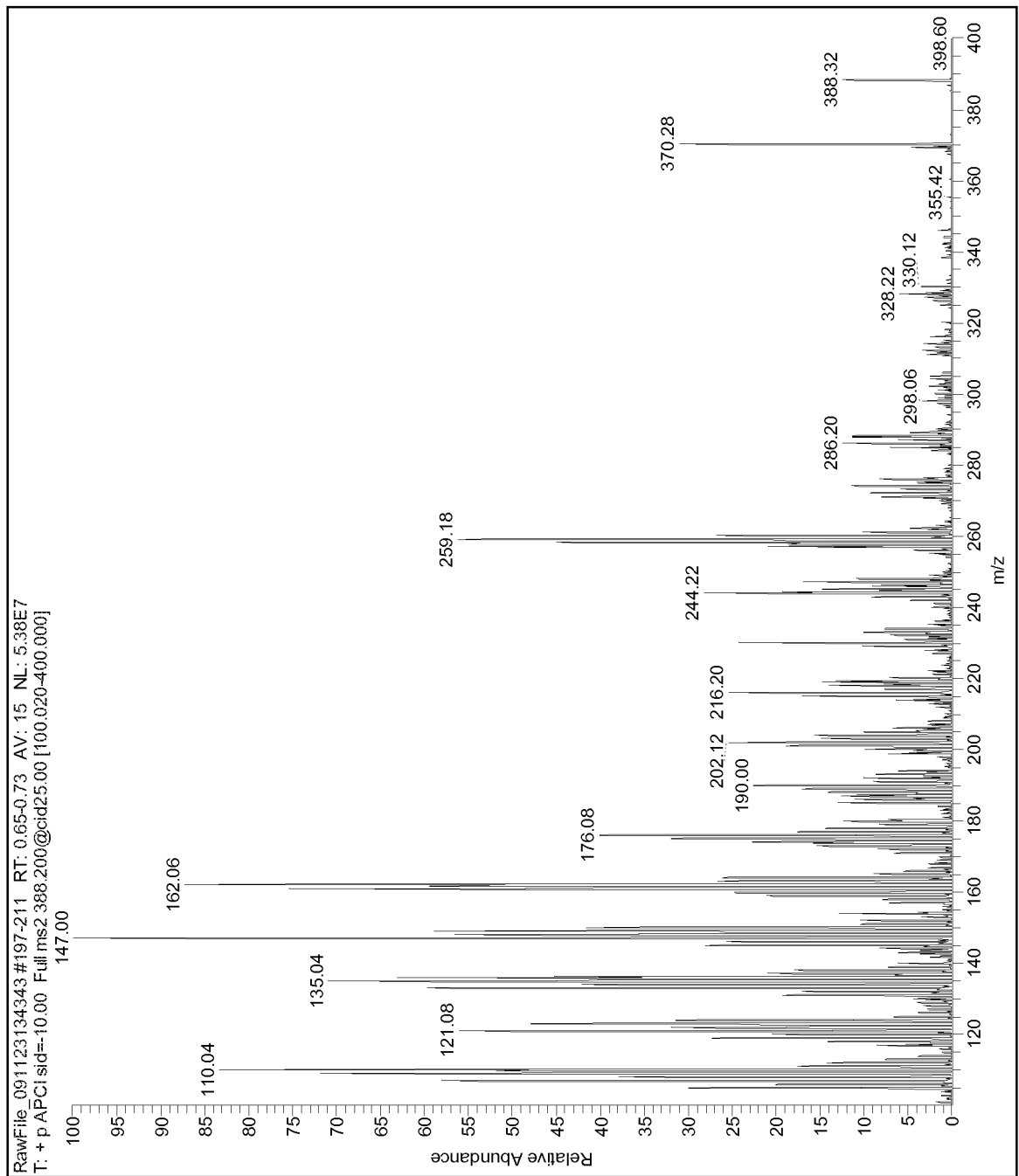


20/35

Figure 9C

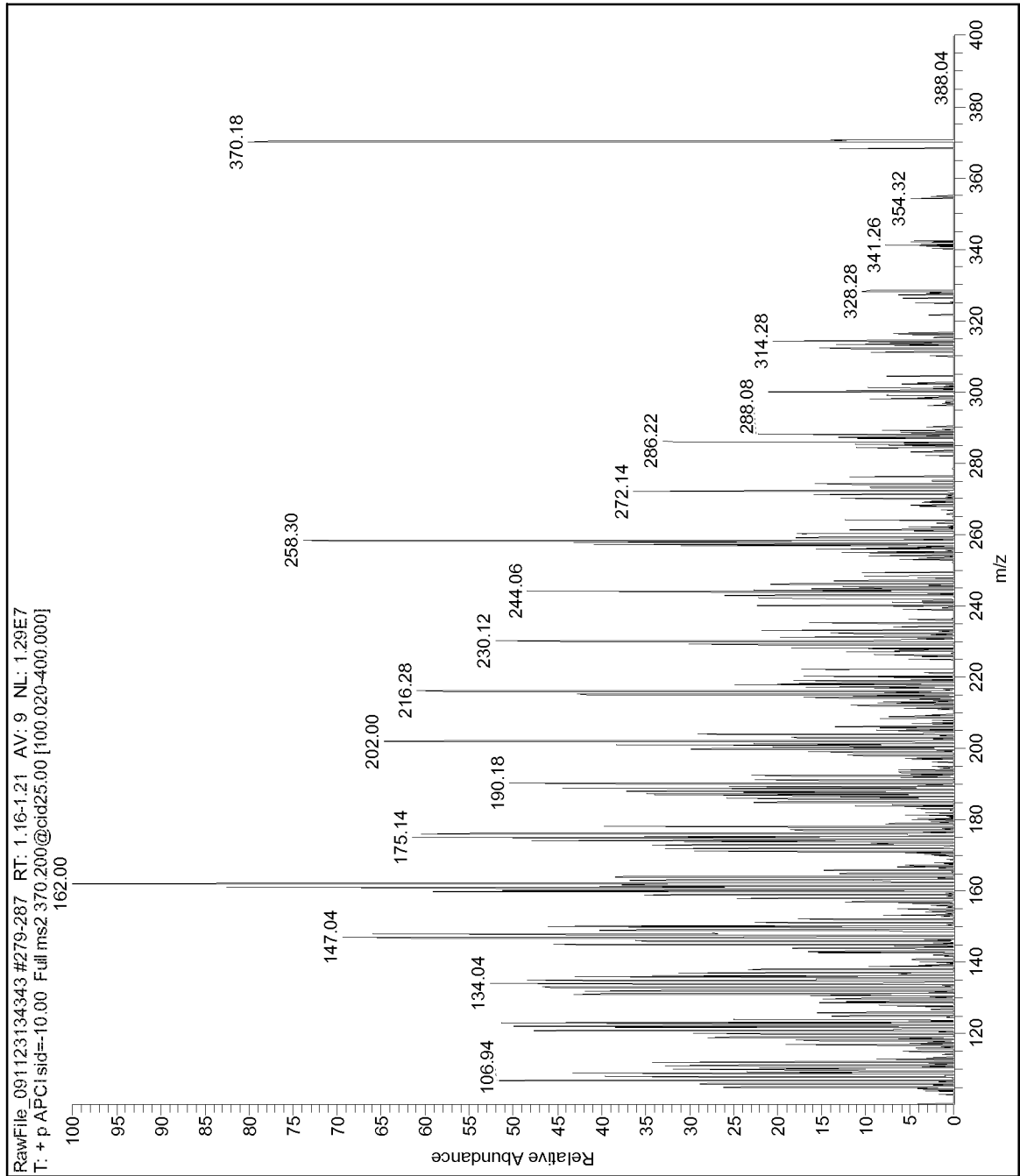


21/35
Figure 10A

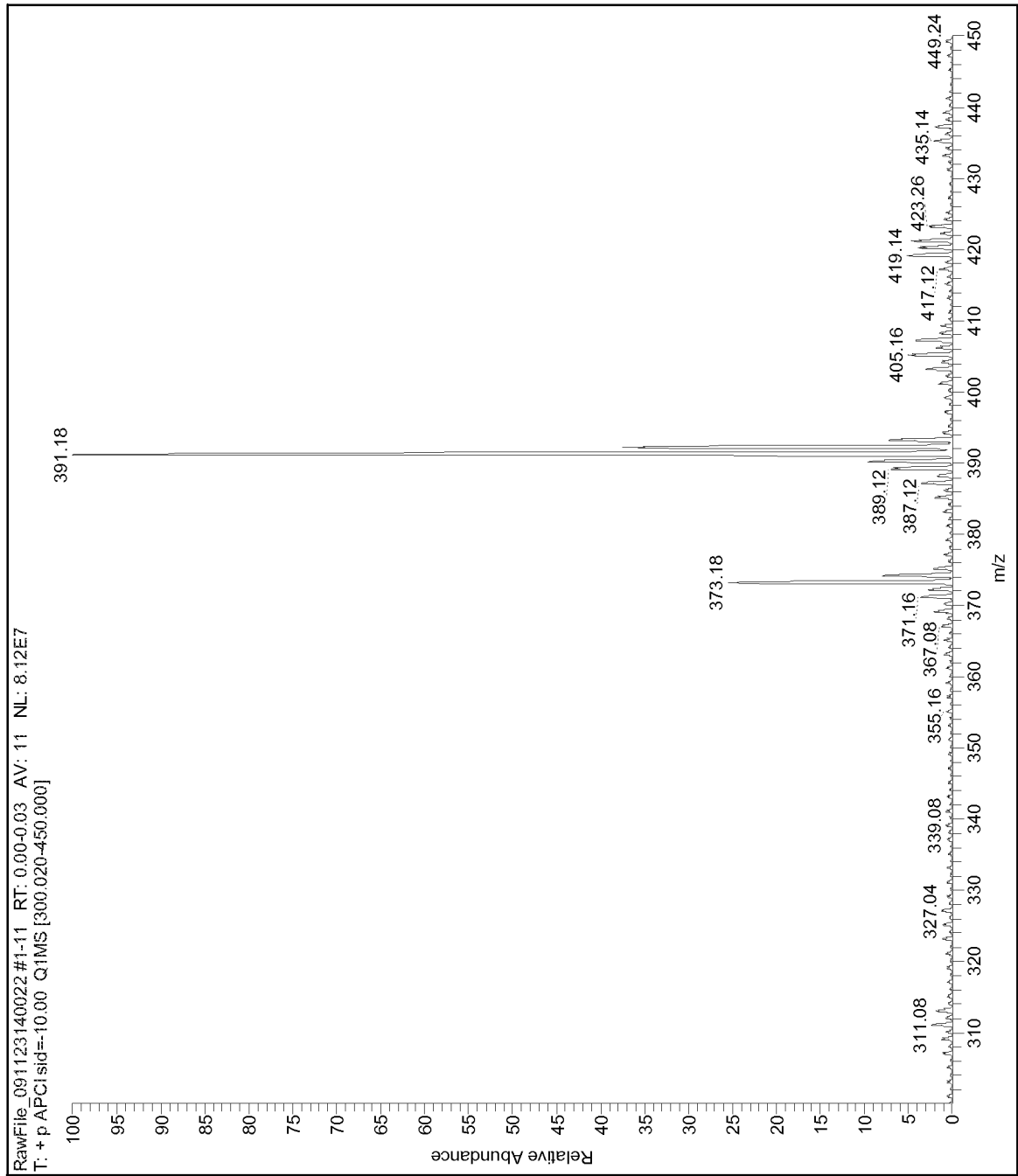


22/35

Figure 10B

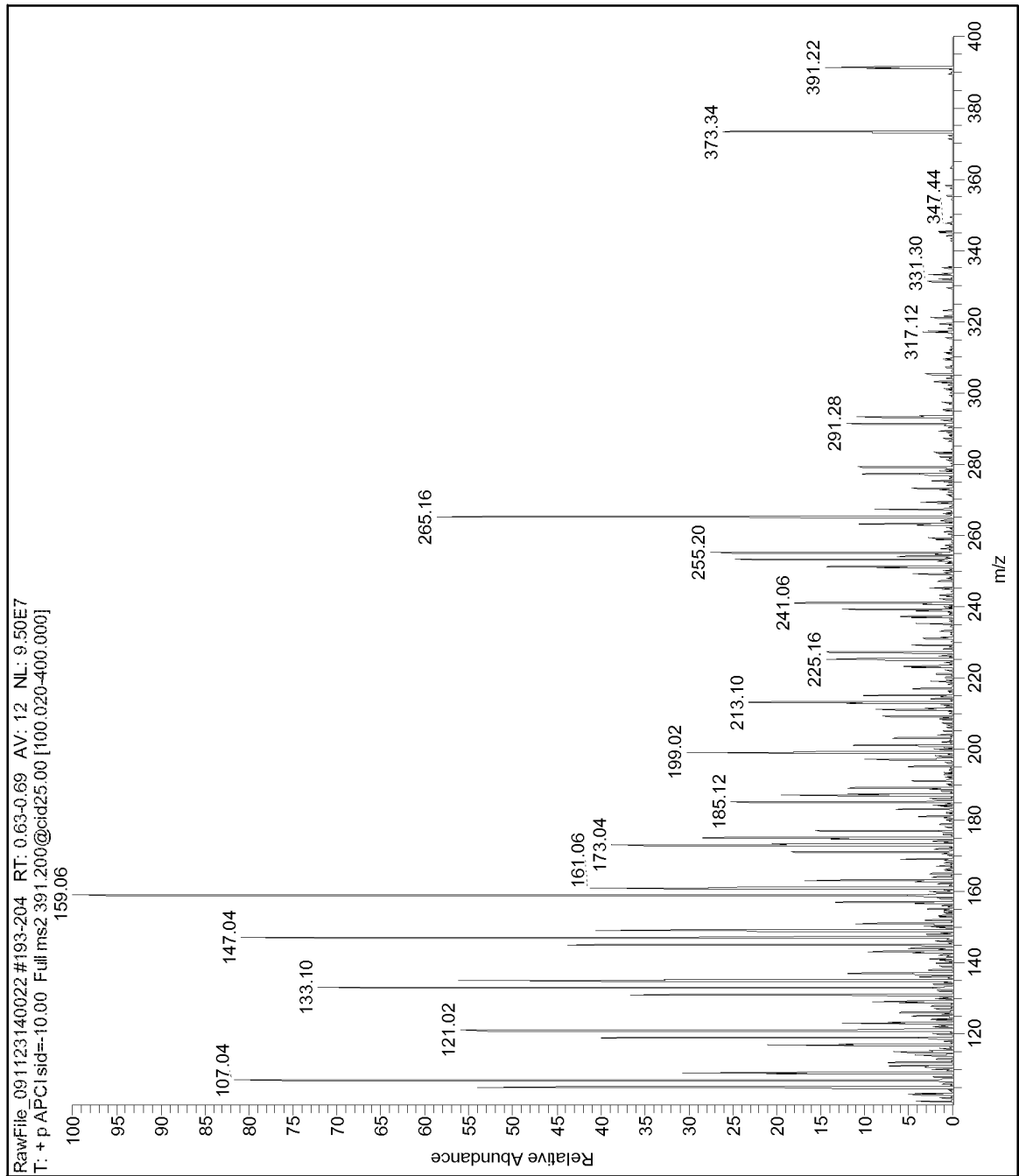


23/35

Figure 10C

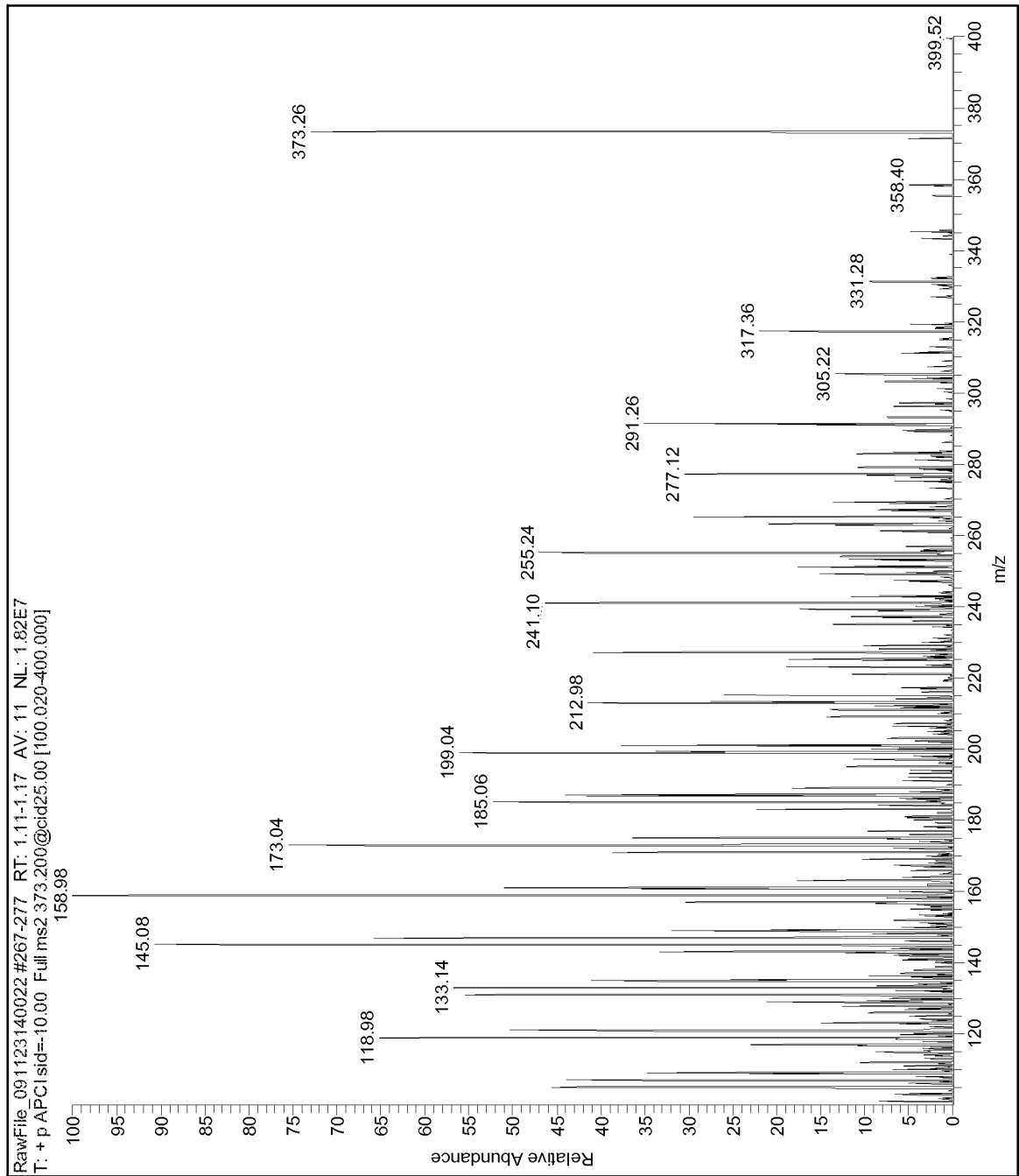


24/35

Figure 11A

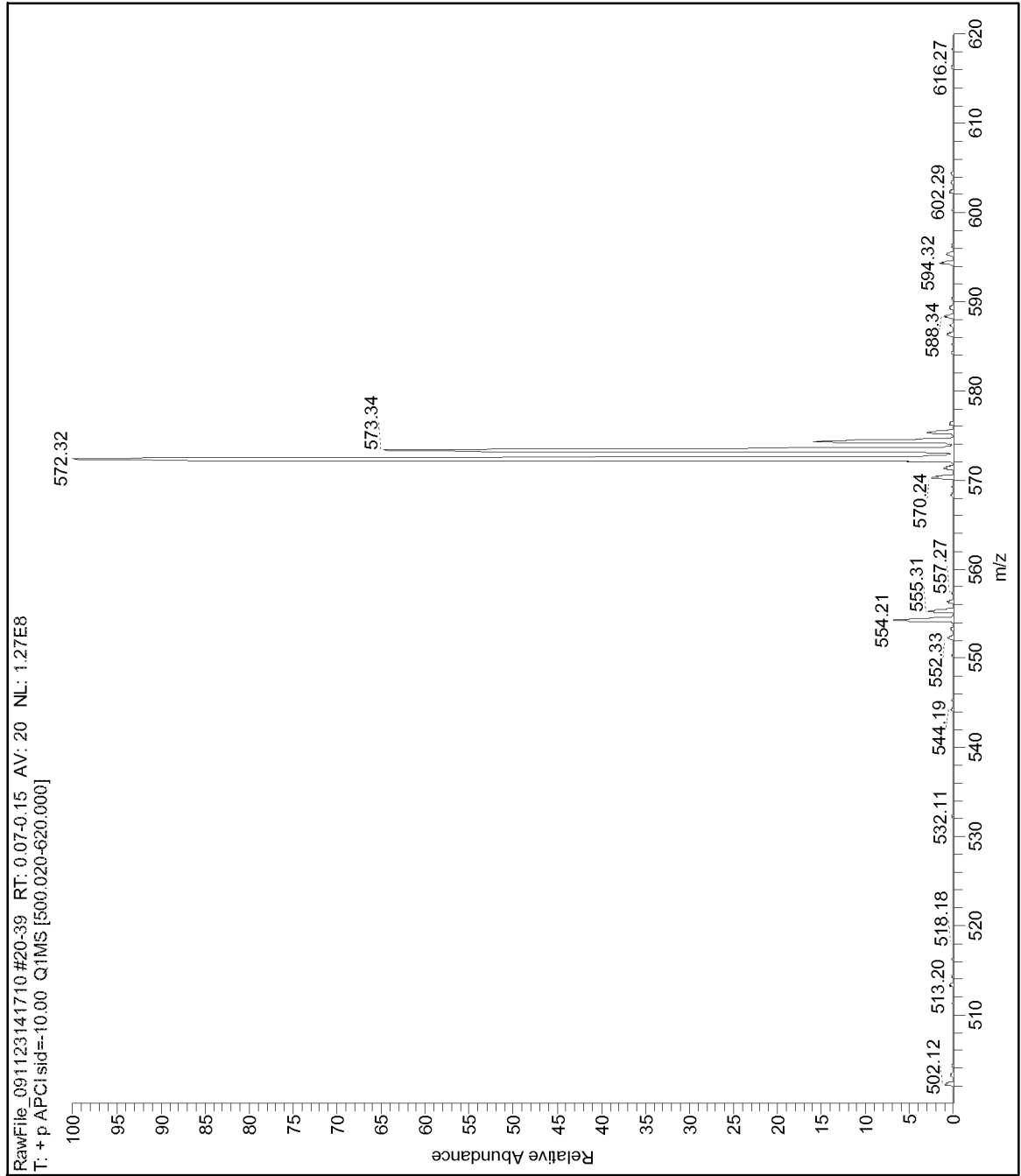


25/35

Figure 11B

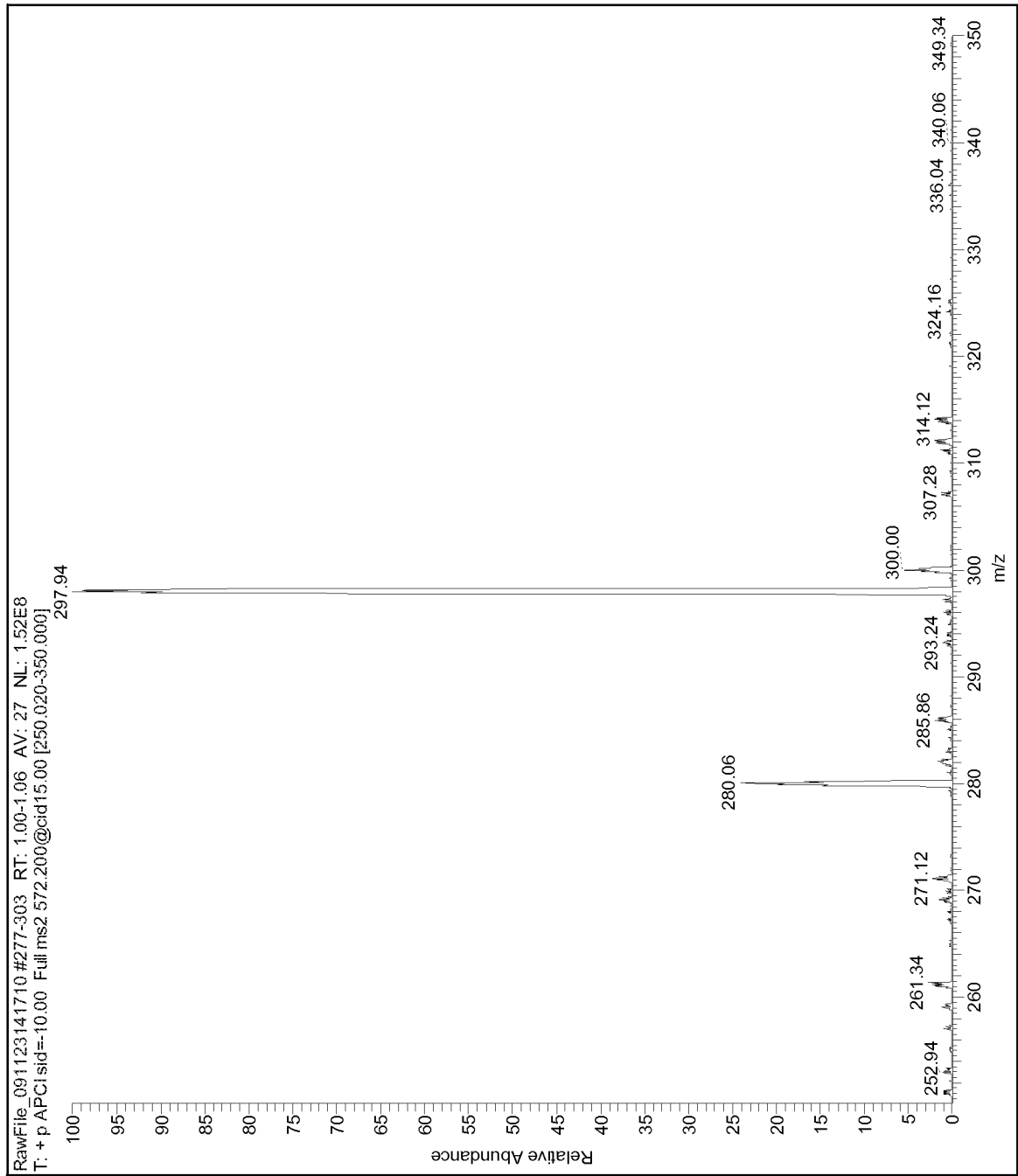


26/35

Figure 12A

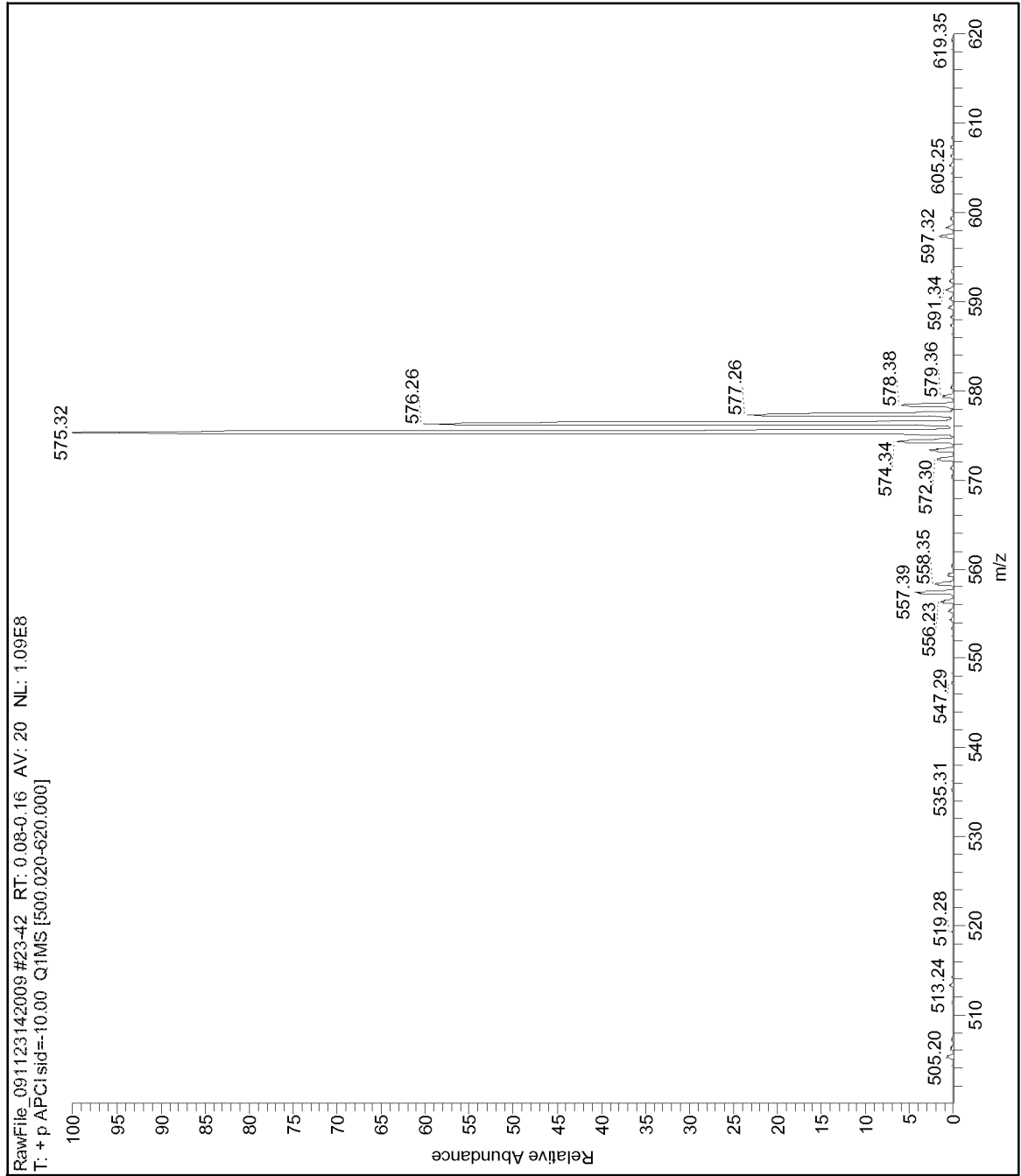


27/35
Figure 12B

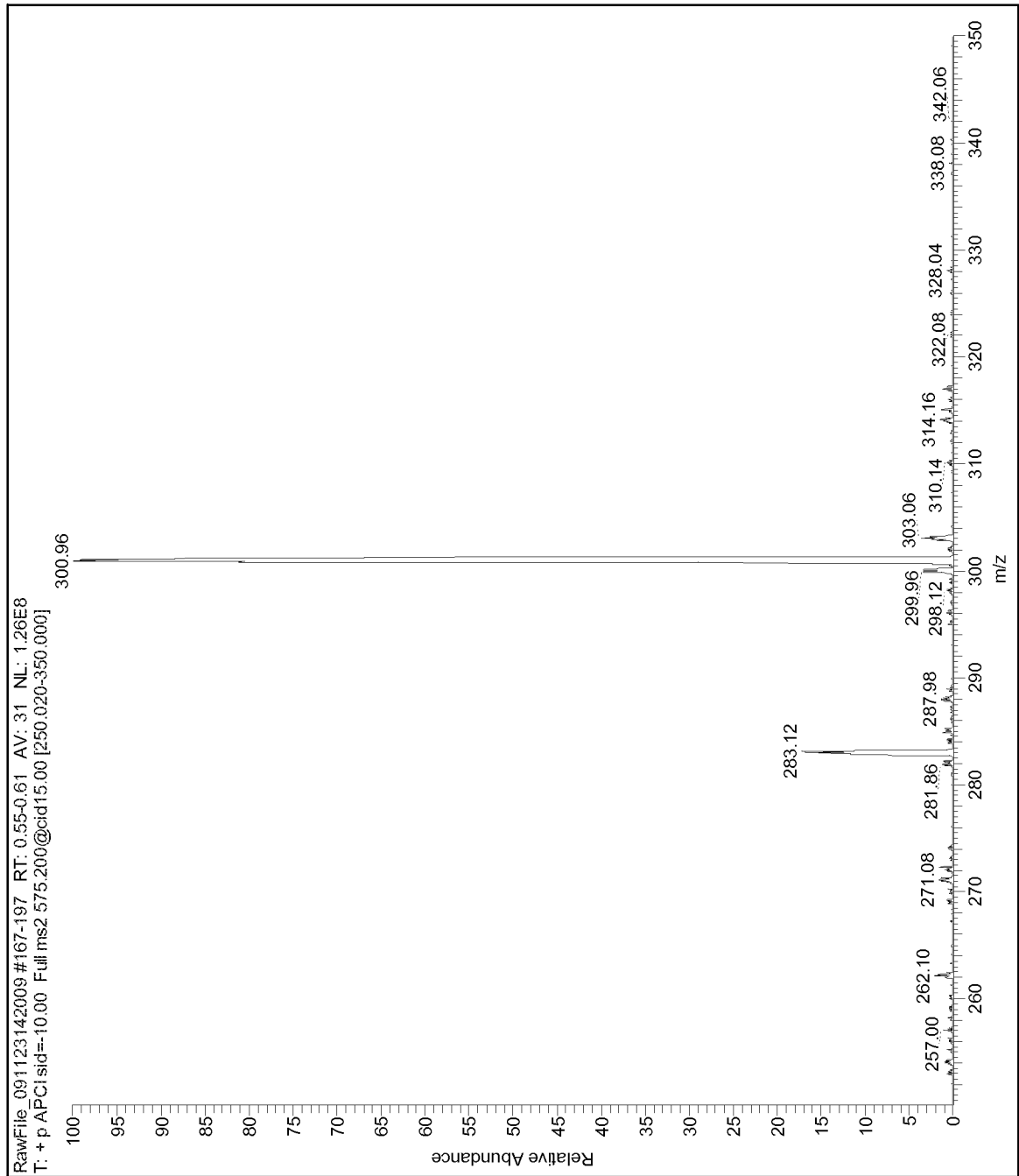


28/35

Figure 13A

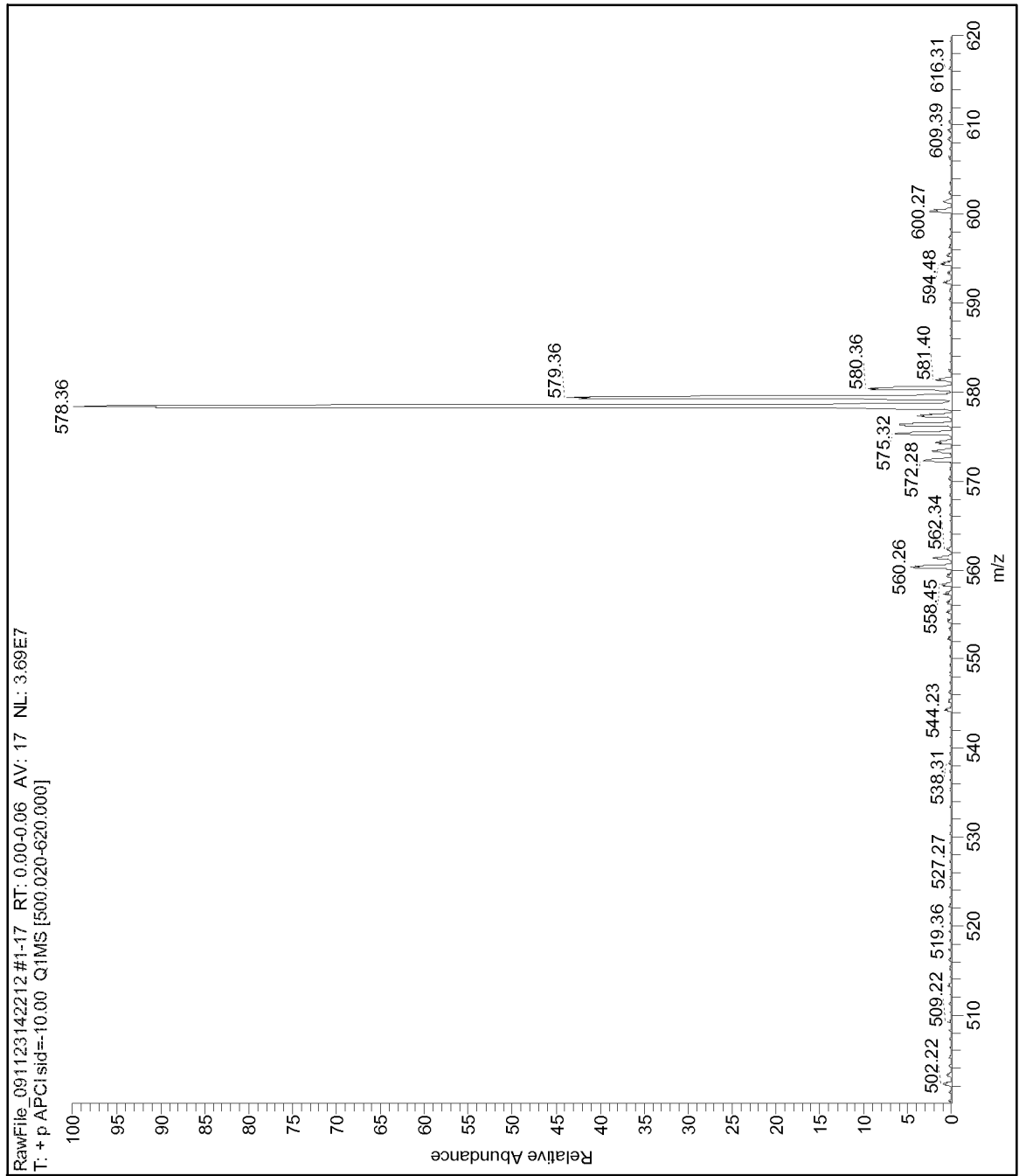


29/35

Figure 13B

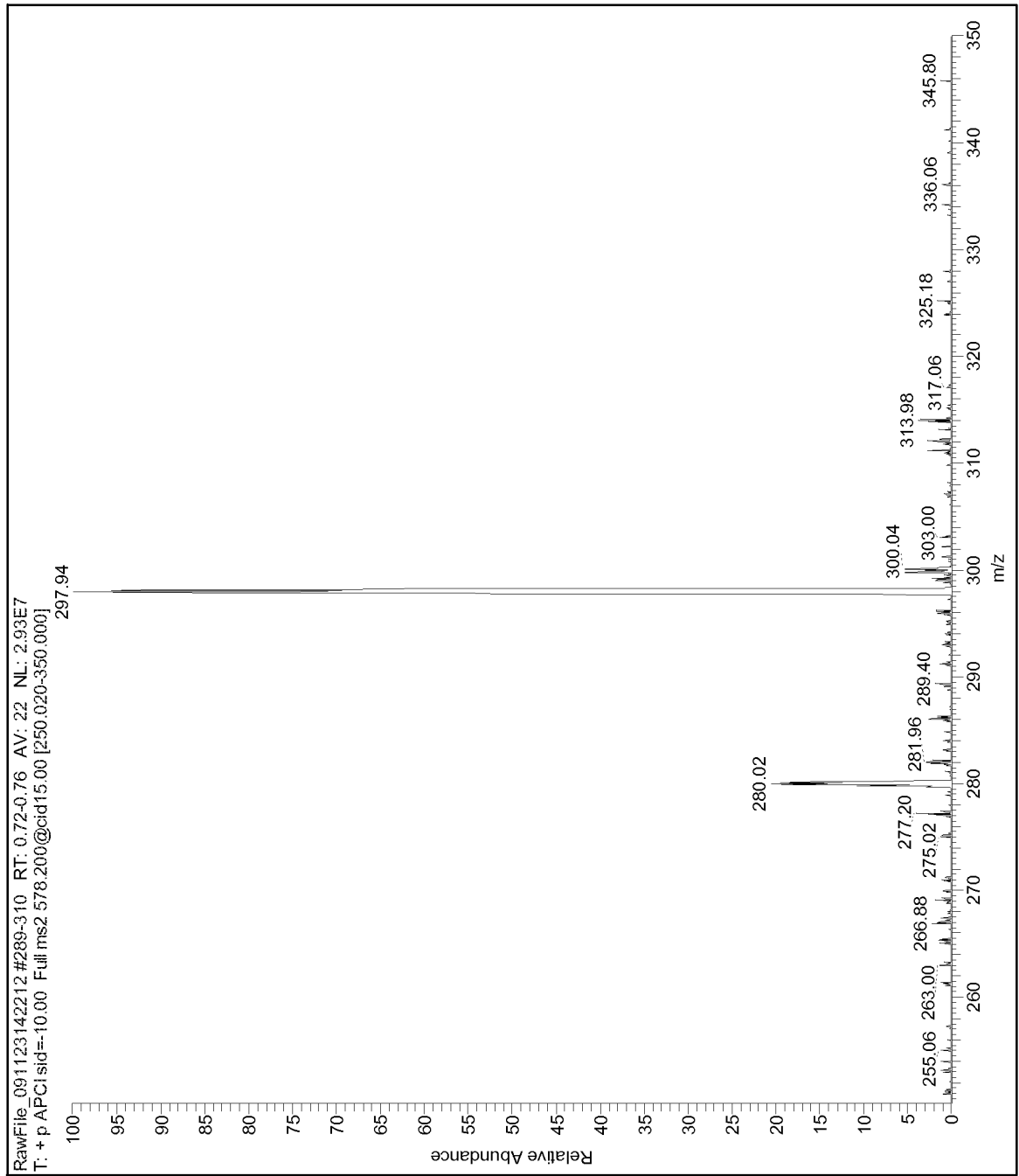


30/35

Figure 14A

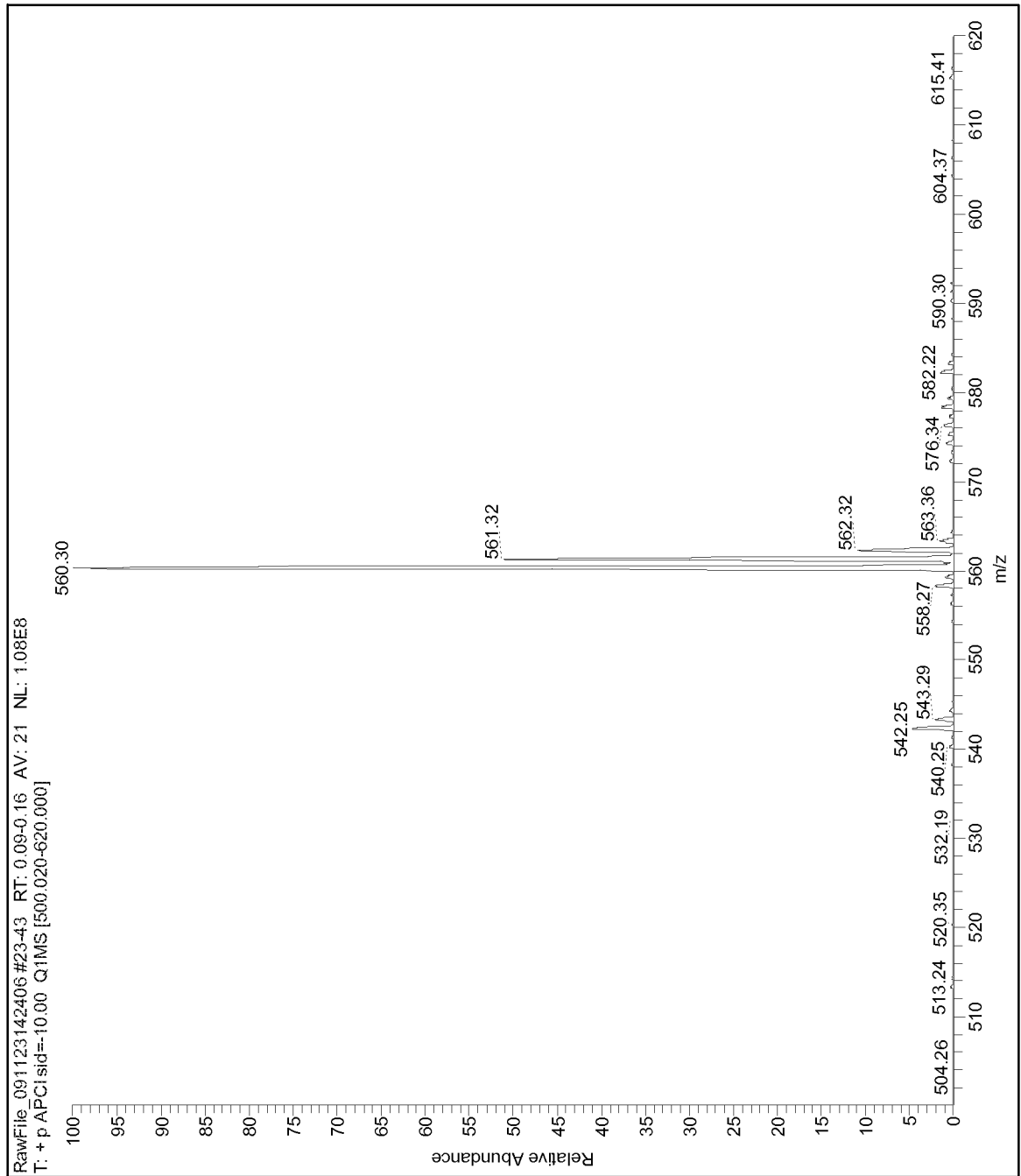


31/35
Figure 14B

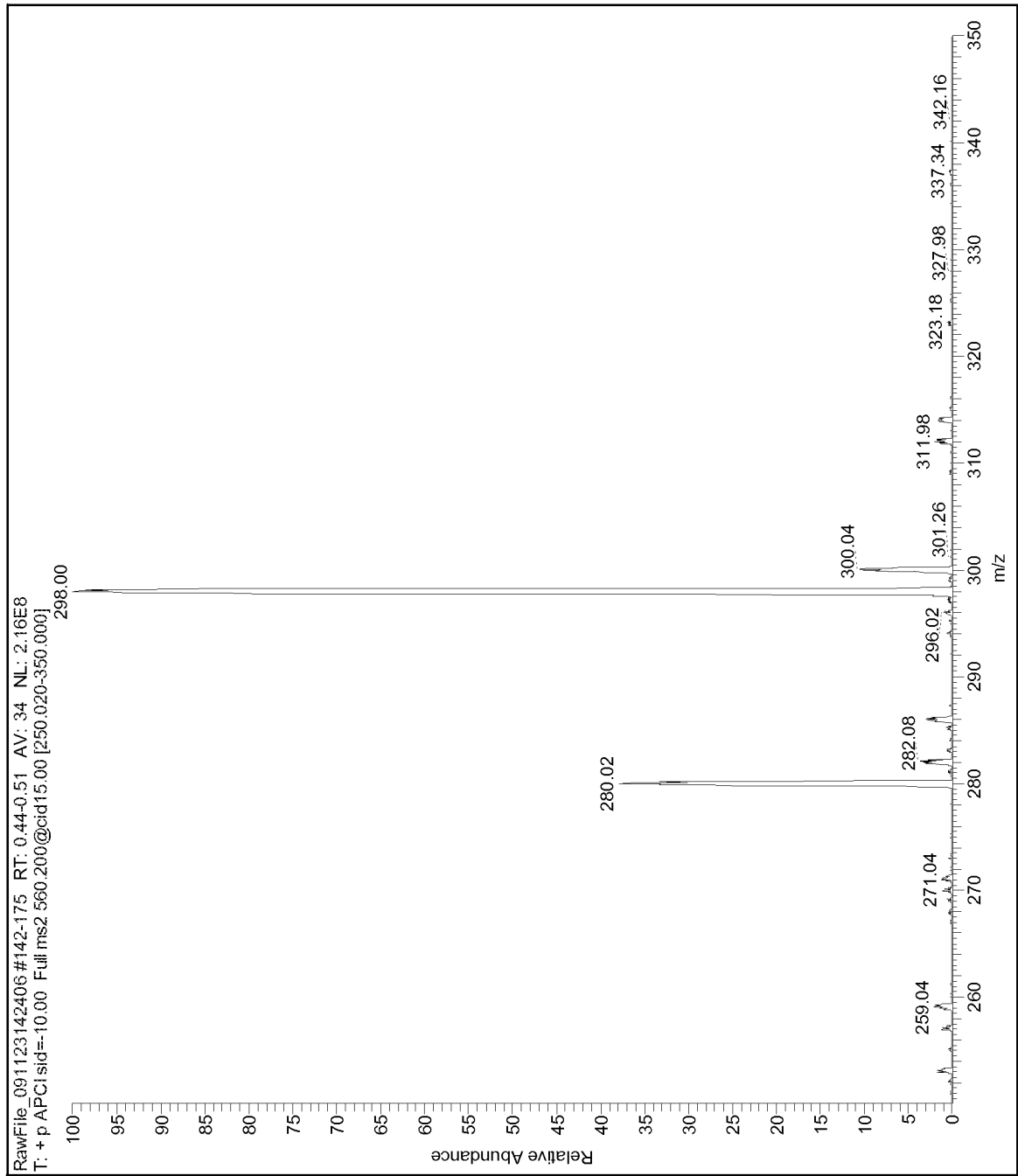


32/35

Figure 15A

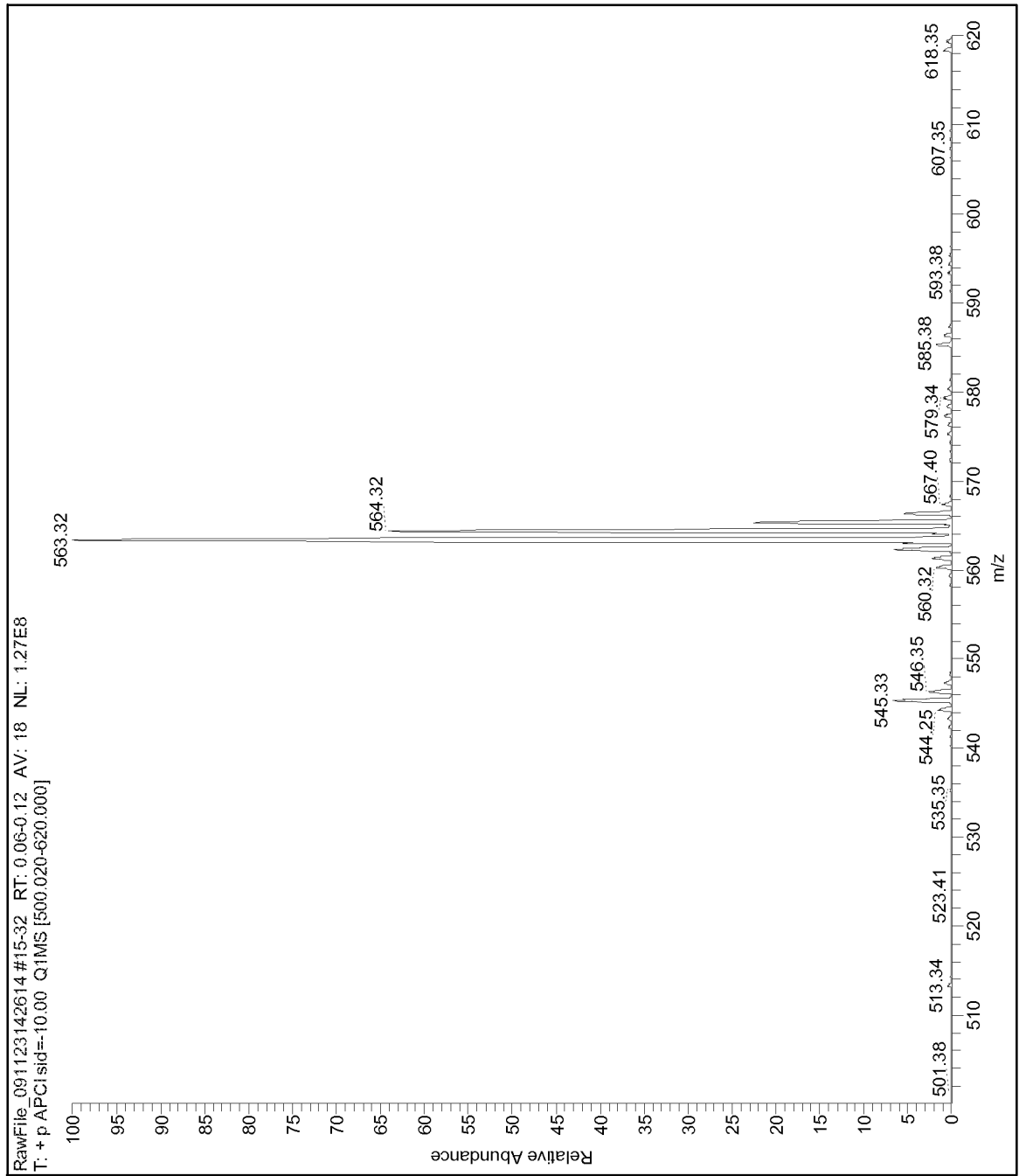


33/35

Figure 15B

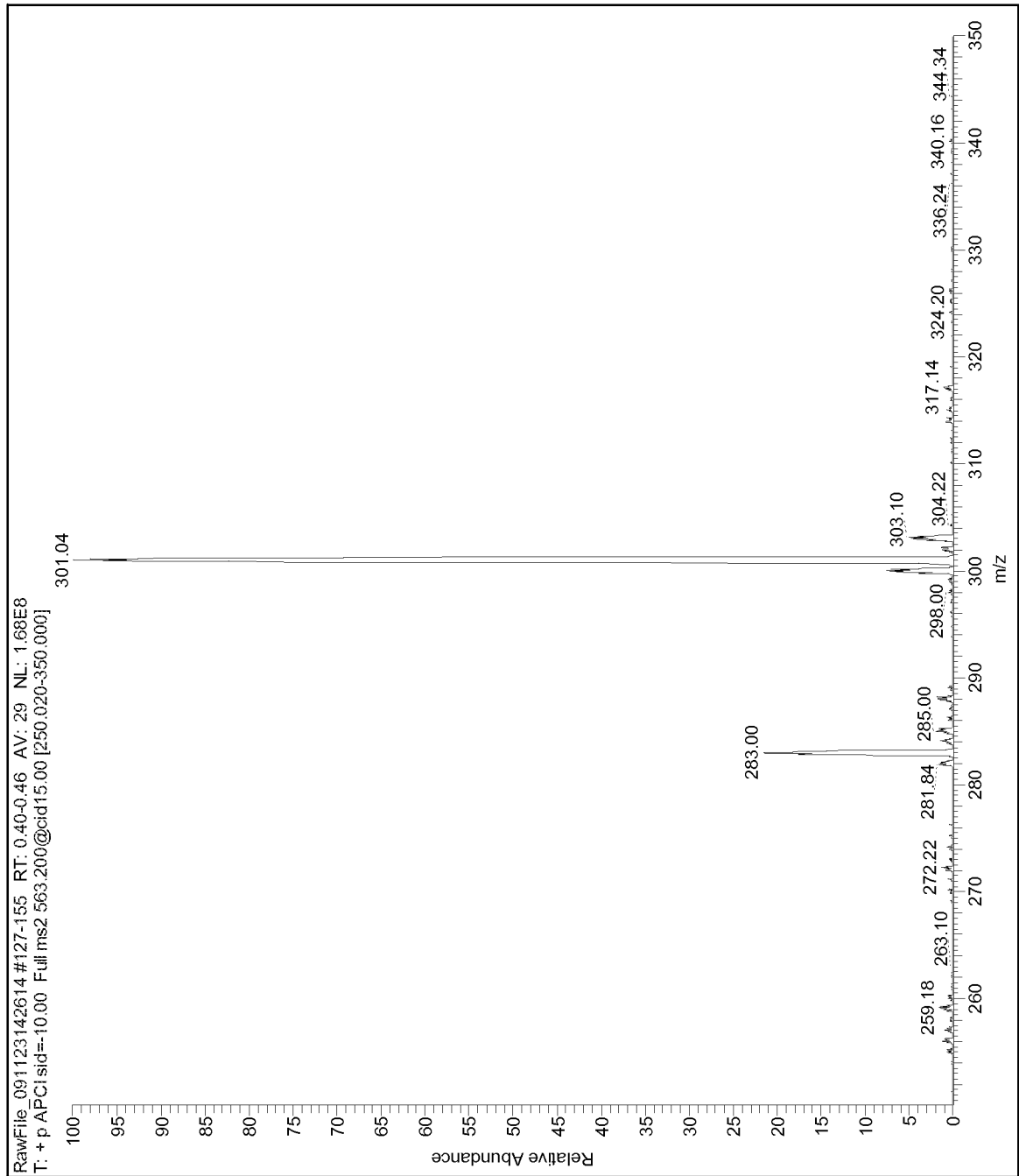


34/35

Figure 16A

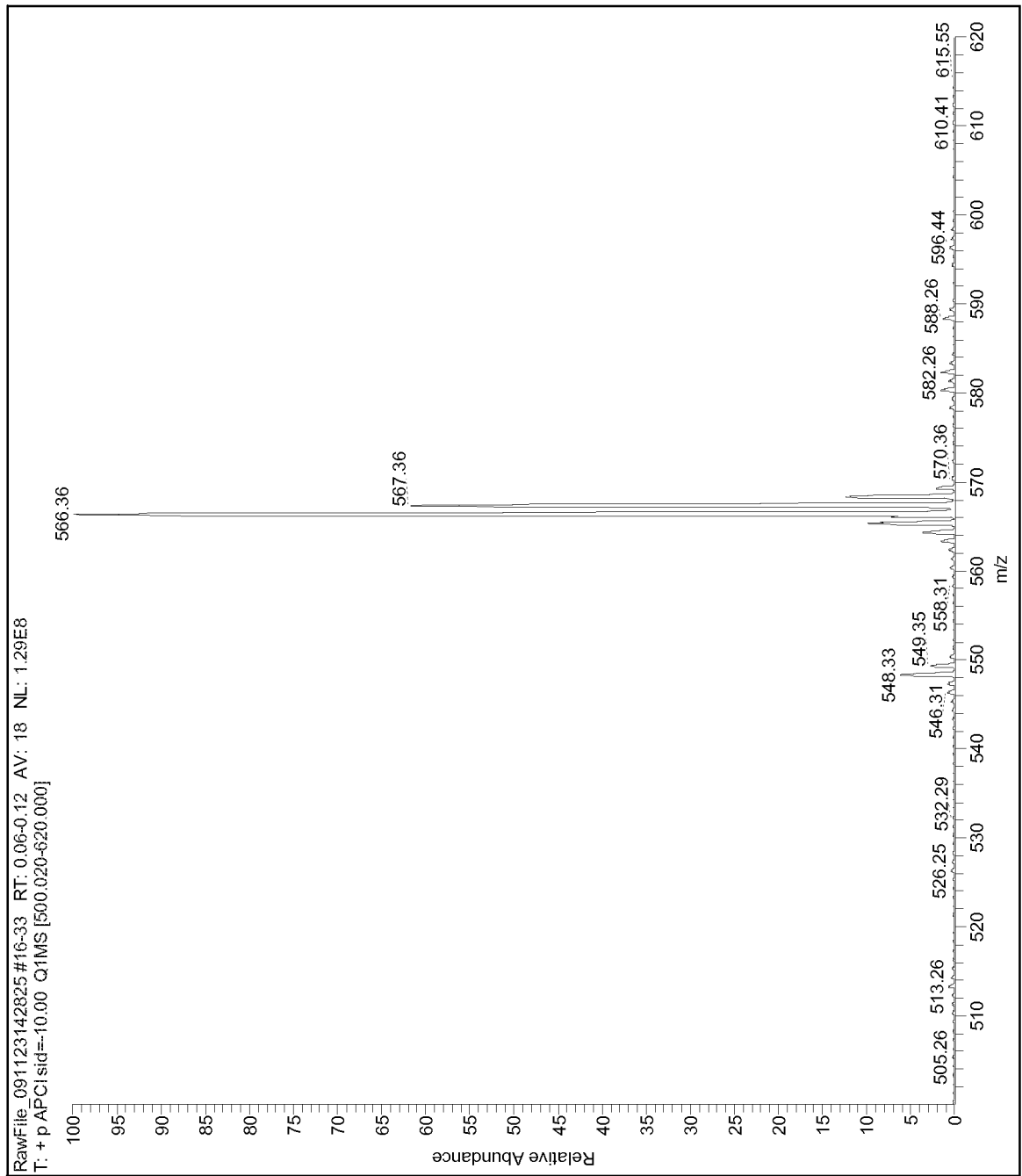

35/35

Figure 16B

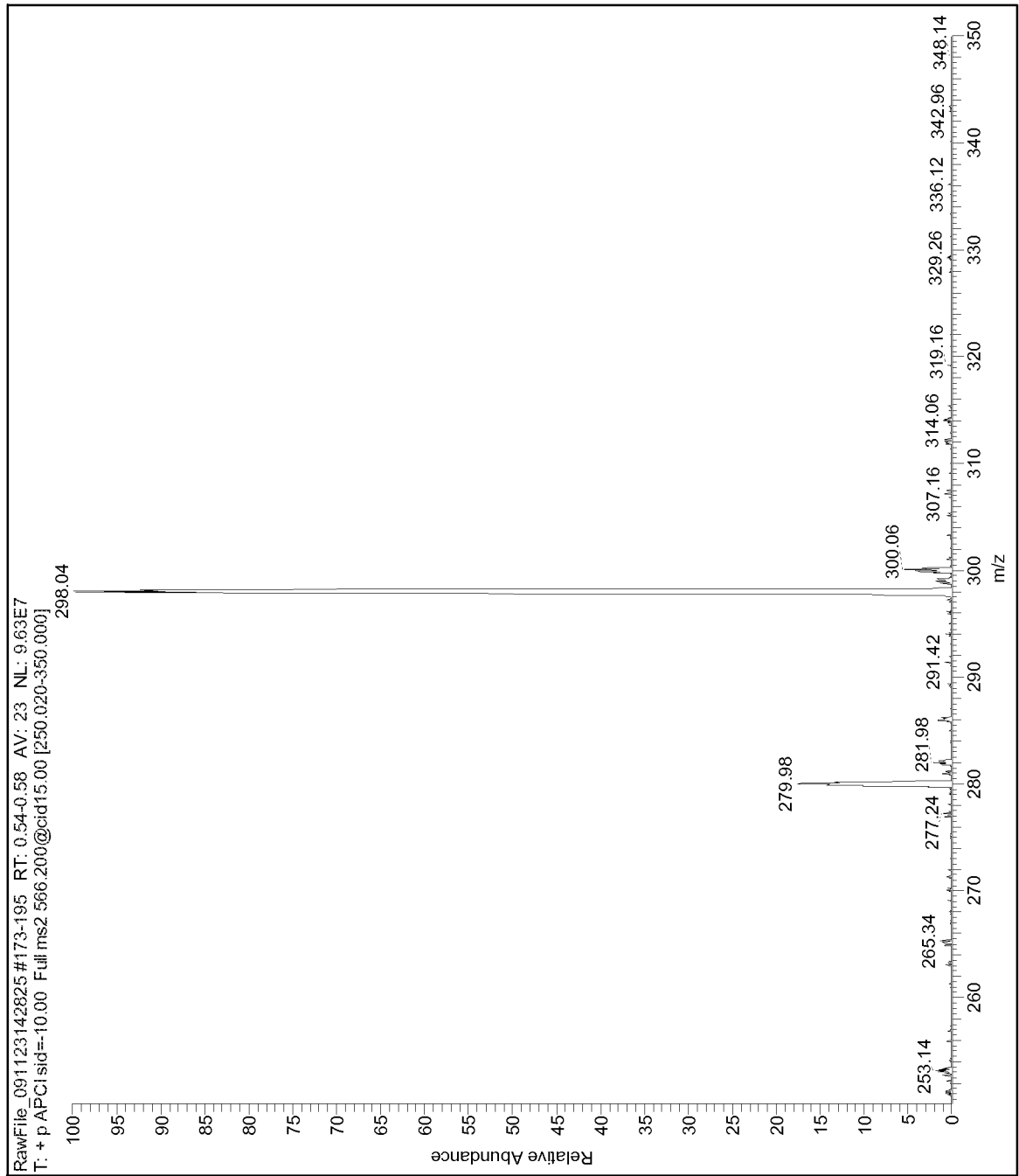

INTERNATIONAL SEARCH REPORT

International application No.

PCT/US 10/59771

A. CLASSIFICATION OF SUBJECT MATTER
 IPC(8) - C12Q 1/68; G06F 19/00 (2011.01)
 USPC - 435/6; 436/86; 702/19
 According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)
 IPC(8) - C12Q 1/68; G06F 19/00 (2011.01)
 USPC -- 435/6; 436/86; 702/19

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched
 IPC(8) - C12Q 1/68; G06F 19/00 (2011.01)
 USPC -- 435/6; 436/86; 702/19 (keyword delimited)

Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)
 PubWEST (USPT, PGPB, EPAB, JPAB), Google
 Search terms used: vitamin D, mass spectrometry, ionization, laser, volatile, heat, liquid chromatography, sample, assay, detect, serum, plasma

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
Y	US 2009/0137056 A1 (Holmquist et al.) 28 May 2009 (28.05.2009), entire document especially Abstract, para [0059], [0009], [0019], [0056], [0057]	1-5, 24-27
Y	US 2006/0228809 A1 (Clarke et al.) 12 October 2006 (12.10.2006), entire document especially Abstract, para [0040]-[0042]	1-5, 24-27
Y	US 2003/0171605 A1 (Reddy et al.) 11 September 2003 (11.09.2003), entire document especially para [0167]	1-5, 24-27

Further documents are listed in the continuation of Box C.

* Special categories of cited documents:

"A" document defining the general state of the art which is not considered to be of particular relevance	"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
"E" earlier application or patent but published on or after the international filing date	"X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)	"Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art
"O" document referring to an oral disclosure, use, exhibition or other means	"&" document member of the same patent family
"P" document published prior to the international filing date but later than the priority date claimed	

Date of the actual completion of the international search 27 January 2011 (27.01.2011)	Date of mailing of the international search report 11 FEB 2011
---	--

Name and mailing address of the ISA/US Mail Stop PCT, Attn: ISA/US, Commissioner for Patents P.O. Box 1450, Alexandria, Virginia 22313-1450 Facsimile No. 571-273-3201	Authorized officer: Lee W. Young PCT Helpdesk: 571-272-4300 PCT OSP: 571-272-7774
---	--

INTERNATIONAL SEARCH REPORT

International application No.

PCT/US 10/59771

Box No. II Observations where certain claims were found unsearchable (Continuation of item 2 of first sheet)

This international search report has not been established in respect of certain claims under Article 17(2)(a) for the following reasons:

1. Claims Nos.:
because they relate to subject matter not required to be searched by this Authority, namely:

2. Claims Nos.:
because they relate to parts of the international application that do not comply with the prescribed requirements to such an extent that no meaningful international search can be carried out, specifically:

3. Claims Nos.: 6-23, 28-44
because they are dependent claims and are not drafted in accordance with the second and third sentences of Rule 6.4(a).

Box No. III Observations where unity of invention is lacking (Continuation of item 3 of first sheet)

This International Searching Authority found multiple inventions in this international application, as follows:

1. As all required additional search fees were timely paid by the applicant, this international search report covers all searchable claims.
2. As all searchable claims could be searched without effort justifying additional fees, this Authority did not invite payment of additional fees.
3. As only some of the required additional search fees were timely paid by the applicant, this international search report covers only those claims for which fees were paid, specifically claims Nos.:

4. No required additional search fees were timely paid by the applicant. Consequently, this international search report is restricted to the invention first mentioned in the claims; it is covered by claims Nos.:

Remark on Protest

- The additional search fees were accompanied by the applicant's protest and, where applicable, the payment of a protest fee.
- The additional search fees were accompanied by the applicant's protest but the applicable protest fee was not paid within the time limit specified in the invitation.
- No protest accompanied the payment of additional search fees.